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EN - INSTRUCTIONS FOR USE

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KYRI CART

Art. CA001TYG

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Carrello mobile riunito

TYGI 350 B

CONTENUTO

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TYGI 350 B

CONTENTS

1	PRODUCT DESCRIPTION	3
2	MAIN DIMENSIONS.....	4
3	MAIN CONTROLS	5
4	TECHNICAL DATA	7
4.1	Unit Label	7
5	PRE-INSTALLATION REQUIREMENTS.....	8
5.1	Environmental conditions	8
5.2	Requirements for the installation of media.....	8
5.3	Floor.....	8
5.4	Environment.....	8
6	ASSEMBLY AND INSTALLATION	8
7	PUTTING THE UNIT INTO OPERATION.....	9
8	PRODUCT OPERATION.....	10
8.1	Control panel with instruments	10
8.1.1	<i>Main screen.....</i>	<i>10</i>
8.1.2	<i>System of changing parameters in individual screens.....</i>	<i>13</i>
8.1.3	<i>Alarm indication</i>	<i>22</i>
8.1.4	<i>Operation of individual instruments of control panel.....</i>	<i>23</i>
8.1.5	<i>Syringe</i>	<i>23</i>
8.1.6	<i>Turbine.....</i>	<i>23</i>
8.1.7	<i>Micromotor BLDC - DX, DX BLUE, DX PRO.....</i>	<i>27</i>
8.1.8	<i>Surgical micromotor DX SRG with peristaltic pump.....</i>	<i>31</i>
8.1.9	<i>Scaler</i>	<i>34</i>
8.2	Foot controller.....	37
8.3	Bottle with distilled water	41
8.4	Instruments cooling - setting of the water amount.....	41
8.5	Adjusting the height of the unit	42
8.6	Completion of work.....	42
9	PRODUCT MAINTENANCE	43
10	CLEANING AND DISINFECTION.....	44
10.1	Elementary disinfection of the cooling waterlines of instruments.....	44
10.2	Cleaning and disinfection of other parts of the dental unit	44
11	DISPOSAL OF THE EQUIPMENT	45
12	REPAIR SERVICE	45
13	CONTENTS OF THE PACKAGING	46
14	WARRANTY.....	46
15	TRANSPORTING	46
16	STORAGE	46
17	REQUIREMENTS ON ELECTROMAGNETIC COMPATIBILITY ACCORDING TO EN 60601-1-2	47

These Instructions for use describe how to use the dental unit **TYGI 350 B**. Please read these Instructions for use thoroughly before using them.

1 PRODUCT DESCRIPTION

The dental unit TYGI 350 B is a special mobile cart version of dentist's control panel with lower hose delivery and with adjustable height.

The use of mobile carts allows quick set up of supplies and equipment used in endodontics, oral surgery and implant placement.

The dentist's element is on wheels. Twin wheel casters are suitable even for hard floor.

The capacitive touch screen display is built in the body of the control panel. Software of control panel touch screen is designed for effective working with a chair KYRI.

Instruments and the chair can be controlled via touch screen display or either by foot controller UNO or NOK (cabled or WIFI version).

Handles are used for moving the unit to desired place and either for changing of the control panel height position. To find your ideal operating position, you can adjust the height of the instrument table (200 mm height range).

On the top of the control panel is located portable stainless tray table plate with dimensions 180 x 280 mm.

Regulation of the cooling water for instruments is optionally based on order configuration:


- Manual water cooling regulation using one mechanical regulator for all instruments
- Electronic regulation of cooling water using proportional valve, which allows separate programming of each instrument

Self-contained water system utilizes 2-liter bottle for distilled water.




The panel can be equipped with 5 instruments. Typically with following instruments: syringe, turbines, micromotors, surgery micromotor, piezo scaler and curing lamp.

The control panel can be fitted with the following instruments:

Condition	Remark
maximum 5 instruments	
minimum 1 syringe	
maximum 4 rotating instruments	Possible combinatios: - max 3x turbines with light - max3x brushless micromotors DENSIM DX or DENSIM DX BLUE - 2x DX PRO + 1x DX SRG - 2x DX PRO BLUE + 1x DX SRG
1 piezo scaler	USS – Ultra Sonic Scaler
1 curing polymerization lamp	PLM

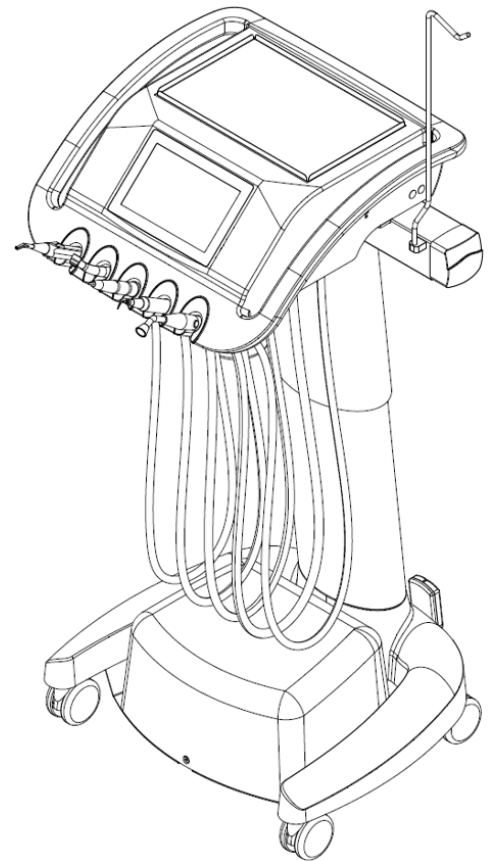
2 MAIN DIMENSIONS

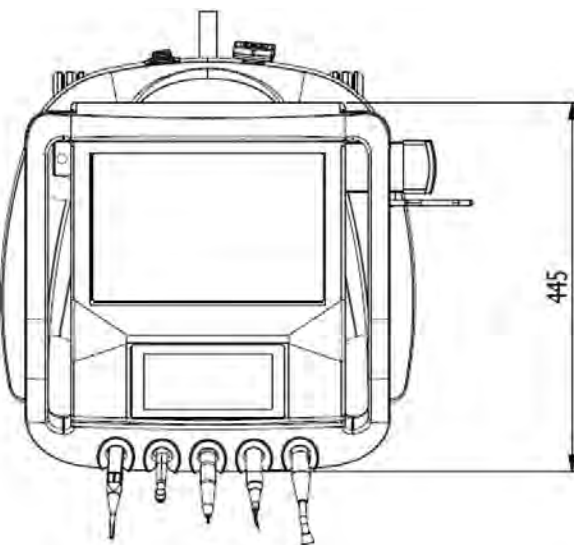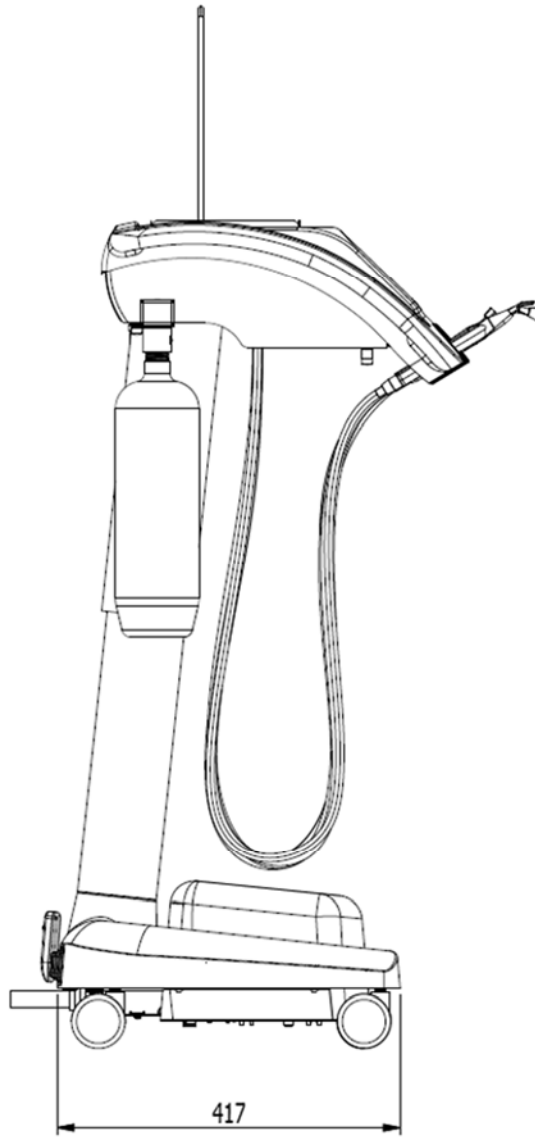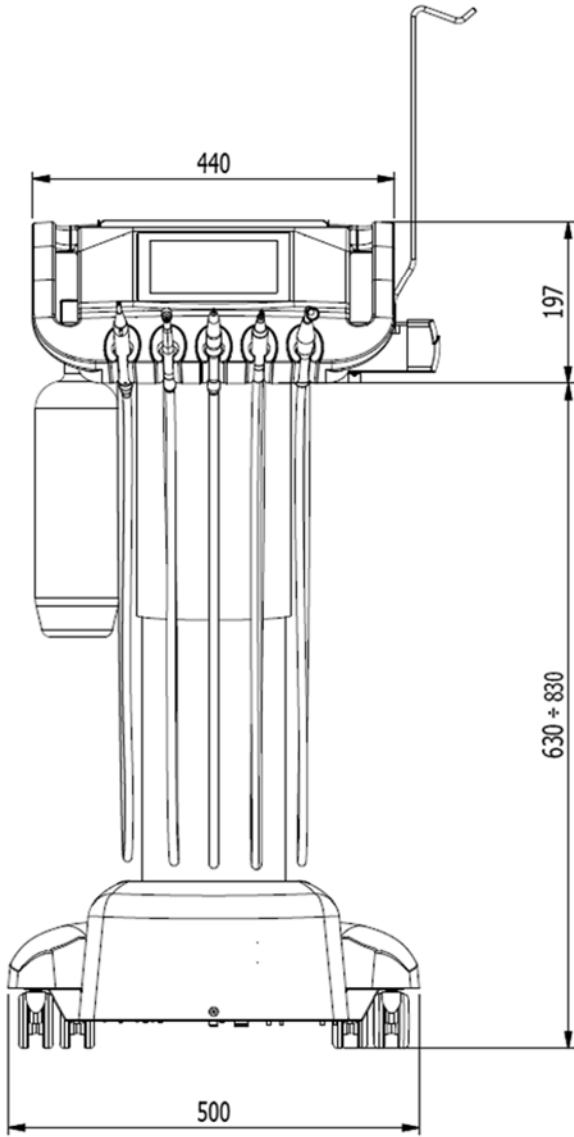3 MAIN CONTROLS

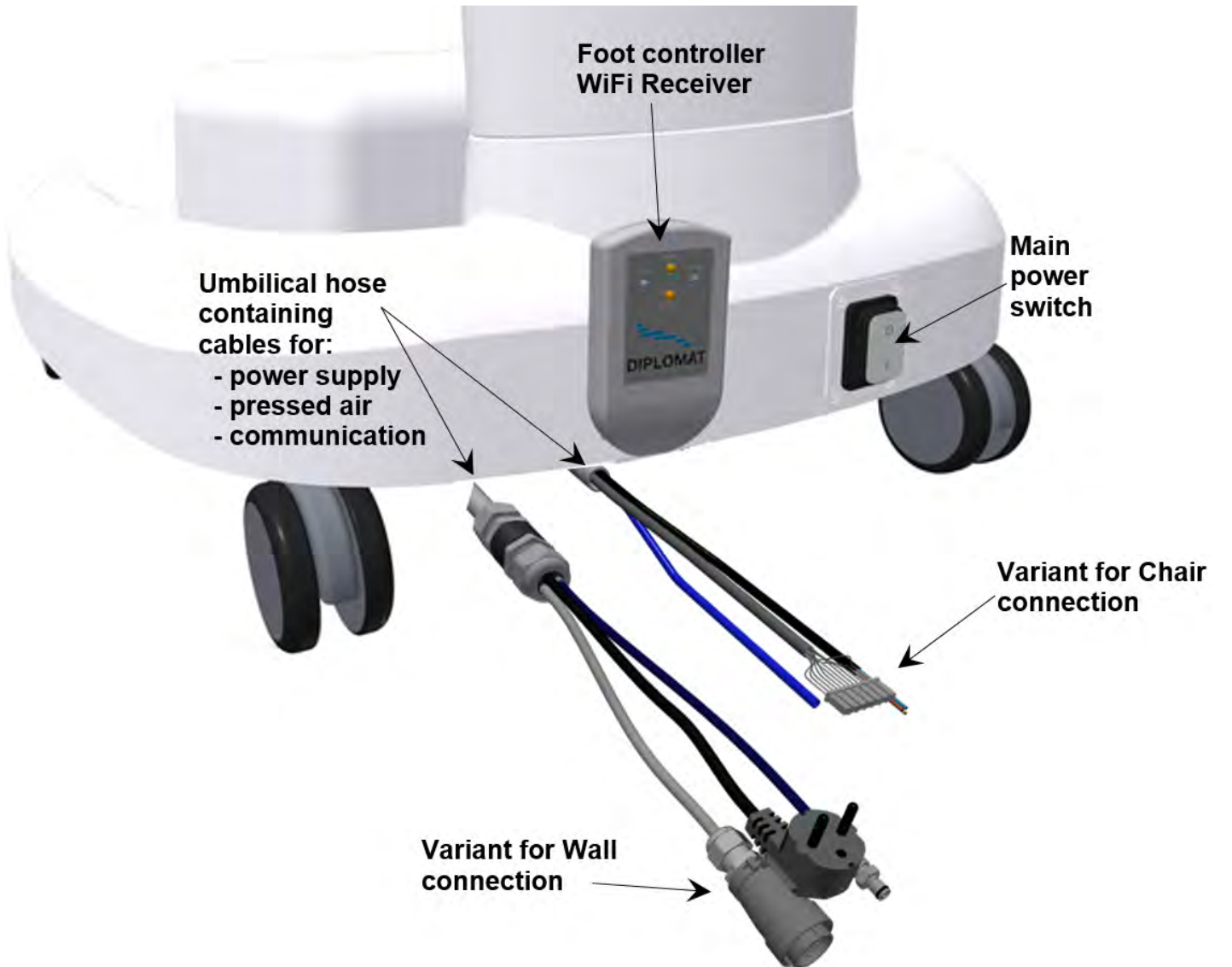
4 TECHNICAL DATA

Supply voltage	230V \pm 10%
Frequency	50 Hz \pm 2 %
Max. power consumption	150 VA
Input pressure of air	from 0,45 to 0,8 MPa
Unit weight netto	27 kg
Type of protection against electric shock	Class I equipment
Degree of protection against electric shock	applied parts of B type
Degree of protection by cover	IP21
Temperature of water for the instruments	35 \pm 5 °C (with heater fitted)
Max. loading capacity of the tray table	3 kg

Parts of the dental unit coming into contact with the patient:

- Instruments located on the control panel

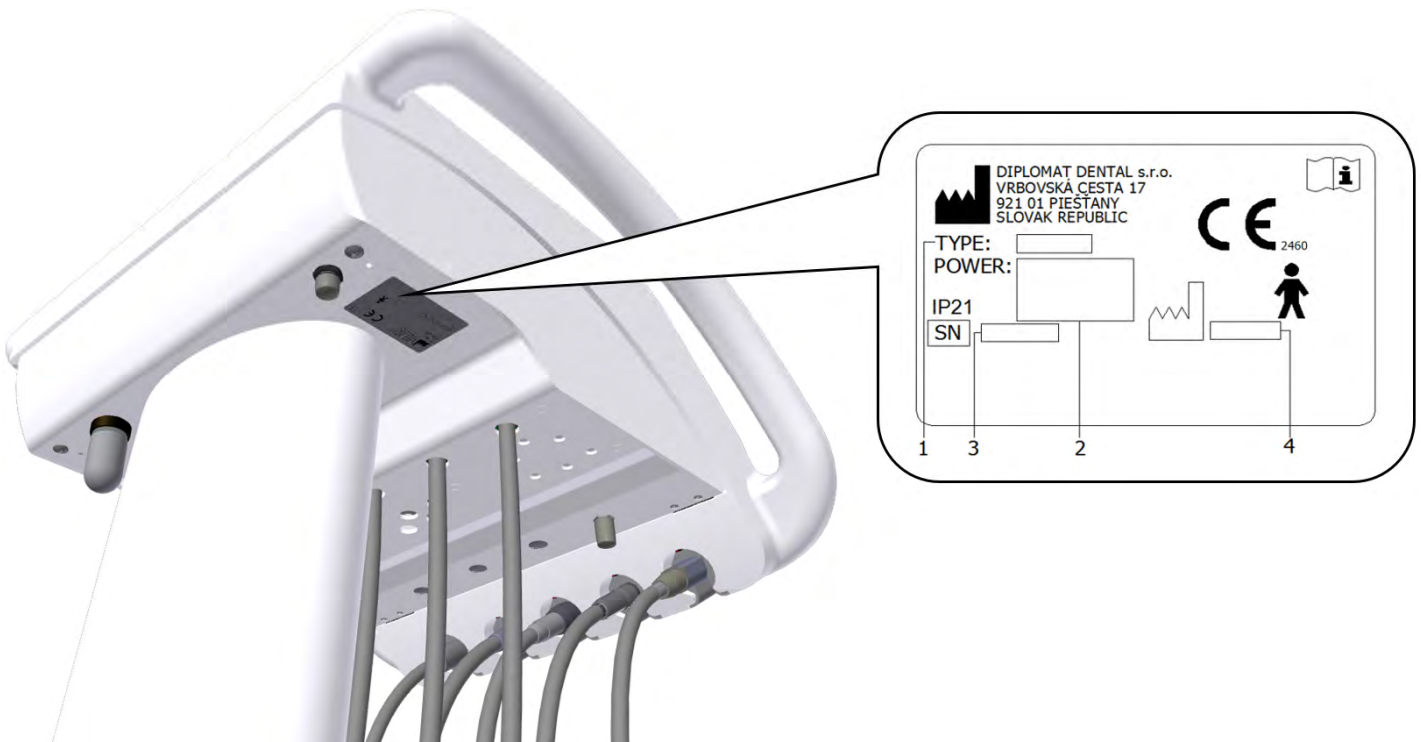4.1 Unit Label

5 PRE-INSTALLATION REQUIREMENTS

Caution

Pre-installation and installation must be performed according to the applicable standards of the particular country and in accordance with the valid documentation of the manufacturer, which is owned by each authorized representative of DIPLOMAT DENTAL s.r.o.

Caution

To avoid the risk of electric shock, device must be connected to the power supply with a protective earthing. Do not modify this device without prior authorization of the manufacturer.

5.1 Environmental conditions

Not to be installed in the premises with explosion hazard!

5.2 Requirements for the installation of media

Water

Distilled water is required to be used.

Air

Supply of at least **55 l/min.** of air at the pressure of **0,45 až 0,8 MPa**, oilless, clean and dry must be ensured. Tubes made of Cu and/or PE are recommended.

Recommended mains fuse rating

Recommended rating of the fuse of the supply main is 16A (in the event of circuit breaker – circuit breaker with switching-off characteristic of C type). No other items of equipment should be connected to the supply main in question!

Maximum power input of the dental unit is 150 VA. The supply main must comply with the respective national standard.

Recommendation

Unless the national standard stipulates otherwise, the manufacturer recommends to use current protective switch with the sensitivity of 30mA and instantaneous time of switching-off.

Pre-installation requirement having been met, assembly and installation of the dental unit is carried out and it is connected to the media.

5.3 Floor

The floor must have a concrete foundation of at least 100 mm thickness. The slope of the floor must not exceed 1%. The use of antistatic floor is recommended.

5.4 Environment

ambient temperature range	from +10° C	to +40° C
relative humidity range	from 30%	to 75 %
atmospheric pressure range	from 700hPa	to 1060hPa

6 ASSEMBLY AND INSTALLATION

The installation must be done by the certified service technician only; otherwise no possible future warranty claims will be accepted. The registration form must be filled out and sent to the manufacturer or the seller.

Unpacking of the unit and inspection of the delivery

Inspect the transport packages for damage. If a defect of the transport package is found, do not open the consignment and report the defect to the forwarding agent or seller immediately. If the consignment is intact, carefully open the package and unpack individual parts of the dental unit. Check the completeness of the delivery according to the

Packing list.

7 PUTTING THE UNIT INTO OPERATION

Warning – disinfection of new dental unit before its first use

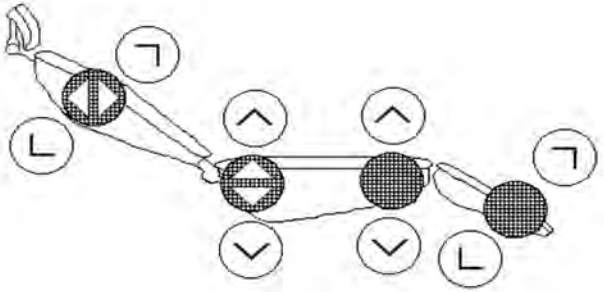
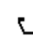
	<p>Before putting the unit into operation, your service technician must carry out disinfection of waterlines of instruments.</p>
--	--

- switch on the compressor and let it get pressurized with air
- check the instruments for their positions
- turn on the main switch

the indicator light of the main switch goes on and the following initial screen is displayed briefly:

After that the main screen is displayed. Grayed options are unavailable.

24.AUG	12:12					
						
1						
2						
3						
4						
DR.SMITH					TYGI 350 B	

The unit is connected to the distribution of air. **After audible signal (short-long tone), the unit is ready to work.**

While the unit is being in process of turning ON:

- no instruments should be taken off its holder,
- the foot controller shall be in the rest position and
- no button on the keyboard should be pressed

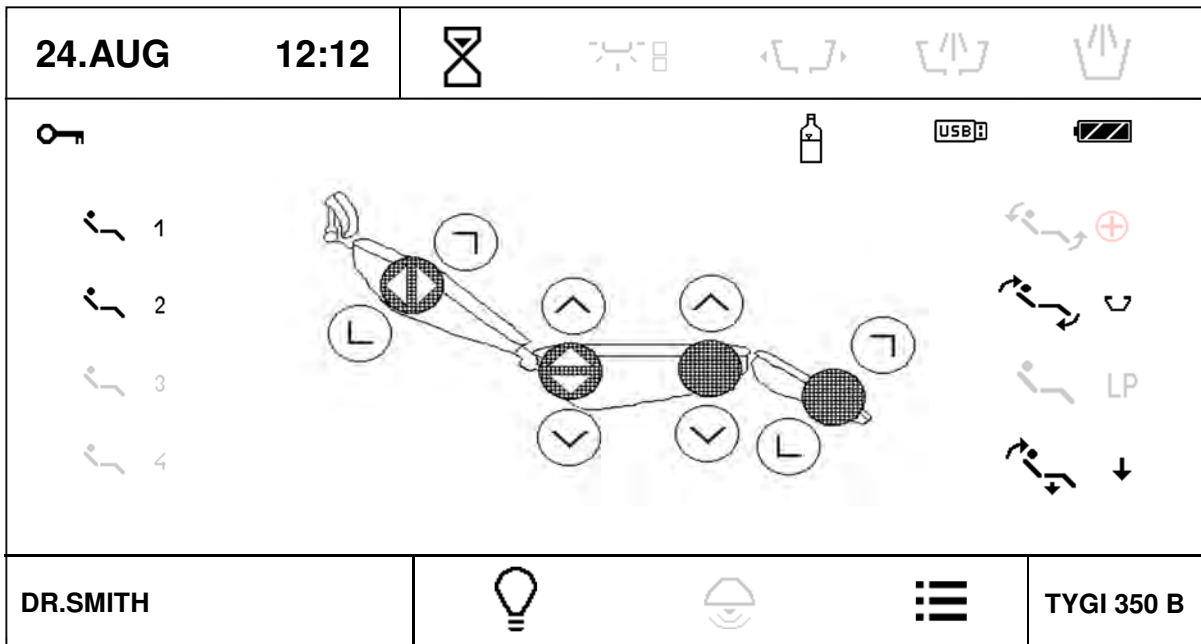
8 PRODUCT OPERATION

8.1 Control panel with instruments

The touch screen display is legible in each working position of the dentist (both sitting and standing).

8.1.1 Main screen

The main screen is displayed after turning the unit on.

Current time and date

24.AUG 12:12

Time and date can be set in user menu

Alarm

It serves to warn the dentist. It has the function of timer.

The time is adjustable. Press the button to display menu with available times. One of them marked as “U” can be set in the user menu. By pressing the button with desired time the time starts to count down to zero value by 1 second. By repeated pressing of the button during countdown the timer function is cancelled. When the time at countdown reaches 00:00, an audible signal sounds.

User

DR. SMITH

Displayed information about the current user. Selection of four users is available, i.e. data for four users is stored in the memory. By pressing the button a pop-up menu is displayed with available users.

DR.SMITH
DR.ANTOV
DR.POLAKOWSKI
DR.SLAVIK

To select the user, press the respective button. User name can be edited in the user menu. After the desired item has been pressed, the display is reset and the parameters are set for the selected user.

Turning the negatoscope on/off

Press the button to increase brightness of the display backlight to 100% and place the negative onto the display. To cancel the function, press the button again and the brightness will return to previous level. The negatoscope can be turned on when an instrument is taken off its holder.

If a picture is selected from inserted and recognized USB key by using USB menu, then the current picture is displayed instead of white background.

Locking of the keyboard

The keyboard can be „locked“ for cleaning. To block the keyboard, press and hold the button for approx. 3 seconds. The condition is indicated as follows:

To unlock the keyboard, press and hold the button for approx. 3 seconds. In both cases the waiting interval is indicated with a small bargraph.

Indication of the water source for instruments

The source is distilled water from the bottle.

Indication of USB key /USB flash drive/

After the USB key has been inserted into the USB connector, the icon  is displayed on the main screen.

If the icon is not displayed within 5 seconds after the key has been inserted, it means that the key has not been activated. If that is the case, pull the USB key out and insert it again quickly. If this fails 3 times, then the USB key probably cannot be recognized. If the icon  is highlighted, we can work with the USB key using the USB menu.


Indication of the accumulator condition in the wireless foot controller

The icon is displayed only when the wireless foot controller is connected. The shape of the icon indicates the condition of the accumulator in the foot controller:


	accumulator is discharged and needs to be connected to charging
	accumulator is charged to ca 33%
	accumulator is charged to ca 66%
	accumulator is charged to ca 100%
	accumulator is being charged

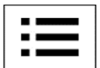Main menu


Press the button to get into the main menu of the unit.

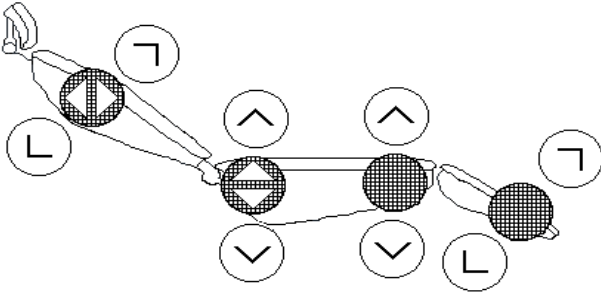Chair movement



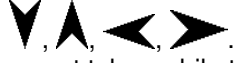
The basic movements of the chair are controlled by using of the buttons with the symbols of the direction of movement. The movement is performed while the button is kept pressed and is indicated by an additional indicator depending on

the orientation of the movement . All the said buttons control the chair directly with the instruments placed in position or with an instrument taken while the pedal of the foot controller is in zero position.

Programmable positions.

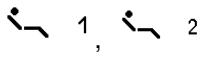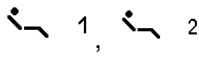These serve to save and recall pre-set positions of the chair.

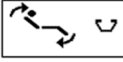The units enables two arbitrary user positions:  1,  2 and two special-purpose positions :

 /  Rinsing position / Return from rinsing position

 Getting-on position

Press the respective button to recall the desired position.

Programming positions

To program the position into the memory, set the chair to desired position and then press and hold the button for the said position for at least 2 seconds.

Getting into program mode is indicated by beeping and changing of the colour of the button frame. After the button has been released the position is written to the memory.

All programmable positions are saved separately for each user and are recalled from the memory automatically when the user has been selected.

8.1.2 System of changing parameters in individual screens

The unit enables the user to set individual parameters.

The parameters can be set using several ways.

1-Setting of the value of a parameter using + / -

24.AUG	12:12					
INSTRUMENT TYPE						40 000
> 1000	> 5000	> 20000	> 40000	X1/min		
	P1	PROGRAM 1				
	1:1		3.5 Ncm			
	NORM	Aft	3.0s			
ENDODONCY			33%			
DR.SMITH					TYGI 350 B	

Press the button with respective parameter. The buttons that enable to change the value are displayed in red. Change the value of the parameter by means of the buttons +/ - .

The new value is saved temporarily (until the change of the program number) by repeated pressing of the button with the parameter under setting or permanently into the program by means of the button

2-Setting of the value of a parameter using keyboard

Press the button to see the keyboard

MIN:100		MAX:40000		
!!!		0		
0	1	2	3	4
5	6	7	8	9
ESC	BS	.		OK


In the keyboard there is displayed the minimum and maximum permissible value for the parameter in question. If the value is outside the range, this is indicated by this symbol „!!!“ and the value will not be accepted when pressing „OK“. The button with decimal point is displayed only if the parameter in question enables so. To return without changing the value, press ESC.

The new value is saved temporarily (until the change of the program number) by repeated pressing of the button with the parameter under setting or permanently into the program by means of the button . Button „BS“ serves to delete the last position.

3-Setting of the value by selection from the options


Some parameters have pre-set values and their selection is made by means of a „pop-up menu“. When the button with the respective parameter is pressed, the menu with pre-set values is displayed, e.g. for the selection of transmission ratio:

	2:1	2,7:1
1:1	4:1	5,4:1
	7:1	7,4:1
1:2	10:1	14,8:1
1:4	16:1	
1:5	20:1	
	30:1	
	100:1	
	128:1	
ESC	U 1.0:1	

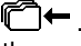To return without changing the value, press the button „ESC“. Press the respective value to set and save it temporarily (until the change of the program number). If you want the change to be saved permanently into the program, first press the button . Saving into the memory of the program is indicated by change in the colour . To cancel the saving into the memory, press the button again.

4-Editing the program name and the user name

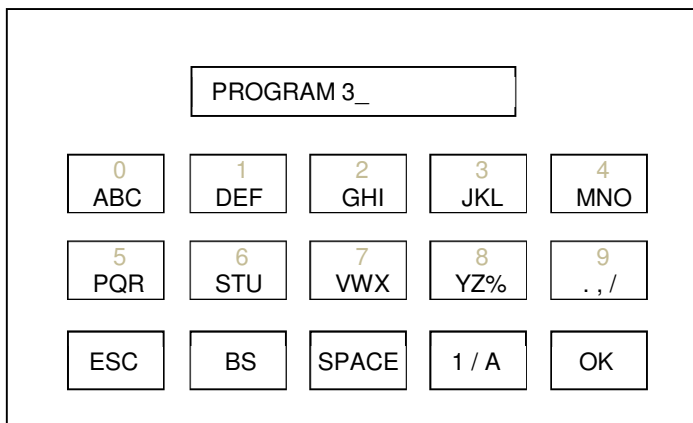Like numerical values, it is also possible to edit the text at program name and the user name. In that case the keyboard looks as follows:

Symbols explained:

ESC	The button serves to return without changing the text
BS	The button serves to delete the last sign
SPACE	The button serves to insert a gap
1/A	The button serves to change-over between numbers / letters
OK	The button serves to end editing and save the changes

Main menu

Press the button to get into the main menu of the unit.
It enables to get access to individual settings of the unit.


SETTINGS	
HYGIENE	ESC
USER	
SERVICE / MANUFACTURER	
TEST	
USB	
INFO	

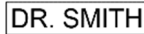USER menu

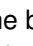
This option enables the user to set some parameters of the unit.

USER				
DR.SMITH		29°C	ESC <hr/>  <hr/>  <hr/>  <hr/> 	
	12:12			5.00:1
	29-APR-2017	 LCD		100%
	01:30			
	50%	 L		
	50%	 R		

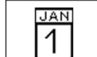
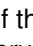
Editing the user menu

 DR. SMITH Press the button to edit the user menu and it is possible to edit the text. Only the text for the current user can be edited.

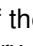
Setting the current time

 Press the button to start the mode of setting the value. Move to hours /minutes (current position displayed in red) by means of the button  and set the desired time by means of the buttons + / - . To save the said time into the memory, press the button again.

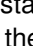
Setting the current date

 Press the button to start the mode of setting the value. Move to day /month / year (current position displayed in red) by means of the button  and set the desired date by means of the buttons + / - . To save the said date into the memory, press the button again.

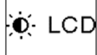
Setting the timer

 Press the button to start the mode of setting the value. Move to minutes / seconds (current position displayed in red) by means of the button  and set the desired time by means of the buttons + / - . To save the said time into the memory, press the button again. At selecting the time for the timer in the main screen this time is marked with „U“ symbol.


Setting the transmission ratio X:1

 Press the button to start the mode of setting the value. Move to the desired position (current position displayed in red) by means of the button  and set the desired value by means of the buttons + / - . To save the said value into the memory, press the button again. At selecting the transmission ratio in the „pop-up menu“ this value is marked with „U“ symbol.

Setting the intensity of display backlight


 LCD Press the button to start the mode of setting the value. Set the desired value of backlight intensity by means of the buttons + / - . To save the said value into the memory, press the button again.

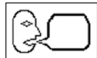Selection of language version

 Press the button to display the „pop-up menu“ with language options.

	CZ		HU
	D		HR
	F		SRB
	GB		FIN
	IT		P
	PL		TR
	PYC		ES
	SK		

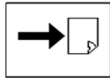Select the desired language by pressing the respective button. Press the button „ESC“ to return without change.

Moving to another page of the user menu.
Press the button to see other parts of the user menu.


USER				
				<p style="text-align: center; color: red;">ESC</p> <hr/> <p style="text-align: center;">+</p> <hr/> <p style="text-align: center;">-</p> <hr/> <p style="text-align: center;">→</p> <hr/> <p style="text-align: center;">←</p>

Turning the audible indication of button depression on/off

Press the button to turn the audible indication of button depression on / off. The condition is indicated by means of the symbol - turned ON / - turned OFF.

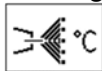Turning on/off the heating of the cooling water for the instruments

Press the button to turn the heating of the water for the instruments on / off. The condition is indicated by the symbol / . Applicable only if the heating of the water for the instrument is fitted.

Choosing of the colour palette

Press the button to display the „pop-up menu“ with the offer of colour palettes.


	Default palette 1
	Default palette 2
	Default palette 3
	User palette

Choose the desired palette by pressing the appropriate button.
Press the button „ESC“ if you do not want the change to take place.
The choose of the user palette is possible only after it has been set up (see below).

Setting the user palette

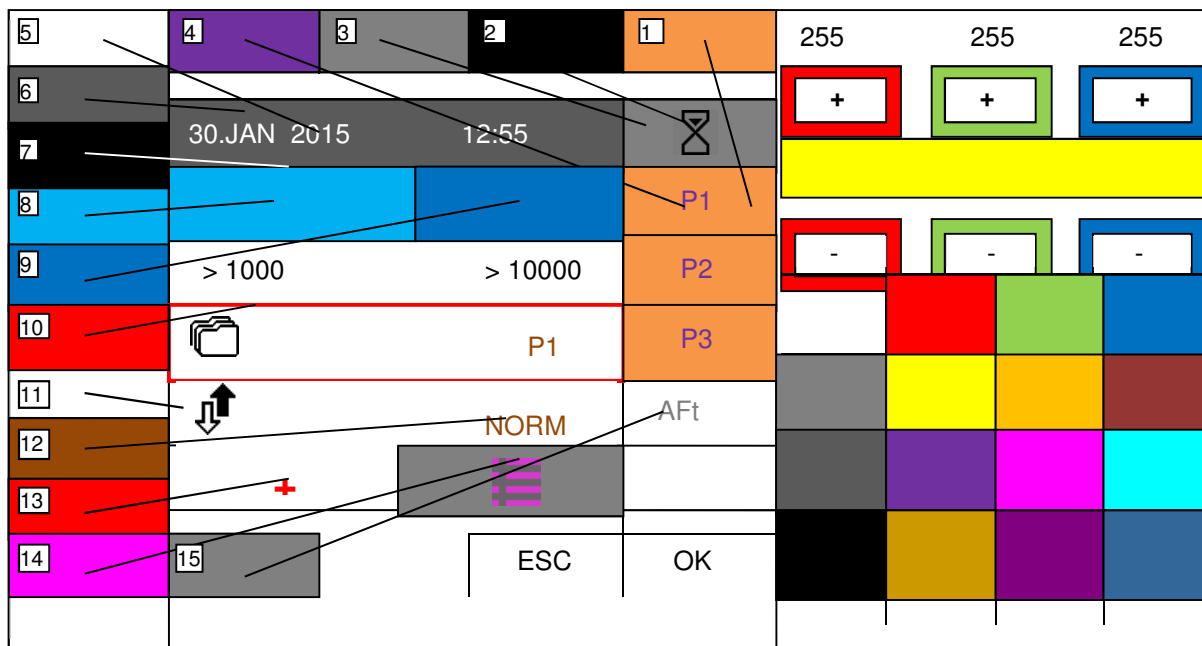Press the button to display the screen of setting of the colours:

Setting of the palette

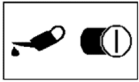
Select the desired colour from pre-set colours at the bottom right part of the screen. If you are not satisfied with the colour, alter it by means of the buttons at the top right part of the screen. The altered colour is displayed in the box between the buttons +/- . After the colour has been set, choose the item to which it should be assigned and press the respective button on the circumference of the left part of the screen to assign it. Repeat the procedure until all the desired colours have been set. Save the palette by pressing the button „OK“. When you press the button „ESC“ the data will not be saved. It is possible to set the following:

- 1- colour of the background of the pop-up menu
- 2- colour of icons
- 3- colour of the background of the active part of the upper and lower row of the display
- 4- colour of the text of the pop-up menu
- 5- colour of the text of the inactive part of the upper and lower row of the display– date, time, type
- 6- colour of the background of the inactive part of the upper and lower line of the display
- 7- colour of lines
- 8- colour of bargraph
- 9- colour of the background of the bargraph
- 10- colour of the frame at pressing the button
- 11- colour of the background of the screen
- 12- colour of adjustable parameters
- 13- colour of controls in active condition / +/-...
- 14- colour of the icon „main menu“
- 15- colour of hidden text and icon

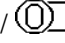
In the event of incorrect setting of the colours (the same colour for more items), some of the items might get illegible and the unit difficult to control. If that is the case, when at least two instruments are taken out of their holders, an error screen is displayed. At the top left corner there is the icon

. Press the icon to set the default colour palette 1.

Button for general enabling of alarm after the lubrication time of the instrument has been reached.

Press the button to enable / disable the indication of the lubrication time of the instrument being reached.

The condition is indicated by means of the symbol  - turned ON /  - turned OFF. The button enables/ disables the indication of the lubrication for all instruments simultaneously. Which instrument should be indicated is set individually for each instrument in the service menu.

Going to the previous page of the user menu.

Press the button to display the previous page of the user menu.

Cancelation of the user menu

Press the button „ESC“ to return to the main screen.

Menu SERVICE / MANUFACTURER

This option serves for service and production setting of the parameters. Particulars can be found in the Instructions for Service.

Menu TEST


This option serves to test some functions of the dental unit.

USB menu



This option serves for the work with USB key /flash drive/. USB key must be formatted to FAT32 exFAT. It is possible to read and display the file with the extension of **.bmp and uncompressed .jpg (max 800x480 pixels)** from the key /drive/ in the main folder. It is also possible to save error messages from SD card into the key /flash drive/ and to move the pictures of icons used from USB key into the SD card in the display. Individual items are accessible only after the USB key has been detected. The condition is indicated by means of the symbol **USB** in the main screen.



USB	
 ← USB	ESC
! → USB	
 ← USB	
 → USB	
 ← USB	

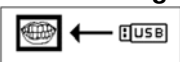Downloading of pictures from the USB key /flash drive/

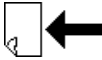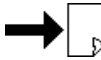
Press the button to download and display the folder with pictures:


	ESC	
00/ picture0.bmp		08/ picture8.jpg
01/ picture 1.bmp		09/ picture9.bmp
02/ picture2.bmp		10/ picture10.bmp
03/ picture3.jpg		11/ picture11.bmp
04/ picture4.bmp		12/ picture12.bmp
05/ picture5.bmp		13/ picture13.bmp
06/ picture6.jpg		14/ picture14.bmp
07/ picture7.bmp		15/ picture15.bmp

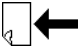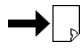Move between the pages in the folder by means of buttons  and . Select the desired picture by pressing the respective button – line. The picture is displayed after pressing the button. To return, press the button „ESC“.

Selected picture is displayed each time the button  is pressed.

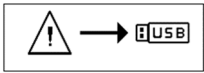Writing of errors into the USB key

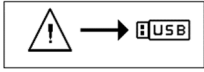Press the button to write the list of detected errors from the terminal memory into the USB key into the file „error.txt“. The said file can be browsed on PC in a text editor. The option is intended for service purposes.

Downloading of icons from the USB key into the SD card

Press the button to copy the files with icons from the USB key into internal SD card. In this way, the SD card do not need to be removed for writing of new icons before they are programmed into the internal FLASH memory. The option is intended for service purposes.

Saving of the user palette into the USB key

Press the button to save the colour palette created by the user into the USB key into the file „palette_usr.hex“. The palette of the current user is copied.

Copying of the user palette from the USB key

Press the button to copy the colour palette created by the user and saved in the USB key in the file „palette_usr.hex“. Copied palette is saved into the settings of the current user.

INFO menu

This menu serves for displaying information about HW and SW versions in each individual electronic module. After pressing button INFO, you will be shown such a screen:


INFO	
MODULE DP LUX:	SW 2.02
MODULE MB LUX:	SW 2.02
MODULE PB LUX:	SW 0.00
MODULE REG LUX:	SW 2.02
ESC	
SERIAL NUMBER:	12345678

8.1.3 Alarm indication

In the event an alarm occurs, the screen with alarm is displayed briefly and audible signal sounds.

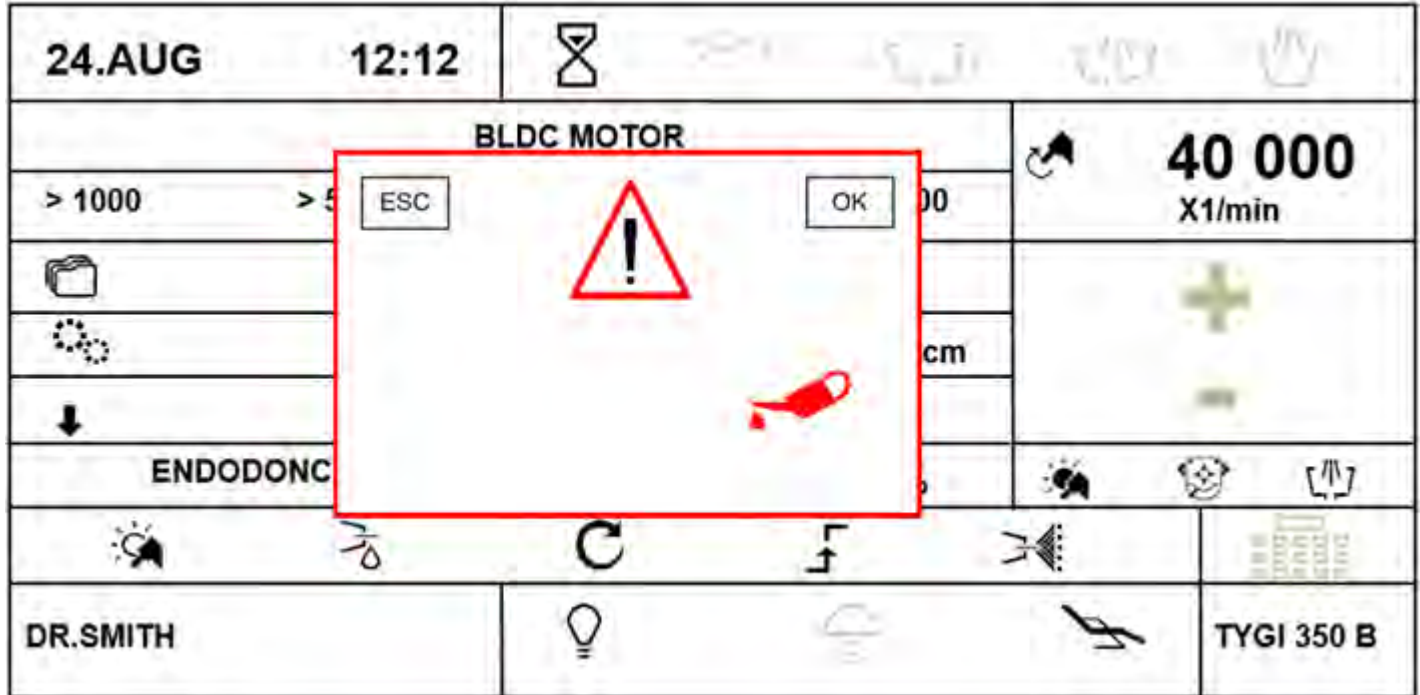Example:

After a short time the screen is cancelled automatically and the alarm is indicated by small icon in the bottom right corner of the display:

Press the button with alarm icon to display the original screen with the alarm type.

If you want to cancel the alarm of instrument lubrication, press the OK button in the screen, thus resetting the time counter for the lubrication alarm and a new interval will start to be counted.

i To be able to receive instrument lubrication period alarming, instrument lubrication has to be enabled by your service technician in service/manufacturer menu.

8.1.4 Operation of individual instruments of control panel

In addition to mentioned below, it is necessary to follow the instructions of the manufacturers of the instruments assembled with the dental unit.

8.1.5 Syringe

It is ready to operate already in the holder:

- For air blowing, press the right lever
- for water rinsing, press the left lever and
- to create water mist (spray), press both levers simultaneously

8.1.6 Turbine

24.AUG	12:12							
TURBINE							100 %	
> 10%		> 50%		> 100%				
	P1	PROGRAM 1				+		
			33%				-	
DR.SMITH						TYGI 350 B		

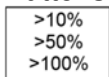
After the instrument has been taken the data on its setting is displayed. Its meaning is as follows:

Instrument type

If the bargraph is inactive, it displays the type of the instrument being used.

Bargraph in active condition displays graphically an informative value of the main parameter (power output).

Buttons of the fast pre-selection of the turbine output

Press the respective button to change the setting of the output to the said value. The data is not saved automatically into the program memory. The buttons are active only if the mode of output control is enabled.

Button for setting the turbine output

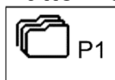
Press the button to start the setting mode in which the buttons whose values can be changed are displayed in red. Change the parameter value by means of the buttons +/- . The new value is saved temporarily (until the change of the program number) by pressing the button of the parameter being set again or permanently

into the program by means of the button

The new value can be entered also directly via numerical keyboard after pressing the button

The button is active only if the output control mode is enabled.

Button for the selection of a program.

Press the button to display a pop-up menu with program options.

PROGRAM 1	
PROGRAM 2	
PROGRAM 3	
PROGRAM 4	
PROGRAM 5	
PROGRAM 6	
PROGRAM 7	
PROGRAM 8	
PROGRAM 9	
ESC	

Press the respective button to select the desired program or cancel selecting a program by means of the button ESC
 To save all changed and unsaved parameters of instrument settings to program memory, just keep the button pressed for a longer time (approximately 2 seconds)
 The save is indicated by changed color of the icon (to red) and by sound signal (2 short beeps).

After the sound beeping, all data from temporal memory are written to permanent memory.

Button to assign program for switching from foot controller

By pressing on appropriate program, the program is included or excluded from programs switching via foot controller side button.

If the foot controller function P# has been chosen, then by repeated pressing of this button, all the included programs are being switched cyclically towards sequentially higher program number. After the included program achieves highest position number, then it again gets position with lowest sequential number.

Buttons for editing the program name.

Press and hold the button for ca 2 seconds. Alphanumeric keyboard is displayed by means of which the text can be edited. Maximum length of the text is 20 characters.

Button for setting the amount of cooling water (if proportional valve installed).

Press the button to start the setting mode in which the buttons whose values can be changed are displayed in red. Change the parameter value by means of the buttons +/- . The new value is saved temporarily (until the change of the program number) by pressing the button of the parameter being set again or permanently into the program by means of the button . The new value can be entered also directly via

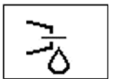
numerical keyboard after pressing the button . The button is active only if the instrument cooling is turned ON.

Button to turn the instrument light on / off.


Press the button to turn the instrument illumination on - / off - .

Button to turn the instrument afterblow on / off .

Press the button to turn the instrument afterblow on - / off - . The button is active only if the cooling is turned on.

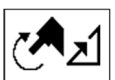Button to change the turbine operation mode.


Press the button to change the mode of the turbine operation. It has two conditions:

- mode of controlled output .
- mode of uncontrolled output , output always 100%.


The button is active only if the unit is equipped with the accessories for controlling the output for the turbine.

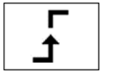Button to change the operation mode of the pedal / foot controller lever.

The button is active only if the output control mode is enabled. It indicates the discrete mode of the pedal, i.e. after pressing the foot controller pedal / lever, the controlled quantity jumps to maximum set values and does not change any more with the pedal position. By pressing the button the following is displayed:

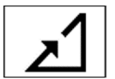Pedal in analogue mode


After pressing the foot controller pedal / lever, the controlled quantity increases proportionally to the deviation of the foot controller pedal / lever up to the maximum value at the full deviation of the foot controller pedal / lever. The button is active only if the instrument is in rest.

Button to turn the instrument cooling on / off.

Press the button briefly to turn the set type of the cooling of the instrument on / off .

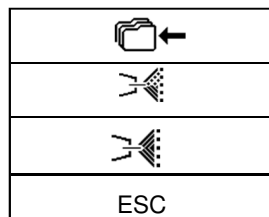By pressing and holding the button for ca 2 seconds the „pop-up “ is displayed with the following options:

Press the respective button to select the desired type of the cooling - water / , or

spray / , or cancel the instrument cooling type by pressing the button „ESC“. To save the changed type of the instrument cooling permanently into the memory, first press the button of saving into the memory . The new condition of saving the data into memory is indicated by change in the colour of the icon . This condition can be cancelled by pressing the button again.

Button to correct the patient position

Press the button briefly to display the screen to control the basic movements of the chair and then the patient position can be corrected. To return, press the button ESC or activate the instrument. It is still possible to correct patient position (while instrument is activated) with the basic moves of the chair via foot controller.

When a parameter is changed, while the instrument is taken out of its holder, without subsequent pressing the button

, the changed parameter will be saved only temporarily. The change will remain saved until the instrument program has been changed (Prog 1 – 9). The change will be maintained ever after the dental unit has been turned on /off. To cancel the temporary setting, press the button  and at the same time select the current program. Data will be over-written with the data from the memory for the said program. To save the change of the parameter permanently into the respective program, press the button  after the change.

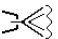
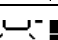
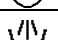
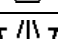
To activate the turbine, take it and move the foot controller lever to the right and/or the press the pedal of the combined foot controller. When moved to the extreme position, the turbine has the output displayed. To stop the activity, move the lever (pedal) of the foot controller to the starting position. Following the termination tool automatically goes blow through for 0.5 seconds. CHIPBLOWER function is recommended to use after the activity has been finished. The turbine output can be set by means of the buttons + / - within the range of 0-100% with the instrument taken in rest or in operation and the lever /pedal of the foot controller is in maximum position. When working with the instrument, the output can be adjusted continuously by pedal or by the lever of the foot controller in the analogue mode from 0 up to maximum set value.

After the operation time of the turbine has reached the lubrication interval, the instrument lubrication alarm  is displayed and audible signal sounds. This will be repeated each time the instrument is used until the information about lubrication need is cancelled by pressing the button OK in the alarm screen.

Assign function to the left / right side button of the foot controller at instrument control

Pressing of left/right side icon of the button opens „pop-up menu“ with the offer of functions. Select desired function by pressing the appropriate button. Press „ESC“ button if you want to cancel. You can setup both: LEFT/RIGHT side buttons of the foot controller. You can choose only those functions, which are available for the current tool. The others are disabled.

	Chipblower function
	Cooling of instrument ON/OFF
	Changing of the direction of motor rotation
	Instrument illumination ON/OFF
	Main light ON/OFF
	Bell
	Cup filling
	Bowl flushing
	No function chosen
P#	Switching the program of the instrument
	Switch color of the instrument light (white/blue) (applies to motors: DENSIM DX BLUE, DENSIM DX PRO BLUE)

Note

 After the turbine has been returned to the holder, the light goes off. Always return the micromotor or turbine to the holder only after the completion of activity (foot controller level (pedal) in the basic position). Used oil is collected in a receptacle – see the illustration below. If the oil exceeds the marked maximum value, the receptacle should be turned in the direction of the arrow and emptied. The filter should be replaced if contaminated.

**Retraction of procedural water**

 Upon the termination of use of dental instrument (after release of the foot controller pedal/lever), dental unit software automatically performs blowing of instrument head (micromotor, turbine) using cooling air of the instrument in a duration of 0.5 seconds.

8.1.7 Micromotor BLDC - DX, DX BLUE, DX PRO

24.AUG	12:12							
BLDC MOTOR							40 000 X1/min	
> 1000	> 5000	> 20000	> 40000					
	P1	PROGRAM 1						
	1:1		3.5 Ncm					
	NORM	Aft	3.0s					
ENDODONCY			33%					
DR.SMITH					TYGI 350 B			

After the instrument has been taken the data on its setting is displayed. Its meaning is as follows:

Instrument type

If the bargraph is inactive, it displays the type of the instrument being used. Bargraph in active condition displays graphically an informative value of the main parameter (power output).

>1000	Buttons of the fast pre-selection of the motor rotations. Press the respective button to change the setting of the output (rotations) to the said value. The data is not saved automatically into the program memory.
>5000	
>20000	
>40000	

Button for setting the motor rotations.

100%	Press the button to start the setting mode in which the buttons whose values can be changed are displayed in red. Change the parameter by means of the buttons +/- . The new value is saved temporarily (until the change of the program number) by pressing the button of the parameter being set again or permanently into the program by means of the button . The new value can be entered also directly via numerical keyboard after pressing the button
------	--

Button for the selection of a program.

P1	Press the button to display a pop-up menu with program options.
----	---

PROGRAM 1	<input checked="" type="checkbox"/>
PROGRAM 2	<input checked="" type="checkbox"/>
PROGRAM 3	<input checked="" type="checkbox"/>
PROGRAM 4	<input checked="" type="checkbox"/>
PROGRAM 5	<input checked="" type="checkbox"/>
PROGRAM 6	<input checked="" type="checkbox"/>
PROGRAM 7	<input checked="" type="checkbox"/>
PROGRAM 8	<input checked="" type="checkbox"/>
PROGRAM 9	<input checked="" type="checkbox"/>
ESC	

Press the respective button to select the desired program or cancel selecting a program by means of the button ESC To save all changed and unsaved parameters of instrument settings to program memory, just keep the button pressed for a longer time (approximately 2 seconds)
The save is indicated by changed color of the icon (to red) and by sound signal (2 short beeps).

After the sound beeping, all data from temporal memory are written to permanent memory.

Button to assign program for switching from foot controller

	By pressing on appropriate program, the program is included <input checked="" type="checkbox"/> or excluded <input type="checkbox"/> from programs switching via foot controller side button. If the foot controller function P# has been chosen, then by repeated pressing of this button, all the included programs are being switched cyclically towards sequentially higher program number. After the included program
--	---

achieves highest position number, then it again gets position with lowest sequential number.

Button for editing the program name.

Press and hold the button for ca 2 seconds. Alphanumeric keyboard is displayed by means of which the text can be edited. Maximum length of the text is 20 characters.

Button for the selection of transmission ratio of the handpiece.

Press the button to display a „pop-up menu“ with transmission ratio options:

	2:1	2,7:1
1:1	4:1	5,4:1
	7:1	7,4:1
1:2	10:1	14,8:1
1:4	16:1	
1:5	20:1	
	30:1	
	100:1	
	128:1	
ESC	U	1,0:1

Press the respective button to select the desired transmission ratio or cancel the selection of transmission ratio by pressing the button „ESC“. If you want to save the new value of the transmission ratio permanently into the memory, first press the button of saving into the memory . The new condition of saving the data into the memory is indicated by change in the colour of the icon . To cancel the condition, press the button again.

Button for setting the motor torque.

Press the button to start the setting mode in which the buttons whose values can be changed are displayed in red. Change the parameter value by means of the buttons +/- .

The new value is saved temporarily (until the change of the program number) by pressing the button of the parameter being set again or permanently into the program by means of the button . The button is active only in the mode ENDODONTIC.

Button for the selection of the motor operation mode NORM / AREV / AFWD

Mod normal – Normal

MX(BLDC) motor is controlled as a standard micromotor, except that the rotations can be set and it is possible to set the torque.

***Pos. Auto-reverse mode - AREV**

After reaching the set limit torque, the micromotor starts to rotate in the opposite direction (counter-clockwise rotations) until the foot controller has been released. Reverse running is indicated by flashing of the instrument light.

***Pos. Auto-forward mode - AFWD**

After reaching the set limit torque, the micromotor starts to rotate in the opposite direction (counter-clockwise rotations), namely until it has been released, then it remains in the said direction for the set time, e.g. AFt = 3,5 s and then return to the original rotations (clockwise rotations). This is repeated until the foot controller has been released. The reverse running is indicated by flashing tool and audible indication during the ATf. In this mode the reverse time can be set within the range of 0,1 – 5 seconds by 0,1 second.

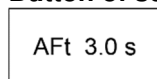
Press the button to display a „pop-up menu“ with the mode options:

↓ NORM
↕ AREV
↕ AFWD
ESC

Press the respective button to select the desired motor operation mode or cancel the selection of the mode by pressing the button „ESC“. If you want to save the new mode permanently into the memory, first press the button of saving into the memory .

The new condition of saving the data into the memory is indicated by change in the colour of the icon . To cancel the condition, press the button again. The button is active only if the instrument is at rest.

Button of setting the time for the mode AUTOFORWARD.

Press the button to start the setting mode in which the buttons whose values can be changed are displayed in red. Change the parameter value by means of the buttons +/- . The new value is saved temporarily (until the change of the program number) by pressing the button of the parameter being set again or

permanently into the program by means of the button .

Button for the working mode of the motor.

 Two working modes are available: ENDODONCY, PREPARATION. In accordance with the chosen mode, the individual items of the setting of the motor parameters are made accessible. Press the button repeatedly to alternate between the two modes.

Mode PREPARATION enables only the normal mode at maximum torque.

Mode ENDODONCY enables all modes – NORMAL/AREV/AFWD. The button is active only if the instrument is at rest.

Button for setting the amount of cooling water (if proportional valve installed).

 33% Press the button to start the setting mode in which the buttons whose values can be changed are displayed in red. Change the parameter value by means of the buttons +/- . The new value is saved temporarily (until the change of the program number) by pressing the button of the parameter being set again or permanently into the program by means of the button . The new value can be entered also directly via numerical keyboard after pressing the button . The button is active only if the instrument cooling is turned ON

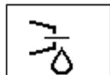
Button to turn the instrument light on / off.

 Press the button to turn the instrument illumination on -  / off - 

If the instrument allows switching two colors of the light (white / blue UV), then

- Quick pressing of the button (time < 1 second) switches color of the light  /  and
- Longer pressing of the button (time >= 1 second) turns the instrument light ON / OFF

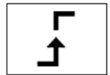
Button to turn the instrument afterblow on / off.

 Press the button to turn the instrument afterblow on -  / off - . The button is active only if the cooling is turned on.

Button to change the direction of the motor rotation.

 Press the button to change the direction of the motor rotation - . The button is active only at the rest position of the lever / pedal of the foot controller.

Button to change the operation mode of the pedal / lever of the foot controller.

 It indicates the discrete mode of the pedal, i.e. after pressing the foot controller pedal / lever, the controlled quantity jumps to maximum set values and does not change any more with the pedal position. By pressing the button the following is displayed:

Pedal in analogue mode

 After pressing the foot controller pedal / lever, the controlled quantity increases proportionally to the deviation of the foot controller pedal / lever up to the maximum value at the full deviation of the foot controller pedal / lever. The button is active only if the instrument is at rest.

Button to turn the instrument cooling on / off.

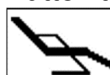
 Press the button briefly to turn the set type of the cooling of the instrument on / off. By pressing and holding the button for ca 2 seconds a „pop-up menu “ is displayed with the following options:

ESC

Press the respective button to select the desired type of the cooling - water / ; or

spray / , or cancel the instrument cooling type by pressing the button „ESC“. To save the changed type of the instrument cooling permanently into the memory, first press the button of saving into the memory . The new condition of saving the data into memory is indicated by change in the colour of the icon . This condition can be cancelled by pressing the button again.

Button to correct the patient position

 Press the button briefly to display the screen to control the basic movements of the chair and then the patient position can be corrected. To return, press the button ESC or activate the instrument. It is still possible to correct patient position (while instrument is activated) with the basic moves of the chair via foot controller.

When a parameter is changed, while the instrument is taken out of its holder, without subsequent pressing the button , the changed parameter will be saved only temporarily. The change will remain saved until the instrument program has been changed (Prog 1 – 9).

The change will be maintained ever after the dental unit has been turned on /off. To cancel the temporary setting, press the button  and select the current program.

Data will be over-written with the data from the memory for the said program. To save the change of the parameter permanently into the respective program, press the button  after the change.

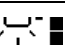
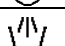
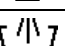
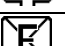
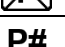
To activate the micromotor, take it and move the foot controller lever to the right and/or the press the pedal of the combined foot controller. When moved to the extreme position, the micromotor has the output displayed. To stop the activity, move the lever (pedal) of the foot controller to the starting position. Following the termination tool automatically goes blow through for 0.5 seconds. CHIPBLOWER function is recommended to use after the activity has been finished. The micromotor output can be set by means of the buttons + / - within the range of 0-100% with the instrument taken in rest or in operation and the lever /pedal of the foot controller is in maximum position. When working with the instrument, the output can be adjusted continuously by pedal or by the lever of the foot controller in the analogue mode from 0 up to maximum set value.

After the operation time of the micromotor has reached the lubrication interval, the instrument lubrication alarm  is displayed and audible signal sounds. This will be repeated each time the instrument is used until the information about lubrication need is cancelled by pressing the button OK in the alarm screen.

Assign function to the left / right side button of the foot controller at instrument control

Pressing of left/right side icon of the button opens „pop-up menu“ with the offer of functions. Select desired function by pressing the appropriate button. Press „ESC“ button if you want to cancel. You can setup both: LEFT/RIGHT side buttons of the foot controller. You can choose only those functions, which are available for the current tool. The others are disabled.

	Chipblower function
	Cooling of instrument ON/OFF
	Changing of the direction of motor rotation
	Instrument illumination ON/OFF
	Main light ON/OFF
	Bell
	Cup filling
	Bowl flushing
	No function chosen
P#	Switching the program of the instrument
	Switch color of the instrument light (white/blue) (applies to motors: DENSIM DX BLUE, DENSIM DX PRO BLUE)

 Availability of individual symbols on the display and options of setting of individual parameters and their limits depends on the type of the motor used. Before using the instrument, it is recommended to make yourself familiar with the performance of the instruments, according to the data from the instrument manufacturer.

Retraction of procedural water

 Upon the termination of use of dental instrument (after release of the foot controller pedal/lever), dental unit software automatically performs blowing of instrument head (micromotor, turbine) using cooling air of the instrument in a duration of 0.5 seconds.

8.1.8 Surgical micromotor DX SRG with peristaltic pump

24.AUG	12:12					
SRG MOTOR						40 000 X1/min
> 1000	> 5000	> 20000	> 40000			
	P1	PROGRAM 1				
	1:1		3.5 Ncm			
				40%		
DR.SMITH				TYGI 350 B		

After the instrument has been taken the data on its setting is displayed. Its meaning is as follows:

Instrument type

If the bargraph is inactive, it displays the type of the instrument being used.
Bargraph in active condition displays graphically an informative value of the main parameter (power output).

- >1000
- >5000
- >20000
- >40000

Buttons of the fast pre-selection of the motor rotations.

Press the respective button to change the setting of the output (rotations) to the said value. The data is not saved automatically into the program memory.

Button for setting the motor rotations.

Press the button to start the setting mode in which the buttons whose values can be changed are displayed in red. Change the parameter by means of the buttons +/- .

The new value is saved temporarily (until the change of the program number) by pressing the button of the parameter being set again or permanently into the program by means of the button . The new value can be entered also directly via numerical keyboard after pressing the button

Button for the selection of a program.


Press the button to display a pop-up menu with program options.

PROGRAM 1	<input checked="" type="checkbox"/>
PROGRAM 2	<input checked="" type="checkbox"/>
PROGRAM 3	<input checked="" type="checkbox"/>
PROGRAM 4	<input checked="" type="checkbox"/>
PROGRAM 5	<input checked="" type="checkbox"/>
PROGRAM 6	<input checked="" type="checkbox"/>
PROGRAM 7	<input checked="" type="checkbox"/>
PROGRAM 8	<input checked="" type="checkbox"/>
PROGRAM 9	<input checked="" type="checkbox"/>
ESC	


Press the respective button to select the desired program or cancel selecting a program by means of the button ESC
To save all changed and unsaved parameters of instrument settings to program memory, just keep the button pressed for a longer time (approximately 2 seconds)

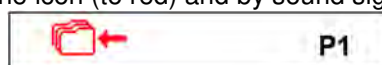The save is indicated by changed color of the icon (to red) and by sound signal (2 short beeps).

After the sound beeping, all data from temporal memory are written to permanent memory.

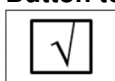Button to assign program for switching from foot controller

By pressing on appropriate program, the program is included or excluded from programs switching via foot controller side button.

If the foot controller function P# has been chosen, then by repeated pressing of this button, all the included programs are being switched cyclically towards sequentially higher program number. After the included program

achieves highest position number, then it again gets position with lowest sequential number.

Button for editing the program name.

Press and hold the button for ca 2 seconds. Alphanumeric keyboard is displayed by means of which the text can be edited. Maximum length of the text is 20 characters.

Button for the selection of transmission ratio of the handpiece.

Press the button to display a „pop-up menu“ with transmission ratio options:

	2:1	2,7:1
1:1	4:1	5,4:1
	7:1	7,4:1
1:2	10:1	14,8:1
1:4	16:1	
1:5	20:1	
	30:1	
	100:1	
	128:1	
ESC	U 1.0:1	

Press the respective button to select the desired transmission ratio or cancel the selection of transmission ratio by pressing the button „ESC“. If you want to save the new value of the transmission ratio permanently into the memory, first press the button of saving into the memory . The new condition of saving the data into the memory is indicated by change in the colour of the icon . To cancel the condition, press the button again.

Button for setting the motor torque.

Press the button to start the setting mode in which the buttons whose values can be changed are displayed in red. Change the parameter value by means of the buttons +/- .

The new value is saved temporarily (until the change of the program number) by pressing the button of the parameter being set again or permanently into the program by means of the button . The button is active only in the mode ENDODONTIC.

Button to active / deactivate the peristaltic pump

Press the button to activate - / deactivate - starting the peristaltic pump. If activated, the peristaltic pump is started together with the instrument when the lever / pedal of the foot controller has been moved.

Button for setting the flowrate and forced starting of the peristaltic pump.

The button is divided into two parts and it is active only if the peristaltic pump has been activated. To set the flowrate of the media through the pump, press the right part of the button - 40%. The setting mode is started, in which the buttons whose values can be changed are displayed in red. Change the parameters by means of the buttons +/- . The new value is saved temporarily (until the change of the program number) by pressing the button of the parameter being set again or permanently into the program by means of the button .

The new value can be entered also directly via numerical keyboard after pressing the button . Pres the button in the left part - to start the forced running of the peristaltic pump with the lever / pedal of the foot controller moved without the motor started. The condition remains until turned off by repeated pressing the button . It is possible to control forced run of peristaltic pump by side button of foot controller by assigning function Chipblower on this side button.

Button to turn the instrument light on / off.

Press the button to turn the instrument illumination on - / off - .

Button to change the direction of the motor rotation.

Press the button to change the direction of the motor rotation - . The button is active only at the rest position of the lever / pedal of the foot controller.

Button to change the operation mode of the pedal / lever of the foot controller.

The button is active only if the output control mode is enabled. It indicates the discrete mode of the pedal, i.e. after pressing the foot controller pedal / lever, the controlled quantity jumps to maximum set values and does not change any more with the pedal position. By pressing the button the following is displayed:

Pedal in analogue mode

After pressing the foot controller pedal / lever, the controlled quantity increases proportionally to the deviation of the foot controller pedal / lever up to the maximum value at the full deviation of the foot controller pedal / lever. The button is active only if the instrument is at rest.

Button to correct the patient position

Press the button briefly to display the screen to control the basic movements of the chair and then the patient position can be corrected. To return, press the button ESC or activate the instrument. It is still possible to correct patient position (while instrument is activated) with the basic moves of the chair via foot controller.

When a parameter is changed, while the instrument is taken out of its holder, without subsequent pressing the button

, the changed parameter will be saved only temporarily. The change will remain saved until the instrument program has been changed (Prog 1 – 9). The change will be maintained ever after the dental unit has been turned on

/off. To cancel the temporary setting, press the button  and select the current program. Data will be over-written with the data from the memory for the said program. To save the change of the parameter permanently into the

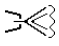
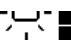
respective program, press the button  after the change.

To activate the micromotor, take it and move the foot controller lever to the right and/or the press the pedal of the combined foot controller. When moved to the extreme position, the micromotor has the output displayed. To stop the activity, move the lever (pedal) of the foot controller to the starting position. Following the termination tool automatically goes blow through for 0.5 seconds. CHIPBLOWER function is recommended to use after the activity has been finished. The micromotor output can be set by means of the buttons + / - within the range of 0-100% with the instrument taken in rest or in operation and the lever /pedal of the foot controller is in maximum position. When working with the instrument, the output can be adjusted continuously by pedal or by the lever of the foot controller in the analogue mode from 0 up to maximum set value.

After the operation time of the micromotor has reached the lubrication interval, the instrument lubrication alarm  is displayed and audible signal sounds. This will be repeated each time the instrument is used until the information about lubrication need is cancelled by pressing the button OK in the alarm screen.

Assign function to the left / right side button of the foot controller at instrument control

Pressing of left/right side icon of the button opens „pop-up menu“ with the offer of functions. Select desired function by pressing the appropriate button. Press „ESC“ button if you want to cancel. You can setup both: LEFT and RIGHT side buttons of the foot controller. You can choose only those functions, which are available for the current tool. The others are disabled.

	Chipblower function
	Cooling of instrument ON/OFF
	Changing of the direction of motor rotation
	Instrument illumination ON/OFF
	Main light ON/OFF
	Bell
	Cup filling
	Bowl flushing
	No function chosen
P#	Switching the program of the instrument
	Switch color of the instrument light (white/blue) (applies to motors: DENSIM DX BLUE, DENSIM DX PRO BLUE)

8.1.9 Scaler

24.AUG	12:12					
SCALER						100%
> 25%	> 50%	> 75%	> 100%			
	P1	PROGRAM 1			+	
ENDO		33%			-	
DR.SMITH				TYGI 350 B		

After the instrument has been taken the data on its setting is displayed. Its meaning is as follows:

Instrument type

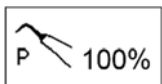
If the bargraph is inactive, it displays the type of the instrument being used. Bargraph in active condition displays graphically an informative value of the main parameter (power output).

Button of the fast pre-selection of the scaler output.

>25%
>50%
>75%
>100%

Press the respective button to change the setting of the output to the said value. The data is not saved automatically into the program memory.

Button for setting the scaler performance output.

Press the button to start the setting mode in which the buttons whose values can be changed are displayed in red. Change the parameter by means of the buttons +/- .

The new value is saved temporarily (until the change of the program number) by pressing the button of the parameter being set again or permanently into the program by means of the button . The new value can be entered also directly via numerical keyboard after pressing the button

Button for the selection of a program.


Press the button to display a pop-up menu with program options.

PROGRAM 1	<input checked="" type="checkbox"/>
PROGRAM 2	<input checked="" type="checkbox"/>
PROGRAM 3	<input checked="" type="checkbox"/>
PROGRAM 4	<input checked="" type="checkbox"/>
PROGRAM 5	<input checked="" type="checkbox"/>
PROGRAM 6	<input checked="" type="checkbox"/>
PROGRAM 7	<input checked="" type="checkbox"/>
PROGRAM 8	<input checked="" type="checkbox"/>
PROGRAM 9	<input checked="" type="checkbox"/>
ESC	

Press the respective button to select the desired program or cancel selecting a program by means of the button ESC To save all changed and unsaved parameters of instrument settings to program memory, just keep the button pressed for a longer time (approximately 2 seconds)

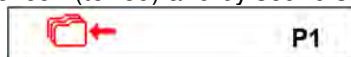The save is indicated by changed color of the icon (to red) and by sound signal (2 short beeps).

After the sound beeping, all data from temporal memory are written to permanent memory.

Button to assign program for switching from foot controller

By pressing on appropriate program, the program is included or excluded from programs switching via foot controller side button.

If the foot controller function P# has been chosen, then by repeated pressing of this button, all the included programs are being switched cyclically towards sequentially higher program number. After the included program achieves highest position number, then it again gets position with lowest sequential number.

Button for editing the program name.

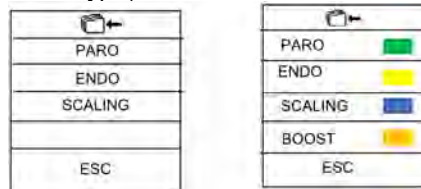
Press and hold the button for ca 2 seconds. Alphanumeric keyboard is displayed by means of which the text can be edited. Maximum length of the text is 20 characters.

Button for the selection of operation mode.

The button is active only for instrument that enables it. Press the button to display a „pop-up“ with the mode options (depending on instrument type):

Press the respective button to select the desired mode or cancel the selection of the mode by pressing the button „ESC“. If you want to save the new mode permanently into the memory, first press the button of saving into the memory . The new condition of saving the data into the memory is indicated by change in the colour of the icon . To cancel the condition, press the button again.

Button for setting the amount of cooling water.

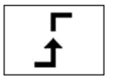
Press the button to start the setting mode in which the buttons whose values can be changed are displayed in red. Change the parameter value by means of the buttons +/- . The new value is saved temporarily (until the change of the program number) by pressing the button of the parameter being set again or permanently into the program by means of the button . The new value can be entered also directly via numerical keyboard after pressing the button . The button is active only if the instrument cooling is turned ON.

Button to turn the instrument light on / off.

Press the button to turn the instrument illumination on - / off - , if it is possible for the instrument in question.

Button to change the operation mode of the pedal / lever of the foot controller.

The button is active only if the output control mode is enabled. It indicates the discrete mode of the pedal, i.e. after pressing the foot controller pedal / lever, the controlled quantity jumps to maximum set values and does not change any more with the pedal position. By pressing the button the following is displayed:

Pedal in analogue mode

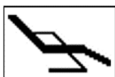
After pressing the foot controller pedal / lever, the controlled quantity increases proportionally to the deviation of the foot controller pedal / lever up to the maximum value at the full deviation of the foot controller pedal / lever. The button is active only if the instrument is at rest.

Button to turn the instrument cooling on / off.

Press the button briefly to turn the instrument cooling on/off. When a parameter is changed, while the instrument is taken out of its holder, without subsequent pressing the button , the changed parameter will be saved only temporarily. The change will remain saved until the instrument program has been changed (Prog 1 – 9). The change will be maintained ever after the dental unit has been turned on /off. To cancel the temporary setting, press the button and select the current program. Data will be over-written with the data from the memory for the said program. To save the change of the parameter permanently into the respective program, press the button after the change .

Button to correct the patient position

Press the button briefly to display the screen to control the basic movements of the chair and then the patient position can be corrected. To return, press the button ESC or activate the instrument. It is still possible to correct patient position (while instrument is activated) with the basic moves of the chair via foot controller.

When a parameter is changed, while the instrument is taken out of its holder, without subsequent pressing the button

, the changed parameter will be saved only temporarily. The change will remain saved until the instrument program has been changed (Prog 1 – 9).

The change will be maintained ever after the dental unit has been turned on /off. To cancel the temporary setting, press the button and select the current program.

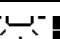
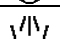
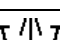
Data will be over-written with the data from the memory for the said program. To save the change of the parameter permanently into the respective program, press the button after the change.

To activate the scaler, take it and move the foot controller lever to the right and/or the press the pedal of the combined foot controller. When moved to the extreme position, the scaler has the output displayed. To stop the activity, move the lever (pedal) of the foot controller to the starting position. CHIPBLOWER function is recommended to use after the activity has been finished. The scaler output can be set by means of the buttons + / - within the range of 0-100% with the instrument taken in rest or in operation and the lever /pedal of the foot controller is in maximum position. When working with the instrument, the output can be adjusted continuously by pedal or by the lever of the foot controller in the analogue mode from 0 up to maximum set value.

Assign function to the left / right side button of the foot controller at instrument control

Pressing of left/right side icon of the button opens „pop-up menu“ with the offer of functions. Select desired function by pressing the appropriate button. Press „ESC“ button if you want to cancel. You can setup both: LEFT/RIGHT side buttons of the foot controller. You can choose only those functions, which are available for the current tool. The others are disabled.

	Chipblower function
	Cooling of instrument ON/OFF
	Changing of the direction of motor rotation
	Instrument illumination ON/OFF
	Main light ON/OFF
	Bell
	Cup filling
	Bowl flushing
	No function chosen
P#	Switching the program of the instrument
	Switch color of the instrument light (white/blue) (applies to motors: DENSIM DX BLUE, DENSIM DX PRO BLUE)


Note

After the scaler has been placed in the holder, the light goes off. Always return the scaler into the holder only after the activity has been finished (lever / pedal of the foot controller in the basic position).

8.2 Foot controller

UNO		
left top button – getting-on position		right top button –recalling of memory position of the chair
left bottom button – optional function		right bottom button - optional function
lever - start with continuous control - starting position (flushing)		cross switch -to control the chair

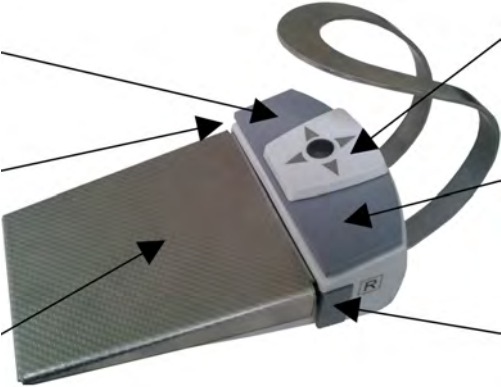NOK		
left top button – getting-on position		cross switch -to control the chair
left bottom button – optional function		right top button –recalling of memory position of the chair
pedal - start with continuous control - starting position (flushing)		right bottom button - optional function

With the micromotor and turbine instruments, when the button with optional function **CHIPBLOWER** is pressed on the foot controller cooling air will blow through the instrument tips.

Foot controller pedal/lever serves to put the instruments into operation, while with the micromotor, turbine it is possible to control the micromotor and turbine rotations by moving the lever of the foot controller (from 0 up to the value set on the display); and with the scaler, cauter, by moving the foot controller lever it is possible to control the output (from 0 up to the value set at the display).

The buttons **GETTING-ON POSITION**, **MEMORY POSITION** and **CROSS SWITCH TO CONTROL THE CHAIR** are intended to control the dental chair.

Memory positions are recalled by pressing the button **MEMORY POSITION + movement of the cross controller –**

upwards =  1 , **downwards** =  2

The second pressing / movement of the cross controller must be made within two seconds after pressing the button **MEMORY POSITION**, else the position will not be recalled.

Wireless foot controller (upon request only)

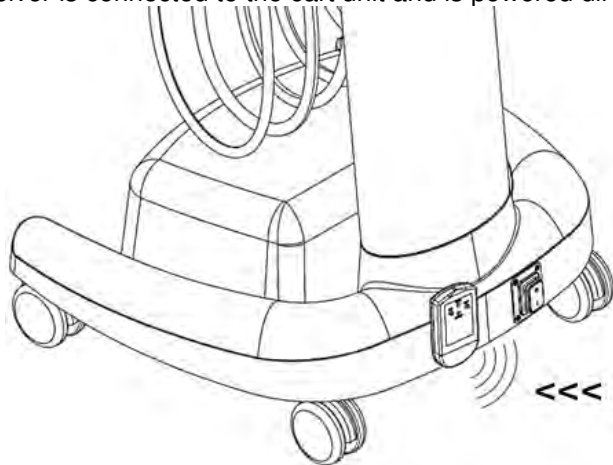
Warning

 WiFi UNO is supplied with a disconnected battery due to unwanted battery drain during periods of transport and storage. Putting into operation must be always done by authorized service technician.


The function of the buttons and the manner of control are identical with those of the cabled foot controller. The difference is in the way of the data transfer between the controller and the dental unit. While with the classic foot controller, the data transfer is made via a connected cable, with the wireless foot controller, the data transfer is carried out by means of radio waves without the need for a cable.

Standard connection of the wireless foot controller with the cart unit

Transmitter is located inside of the foot controller and is charged by battery. Receiver is connected to the cart unit and is powered directly.

<<< Receiver
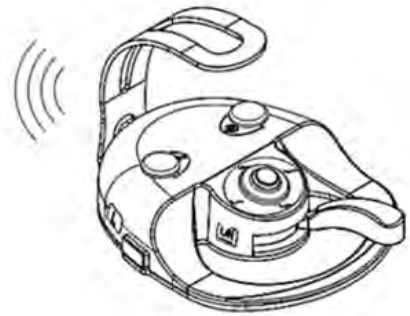
<<< Transmitter

Transmitter:

The foot controller is supplied from two accumulators NiMH 1,2V / 1800-2400 mAh.

Accumulators capacity is at least 2 month of standard operation. **The need for charging is indicated on the RECEIVER by regular flashing of the red LED diode or on the display by the symbol .**

After the moment of need for charging the accumulators has been indicated, it is still possible to work with the foot controller for a few more hours without connection to the charging.

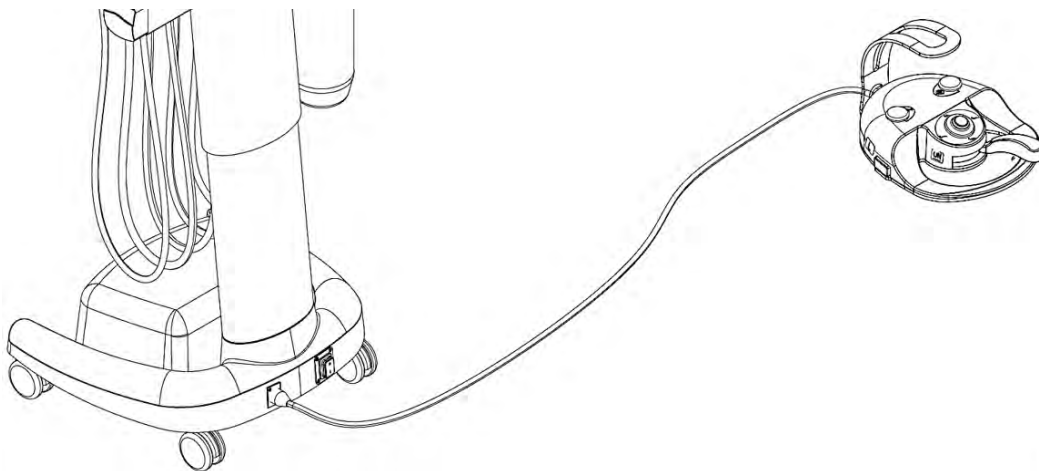
Charging of the foot controller:

The charging is activated by connecting the interconnecting cable between the unit and the foot controller.

During the charging it is still possible to use the WiFi foot controller – it functions as the classic cable foot controller.


The charging  takes about 3 hours and is stopped automatically once the accumulator has been charged.

Completion of the charging is indicated on the display by the symbol .

 To ensure the correct function of the foot controller and the unit, it is necessary to connect and disconnect the cable of the foot controller while the unit is turned off.

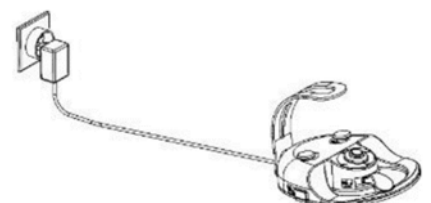 Optionally it is possible to charge the foot controller either via standalone outlet power charger. Ordercode of the power adapter is: ME374000.292.

Receiver

Button PAIR – serves for pairing transmitter with receiver

Button RESET – initialisation of device after settings change or malfunction of the foot controller.

LED1 (ON/RCV) – indicates connection to the power voltage and either detection of transmitting commands from the transmitter

LED2 (BATTERY) - indicates the charge status of the transmitter accumulators

The possible states of charge of the accumulators (LED2):

- Solid green – battery is fully charged
- Steady orange – battery is partially discharged
- Steady red – battery is nearly discharged
- Blinking red – it is necessary to recharge the battery

Procedure for charging the battery of the foot controller:

1. Pull out the receiver from the terminal.
2. Connect the battery charging cable of the foot controller to the terminal instead of the signal receiver
3. Connect the other end of the cable to the terminal on the back of the foot controller
4. Charging is indicated by the flash symbol on the battery indicator on the doctor's control panel display
5. The battery is fully charged when the flash symbol disappears from the battery indicator of the display
6. Disconnect the cable from both terminals
7. Connect the receiver to the terminal

If an error occurs on the transmitter during operation, the signal disappears and the receiver does not receive a new or confirming signal within 2 seconds, the receiver automatically sends the signal as a pedal in the zero position and interrupts the operation of the currently used instrument.

Pairing the receiver and the transmitter:

Each transmitter has a unique address set by the manufacturer. The transmitter comes from the production paired with a receiver. Should it happened that the transmitter is not paired with the receiver (when the instruction is sent from the foot controller, LED1 is not blinking) or in case of the foot controller or receiver replacement, it is necessary to pair the foot controller with the receiver first.

The setting must be performed by the service technician in accordance with the Instructions for Service.

List of possible error conditions

Failure	Possible cause	Solution
Receiver not working - LED 1 is not shining	The receiver is without power	Check if the dental unit is turned on Check if receiver is properly connected with the dental unit
	The receiver is defective	Connect the foot controller and the unit via cable and contact an authorized service or reseller
	The receiver works – LED1 is lit, but the receiver does not accept commands from the foot controller	Discharged battery of the foot controller
Charge the battery the foot of the driver		
The receiver, or the transmitter is defective		Contact an authorized service or reseller
Transmission barriers		Provide direct visibility between transmitter and receiver
At the maximum of pedal lever, the performance of the instruments does not reach the maximum	Loose / damaged potentiometer	Contact an authorized service or reseller
When the pedal is in the starting position the instruments continue to rotate	Loose switch for start	Contact an authorized service or reseller
The pedal does not return to its base position	Out of position or broken spring of pedal	Contact an authorized service or reseller
Button for control of the chair movement remains still pressed	Damaged foam spring button	Contact an authorized service or reseller

Technical data:

Transmitter	
Supply voltage – via cable – battery	10-36 V DC (rated voltage 24 V) 2,4 V DC 2xAA NiMH 1200 mAh (GP Recyko+ or sanyo eneloop)
Frequency	2400-2464 MHz
Battery life on a single charge	2 months
Transmitter reach	minimum 2 m
The transmitting power	-18 dBm according to SW settings till 4 dBm
Mudulation type	GFSK

Receiver	
Supply voltage	24 V DC
Frequency	2400-2464 MHz

8.3 Bottle with distilled water

The bottle with distilled water is located at the bottom of control panel. Distilled water from the bottle is used for the syringe, micromotor, turbine, scaler on the dentist panel.

Replenishing distilled water:

- turn the unit off
- unscrew the bottle
- replenish distilled water in the bottle
- screw the bottle firmly to avoid compressed air leakage during the work
- turn the unit on
- check the bottle for air leakage

In the event that distilled water has been spent and there is air in the water distribution, it is recommended to deaerate the water paths of instruments which use water by letting water flow until there are no bubbles in the water from instruments.

	<p>Warnings</p> <ul style="list-style-type: none"> • When replenishing distilled water, care must be taken, not to allow foreign substances penetrate into the water or to cause changes in its quality or composition. • Distilled water for medical purposes must be used with maximal electrical conductivity of the water up to 2000 µS/cm. • Do not use de-mineralized water for industrial purposes. • The manufacturer recommends to replace the bottle once a year.
---	--

8.4 Instruments cooling - setting of the water amount

To adjust water amount for instruments cooling, use solution which is applied on your unit:

Manual adjusting version

The amount of cooling water for all instruments, except for the syringe, can be adjusted by using of the needle valve located on the bottom side of the control panel. When the valve is turned fully clockwise, the water supply stops.

Electronic adjusting version

If the dental unit is fitted with a proportional valve, the amount of the cooling water can be adjusted from the touch panel for each instrument (except for the syringe). See illustration below.

12.09.2015	12:12					
INSTRUMENT TYPE						40 000 X1/min
> 1000	> 5000	> 20000	> 40000			
	P1	PROGRAM 1				
	1:1		3.5 Ncm			
	NORM					
ENDODONCY			33%			
DR.SMITH					UNIT MODEL	

Press this button to adjust amount of the cooling water

8.5 Adjusting the height of the unit

The working height of the control panel can be changed within the range of 200 mm. Grip the control panel (preferably while standing behind it) from the sides with both hands, while pressing (and holding) with a finger the button located on the bottom side of the panel, thus unblocking the mechanism of height adjustment. Re-position the panel gripped in the described way to the desired working height. Block the panel in the set height by releasing the mechanism button.

When moving the cart, ensure that the connection cable of the foot controller or any other object does not get in touch with the wheels of the cart.

8.6 Completion of work

The following is important:

- turn off the main switch – position 0
Thus disconnecting the supply of electric power, water and air and the whole unit is depressurized
- turn off /close/ the main supply of water (at the work place) to the dental unit
- turn off the compressor – open the sludge valve

9 PRODUCT MAINTENANCE

Maintenance of the instruments and handpieces should be carried out according to the instructions of their manufacturer. In case of the execution of the spittoon block with the connection to the central distribution, check the cleanliness of the advance filter (strainer) and the operation of the water hardness treatment device (according to the instructions of the manufacturer).

Inspections during the guarantee period

During the warranty period, the user is recommended regularly invite an authorized service technician for preventive periodical inspection in **3-month intervals**.

The inspection is focused on the following:

- check of input filters (with regard to the cleanliness of input media)
- providing additional information and practice advice concerning the dental unit
- check of the proper use and maintenance of the dental unit and its instruments (according to the Instructions for Use and training)
- check and/or adjustment of all media (input, setting of turbine pressures and the like)
- the manufacturer determines the range of above-mentioned work to be **1 to 1,5 hours**
- authorized service technician is obliged to confirm the periodical inspection carried out in the guarantee card

Inspection and audit /revision/ after the guarantee period:

The manufacturer recommends the periodical inspection of the product to be performed in **6-month intervals** by an authorized service technician, who is **to carry out** the following:

- general inspection of the dental unit and functional parts thereof,
- check and adjustment of the working pressures of water and air,
- check of air filter in the power block,
- check of the integrity of the electropart and electric wiring (electrical safety).

Audit /revision/ of electrical safety

It is carried out according to the regulations of the country in which the dental unit is installed.

10 CLEANING AND DISINFECTION

10.1 Elementary disinfection of the cooling waterlines of instruments

It is recommended to use Alpron or Sanosil S003 or Dentosept P agent in 1% concentration with distilled water. The solution is to be poured into the bottle for distilled water and can be used permanently. The agents in 1% concentration are harmless to the health of the patient. With regular use, the cooling system is maintained in cleanliness and it is not necessary to use other disinfecting agents. Contact your dealer for details about the possibilities of the purchase and use of the said products.

10.2 Cleaning and disinfection of other parts of the dental unit

Clean the external surfaces of the unit including the glass surface of the touch screen keyboard with a wet cloth. Use the INCIDIN FOAM – spray (HENKEL- ECOLAB) agent according to the instructions for use at least once a day or when the surface of the equipment is incidentally contaminated with biological material.

Caution

When cleaning the floor (PVC floor covering) with a disinfecting agent, it is forbidden to place the foot controller onto the still wet floor covering. It is forbidden to use agents with detrimental effects on the structure of the varnish or plastic materials (agents based on phenols and aldehydes).

The manufacturer shall not be taken responsible for damages arisen due to the use of disinfecting and cleaning agents others than those recommended.

Instruments and handpieces

Cleaning, disinfection and sterilization of the instruments and their handpieces should be carried out according to the instructions of their manufacturer, furnished with the instrument.

11 DISPOSAL OF THE EQUIPMENT

Part	Basic material	Recycleable material	Material to be placed in waste dumping grounds	Dangerous material
Metal	Steel	X		
	Aluminium	X		
Plastic	PUR		X	
	PVC			X
	PA, ABS	X		
	Laminated glass		X	
	Other plastic	X		
Rubber			X	
Glass		X		
Instruments			X	
Electronic		X		
Cable	Copper	X		
Transformer		X		
Amalgam separator	Filters			X
	Collecting vessel with amalgam			X
Packaging	Wood	X		
	Cardboard	X		
	Paper	X		
	PUR		X	

Note

In the disposal of the dental unit, the legislation of particular country should be observed. It is recommended to appoint an expert company with the disposal of the unit.

Not to be disposed of within the communal waste.
Waste material can be handed at destined places, e.g. electrical waste!

12 REPAIR SERVICE

In the event of a failure of the dental unit, contact the nearest service centre or your dealer who will give you information about the service network.

13 CONTENTS OF THE PACKAGING

Standard equipment

	TYGI 350 B
Mobile cart	1
Built-in Powerblock	1
Umbilical hose 5m	1
Foot controller	1
Tray table	1
Bottle for distilled water 2l	1
Instruments, accessories, small parts and completion sheet, sealed in a paper carton	1

Accompanying documentation

- These Instructions for Use
- Guarantee Card
- Manuals from subcontractors
- Completion sheet (placed in the sealed carton with the instruments)
- Registration form

14 WARRANTY

The manufacturer gives guarantee on the product according to the Gurantee card.

The danger of damage to the goods passes from the seller to the buyer by the moment of passing over the goods to the first forwarding agent for the transport for the buyer and/or by the moment of taking over the goods directly by the buyer.

The manufacturer reserves the right for changes within the innovation of the product.

When taking the product over into use the buyer is obliged to complete the guarantee form and to send it back to the manufacturer.

Caution

Defects /failures/ caused by negligent operation and/or by not observing the instructions in the Instructions for Use shall not be accepted as the subject of the claim under guarantee.

15 TRANSPORTING

Packages should be transported in covered vehicles in max. 3 plies and secured to prevent movement. Do not drop or tilt the package while loading or unloading.

16 STORAGE

Dental units can be stored in dry indoor warehouses in max. 3 plies with no dramatic temperature changes. Protect non-varnished parts from corrosion with rust-preventive treatment. Do not store together with chemicals!

Required environmental conditions:

relative humidity at max. 75%, temperature range for transporting and storage: -25°C to +50°C.

17 REQUIREMENTS ON ELECTROMAGNETIC COMPATIBILITY ACCORDING TO EN 60601-1-2

Caution

The use of accessories other than those stated in the Instructions for Use of the dental unit may result in increased electromagnetic emissions or decrease electromagnetic immunity and invoke disfunction of the dental unit.

Portable RF communications equipment must not be used at the distance of less than 30 cm from any part of the dental unit. Otherwise the functionality of the dental unit may be impaired.

17.1 Electromagnetic emissions

Emissions test	Compliance	Electromagnetic environment
RF emissions CISPR11	Group 1	The dental unit uses RF energy only for its function. Its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR11	Class B	The dental unit is designed to be used in all environments, including residential area and it can be directly connected to the public mains network.
Harmonic emissions EN 61000-3-2	Class A	
Voltage fluctuations/flicker emissions EN 61000-3-3	Complies	

17.2 Electromagnetic immunity

The dental unit is intended for use in the electromagnetic environment meeting the requirements in Table 17.2. The customer and/or the user of the dental unit must assure that the dental unit is used in such an environment. Table 17.2

Immunity test	EN 60601 test level	Compliance	Electromagnetic environment
Immunity to electrostatic discharge according to EN 61000-4-2	Contact discharge ± 6 kV Air discharge ± 8 kV	Contact discharge ± 6 kV Air discharge ± 8 kV	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Immunity to electrical fast transients and bursts according to EN 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	± 2 kV for power supply lines ± 1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Immunity to surges according to EN 61000-4-5	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth	Mains power quality should be that of a typical commercial or hospital environment.
Immunity to power frequency magnetic field according to EN 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should not exceed values typical for commercial or hospital environment.
Immunity to voltage dips, short interruptions and voltage variations on power supply input lines according to EN 61000-4-11	< 5% UT (> 95% dip in UT for 0,5 cycle) 40% UT (60% dip in UT for 5 cycles) 70% UT (30% dip in UT for 25 cycles) < 5% UT (> 95% dip in UT for 5 cycles)	< 5% UT (> 95% dip in UT for 0,5 cycle) 40% UT (60% dip in UT for 5 cycles) 70% UT (30% dip in UT for 25 cycles) < 5% UT (> 95% dip in UT for 5 cycles)	Mains power quality should be that of a typical commercial or hospital environment. If the user of the dental unit requires continued operation during power mains interruption, it is recommended that the dental unit be connected to a backup source.

17.3 Electromagnetic immunity

The dental unit is intended for use in the electromagnetic environment meeting the requirements in Table 17.3. The customer or the user must assure that the dental unit is used in such an environment.

Table 17.3

Immunity test	EN 60601 Test level	Compliance	Electromagnetic environment
Immunity to conducted disturbances, induced by radio-frequency fields according to EN 61000-4-6	3 V _{eff} 150 kHz to 80 MHz	3 V _{eff}	Portable and mobile RF communications equipment should be used no closer than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.
Immunity to radiated, radio-frequency electromagnetic field according to EN 61000-4-3	3 V/m 80 MHz to 2,5 GHz	3 V/m	<p>Recommended minimum distances:</p> $d = 1,167 \sqrt{P} \quad 150 \text{ kHz to } 80 \text{ MHz}$ $d = 1,167 \sqrt{P} \quad 80 \text{ MHz to } 800 \text{ MHz}$ $d = 2,333 \sqrt{P} \quad 800 \text{ MHz to } 2,5 \text{ GHz}$ <p>P[W]- rated maximum output power d[m] – recommended separation distance</p> <p>Field strengths from fixed RF transmitters should be less than the compliance level in each frequency range.</p> <p>Interference may occur in the vicinity of equipment marked with the symbol .</p>

Note 1

For frequency of 80 MHz the frequency range of 80 MHz to 800 MHz applies and for the frequency of 800 MHz the frequency range of 800 MHz to 2,5 GHz applies.

Note2

These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones, mobile radios, amateur radio, AM and FM radio and TV broadcast and the like cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to the fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the dental unit is used exceeds the applicable RF compliance level above, its normal operation must be verified. If abnormal performance is observed, the dental unit must be relocated.

Over the frequency range of 150 KHz to 80 MHz, field strengths should be less than 3 V/m.

17.4 Recommended separation distances between portable and mobile RF communications equipment and the dental unit

The dental unit is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the dental unit can help prevent electromagnetic interference by maintaining minimum distances between portable and mobile RF communications equipment and the dental unit according to Table 17.4.

Table 17.4

Rated maximum output power of transmitter P[W]	Separation distance according to frequency of transmitter d[m]		
	150 kHz to 80 MHz $d = 1,167 \sqrt{P}$	80 MHz to 800 MHz $d = 1,167 \sqrt{P}$	800 MHz to 2,5 GHz $d = 2,333 \sqrt{P}$
0,01	0,117	0,117	0,233
0,1	0,369	0,369	0,738
1	1,167	1,167	2,333
10	3,69	3,69	7,377
100	11,67	11,67	23,33

P[W]- rated maximum output power

d[m] – recommended separation distance

For transmitters rated at a maximum output power not listed above, the distance is calculated using the equation applicable to the respective frequency.

Note 1

For frequency of 80 MHz the frequency range of 80 MHz to 800 MHz applies and for the frequency of 800 MHz the frequency range of 800 MHz to 2,5 GHz applies.

Note2

These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

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