## VALO™ X Family Curing Lights

### Product Description

With its broadband spectrum, the VALO™ X family of curing lights are designed to polymerize all light-cured products in the wavelength range of 380–515 nm per ISO 10650.

The VALO X family of curing lights can be used in a corded or cordless configuration using the Ultradent VALO X rechargeable batteries or VALO X light cord adapter. These curing lights are designed to rest in a standard dental unit bracket or can be custom-mounted using the VALO surface mounting bracket.

**Note:** The VALO X family of curing lights include VALO X, VALO XG and VALO XS.

Overview of Controls:

|                                                                                                                                                                                                         |                                                                                                         |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------|
| Top button                                                                                                                                                                                              | Activates LED light (single button press)<br>Cycles between Curing Modes (1-second button hold)         |
| Bottom button                                                                                                                                                                                           | Activates LED light (single button press)<br>Cycles between Diagnostic Aid Modes (1 second button hold) |
| To cycle to Diagnostic Aid Modes from Curing Modes, press and hold the bottom button for 1 second. To cycle back to Curing Modes from Diagnostic Aid Modes, press and hold the top button for 1 second. |                                                                                                         |
| Alternatively, to cycle between Diagnostic Aid Modes and Curing Modes using the AF (accelerometer function), see section 8 (stepwise instructions).                                                     |                                                                                                         |

**Note:** Button location is the same for VALO X, VALO XS and VALO XG curing lights

Quick Start Guide for Indicators and Alerts (see stepwise instructions for modes relevant to each device):

| Mode Indicators on VALO X Family of Curing Lights | State or function    |
|---------------------------------------------------|----------------------|
| Blue and Green                                    | Curing Modes         |
| White and Violet                                  | Diagnostic Aid Modes |
| Flashing Green                                    | Sleep Mode           |

| Alert Indicators on VALO X Family of Curing Lights | State or function                         |
|----------------------------------------------------|-------------------------------------------|
| Red Flashing                                       | Low battery                               |
| Red Steady                                         | Replace/recharge battery                  |
| Amber                                              | Temporary LED over temperature protection |
| Teal Blue Chasing                                  | Mode actively running and light is busy   |

| Other equipment: Battery Charger | State           |
|----------------------------------|-----------------|
| Solid Amber                      | Charging        |
| Solid Green                      | Charge complete |
| Both off                         | Standby         |

For all products described, carefully read and understand all instructions and SDS information prior to use.

### Indications for Use/Intended Purpose

The VALO X family of curing lights are a source of illumination for curing photo-activated dental restorative materials and adhesives. It is also intended to provide illumination to aid in visualization during oral procedures. VALO X family of curing lights accessory/diffuser lenses are not intended for complete cure of photo-activated materials and adhesives.

### Contraindications

- No known contraindications have been identified for this product.
- For patients or users with allergy concerns, refer to product allergen document available at ultradent.com. If allergic reaction is observed, rinse exposed area thoroughly with water and have the patient consult their physician.

### Warnings and Precautions

| Risk Group 2                                                                                                                        |
|-------------------------------------------------------------------------------------------------------------------------------------|
| CAUTION Possibly hazardous optical radiation emitted from this product. DO NOT stare at operating lamp. May be harmful to the eyes. |
| Product tested against IEC62471                                                                                                     |

**Note:** This applies to VALO XG, and VALO XS only since they are risk group 2 under IEC62471. VALO X is risk group 1.

- CAUTION Possibly hazardous optical radiation is emitted from these products. DO NOT look directly into the curing lights' output. Patient, clinician, and assistants should always use amber-colored eye protection when curing lights are in use.
- To prevent the risk of electric shock, no modification of this equipment is allowed. Use only the included Ultradent VALO X light power supply and plug adapters. If these components are damaged, do not use and call Ultradent Customer Service to order a replacement.
- The VALO X family of curing lights are intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. Use of the VALO X family of curing lights adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, the VALO X family of curing lights and the other equipment should be observed to verify that they are operating normally. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the VALO X family of curing lights.
- Use only authorized accessories, cables, chargers, batteries, and power supplies to prevent improper operation, increased electromagnetic emissions, or decreased electromagnetic immunity
- Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 in.) to any part of the VALO X family of curing lights. Otherwise, degradation of the performance of this equipment could result.
- To avoid the risk of electrical fire associated with handling of batteries:
  - DO NOT autoclave or spray batteries, battery contacts, chargers, or AC power supply with liquid of any kind. If corrosion appears on the contacts of a charger, call Ultradent Customer Service to order a replacement.
  - DO NOT charge batteries around flammable materials
  - DO NOT keep charger in clinical operator or in an area where they will come into contact with disinfectant spray or water
  - DO NOT charge non-rechargeable batteries/cells. Charging non-rechargeable batteries can result in injury or property damage.
  - DO NOT recharge Ultradent batteries with a non-Ultradent charger
- To avoid the risk of injury, DO NOT use batteries that are corroded (rust), dented, emit an odor or fluids, have a torn or missing wrapping, or are otherwise damaged. Call Ultradent Customer Service to order replacement batteries/cells.
- Only use recommended batteries. Alternative batteries can cause malfunctions.
- To prevent the risk of thermal irritation or injury, avoid back-to-back curing cycles and do not expose oral soft tissues at close proximity for more than 10 seconds in Standard Power Mode, 5 seconds in Xtra Power Mode, or 3 seconds in Ortho Quadrant Mode. If longer curing times are required, use multiple curing cycles with periods of rest between cycles or use a dual-cure product to avoid heating soft tissue.
- Use caution when treating patients who suffer from adverse photobiological reactions or sensitivities, patients who are undergoing chemotherapy treatment, or patients being treated with photosensitizing medication
- These units may be susceptible to strong magnetic or static electric fields, which could disrupt the programming. If you suspect this has occurred, unplug the affected unit momentarily and then re-plug it into the outlet.
- DO NOT wipe down the curing lights with caustic or abrasive cleaners, autoclave, or immerse in any kind of ultrasonic bath, disinfectant, cleaning solution, or liquid. Failure to follow included processing instructions may render the curing lights inoperable.
- To avoid damaging the equipment, DO NOT insert fingers, instruments, or other objects into the battery/cell compartment of the curing light
- To avoid damaging the equipment, DO NOT attempt to clean the electrical contacts, or any part of the battery/cell compartment. Call Ultradent Customer Service if there is a concern.
- To help prevent cross-contamination and help keep dental composite material from adhering to the surface of the lens and wand body, an Ultradent approved barrier sleeve must be used over the VALO X family of curing lights with each use. Barrier sleeves are intended for single-patient use.
- To reduce the risk of corrosion, remove barrier sleeve after use.
- To reduce the risk of under-cured resins, do not use curing lights if lens is damaged
- To prevent the risk of product damage, do not autoclave or dry heat sterilize the VALO X family of curing lights accessory lenses

- DO NOT use any of the VALO X family of curing lights accessory lenses for complete/final cure. Remove all curing accessory lenses after use and perform a complete/final cure without the use of accessory lenses.
- DO NOT use illumination/visualization accessory lenses for curing dental restorative materials and adhesives

### Stepwise Instructions

#### Preparation

- Charge batteries before using the curing light.
- Prior to each use, place a new barrier sleeve over the curing light, minimizing wrinkles over the lens for best results.
- Select desired Curing Mode or Diagnostic Aid Mode. (See below for information on modes and how to cycle through each mode)
- For optimal results when curing, position curing light centered over and as close to resin as possible without contacting the surface.

**Note (VALO XS only):** Ortho mode is intended as a tacking mode with one pulse per bracket on a quadrant. As this is intended for tacking, an additional curing cycle would typically be performed to fully polymerize the composite.

Use

Curing Modes

VALO X:

| Curing Mode    | Power (mW) | Irradiance (mW/cm <sup>2</sup> ) | Total Exposure Time (Seconds) | Energy (Joules) | Button Indicator |
|----------------|------------|----------------------------------|-------------------------------|-----------------|------------------|
| Standard Power | 1,350      | 1,100                            | 10                            | 13.5            | Steady Blue      |
| Xtra Power *   | 2,700      | 2,200                            | 5                             | 13.5            | Pulsing Blue     |

VALO XG:

| Curing Mode    | Power (mW) | Irradiance (mW/cm <sup>2</sup> ) | Total Exposure Time (Seconds) | Energy (Joules) | Button Indicator |
|----------------|------------|----------------------------------|-------------------------------|-----------------|------------------|
| Standard Power | 1,077      | 1,000                            | 10                            | 10.8            | Steady Blue      |
| Xtra Power *   | 2,154      | 2,000                            | 5                             | 10.8            | Pulsing Blue     |

VALO XS:

| Curing Mode      | Power (mW) | Irradiance (mW/cm <sup>2</sup> ) | Total Exposure Time (Seconds) | Energy (Joules) | Button Indicator |
|------------------|------------|----------------------------------|-------------------------------|-----------------|------------------|
| Standard Power   | 747        | 1,000                            | 10                            | 7.5             | Steady Blue      |
| Xtra Power *     | 1,493      | 2,000                            | 5                             | 7.5             | Pulsing Blue     |
| Ortho Quadrant** | 1,493      | 2,000                            | 3 x 5                         | 4.48 x 5        | Pulsing Green    |

\* Nominal value.

\*Xtra Power has a triple pulse with a 0.5 second auto rest between pulses. Please note, there is a 2-second cooling period after using Xtra Power mode.

\*\*Ortho Quadrant consists of 5 consecutive bursts at 3 seconds each with a 2-second pause between pulses for repositioning. Please note, there is a 2-second cooling period after using Ortho Quadrant Mode.

**Changing between curing modes using "drum tap" motion or button push:**

**Drum Tap:** the VALO X family of curing lights are equipped with an accelerometer function allowing the user to cycle through Curing Modes via a "drum tap" motion.

Use a downward tap motion to activate/change Curing Mode (see visual aid below).

- Steady Blue Button Indicator = Standard Power
- Pulsing Blue Button Indicator = Xtra Power
- Pulsing Green Button Indicator = Ortho Quadrant (VALO XS)

**Button Push** – hold **top** button for 1 second to activate/change Curing Mode.

- Steady Blue Button Indicator = Standard Power
- Pulsing Blue Button Indicator = Xtra Power
- Pulsing Green Button Indicator = Ortho Quadrant (VALO XS)

Diagnostic Aid Modes

**VALO X**

| Diagnostic Aid Mode | Power <sup>2</sup> (mW) | Irradiance <sup>2</sup> (mW/cm <sup>2</sup> ) | Total Exposure Time (Seconds) | Button Indicator |
|---------------------|-------------------------|-----------------------------------------------|-------------------------------|------------------|
| Black Light         | 522                     | 425                                           | 60                            | Pulsing Violet   |
| White Light         | 282                     | 230                                           | 60                            | Pulsing White    |

**VALO XG**

| Diagnostic Aid Mode | Power <sup>2</sup> (mW) | Irradiance <sup>2</sup> (mW/cm <sup>2</sup> ) | Total Exposure Time (Seconds) | Button Indicator |
|---------------------|-------------------------|-----------------------------------------------|-------------------------------|------------------|
| Black Light         | 405                     | 376                                           | 60                            | Pulsing Violet   |

**VALO XS**

| Diagnostic Aid Mode | Power <sup>2</sup> (mW) | Irradiance <sup>2</sup> (mW/cm <sup>2</sup> ) | Total Exposure Time (Seconds) | Button Indicator |
|---------------------|-------------------------|-----------------------------------------------|-------------------------------|------------------|
| Black Light         | 287                     | 385                                           | 60                            | Pulsing Violet   |

<sup>2</sup> Nominal value, without accessory lens attached.

**Changing between Diagnostic Aid modes using "wave" motion or button push:**

**Wave:** the VALO X family of curing lights are equipped with an accelerometer function allowing the user to cycle through Diagnostic Aid modes via a "wave" motion.

Use a sideways tap motion to activate/change Diagnostic Aid Mode (see visual aid below)

- Violet Button Indicator = Black Light Diagnostic Aid Mode
- White Button Indicator (VALO X) = White Light Diagnostic Aid Mode


**Button Push:** Push and hold **bottom** button for 1 second to cycle through 2 Diagnostic Aid Modes.

Violet Button Indicator = Black Light Diagnostic Aid Mode – 1 cycle at 1 minute - Turn off with either button press or wave tap motion

White Button Indicator (VALO X) = White Light Diagnostic Aid Mode – 1 cycle at 1 minute - Turn off with either button press or wave tap motion

**Accessory/diffuser** lenses are attached magnetically over the curing light lens. Place barrier sleeve on curing light prior to attaching accessory lens.

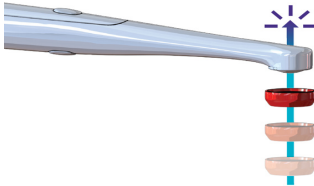**NOTE:** Upon holding either button, you will get an audible indication at 1 second.

**Sleep Mode**

Push and hold both buttons for 1 second to put in sleep mode. Upon holding both top and bottom buttons you will get an audible indication when entering sleep mode. When entering sleep mode, you must continue holding until the audible tone is complete. When you enter the sleep mode by pressing both buttons, the VALO X family of curing lights will not wake until either button is pressed.

**VALO X Family Accessories**

| Accessory                                | Mode                            | Description                                                                                                                      |
|------------------------------------------|---------------------------------|----------------------------------------------------------------------------------------------------------------------------------|
| PointCure™ Lens                          | Curing Modes                    | Augments the curing lights to polymerize composite through a translucent prosthetic.                                             |
| ProxiCure™ Ball Lens                     | Curing Mode                     | Augments the curing lights to polymerize composite and help shape contact area matrix of an interproximal restoration.           |
| Diffuser Lens (VALO X)                   | White Light Diagnostic Aid Mode | Augments the curing light to provide a visual aid for accurate color/shade comparison or whenever natural light is needed.       |
|                                          | Black Light Diagnostic Aid Mode | Augments the curing light to provide visualization of fluorescing chemicals in dental resins.                                    |
| Interproximal Lens (VALO X)              | White Light Diagnostic Aid Mode | Augments the curing lights in visualization of teeth and dental prostheses.                                                      |
| Interproximal Lens (VALO XG and VALO XS) | Standard Power Curing Mode      | Augments the curing lights in visualization of teeth and dental prostheses.                                                      |
| White light Accessory lens (VALO XG)     | Curing Mode                     | Augments the curing light to provide a visual aid for accurate color/shade comparison or whenever natural light is needed.       |
| Translume™ Lens                          | Curing or Diagnostic Aid Modes  | Augments the curing lights in visualization by providing longer wavelength light to transilluminate teeth and dental prostheses. |

**Preparation**

- Position accessory lens near the curing light lens allowing the magnet to snap the accessory lens into place over the barrier sleeve.

Use – PointCure Lens

- Curing Mode - The PointCure Lens concentrates the light into a 2.5 mm aperture. This is ideal for point curing (tacking) veneers and all porcelain crowns.
- Curing Mode - Point curing (tacking) translucent prosthesis: point cure prosthesis according to manufacturer guidelines of the luting material, material thickness, shade, opacity of material, and professional judgement. Clean up the uncured excess luting material around the margins, remove the accessory, then cure the entire restoration. Recommendations: Veneer - point cure the facial/buccal center of prosthesis. Crown - point cure facial/buccal and lingual center of prosthesis.

Use – ProxiCure Ball Lens

- Curing Modes - Supports the matrix and aids in convex proximal contact before and during partial light polymerization. Avoid entrapping ProxiCure Ball Lens in polymerized material. Remove accessory lens and cure the entire restoration for complete cure. Recommendation: Determine curing mode and time as per composite manufacturer's instructions.

Use – Diffuser Lens (VALO X)

- Black Light Diagnostic Aid Mode - Use to aid in visualization of fluorescent particles in various dental resins. Gives 60 seconds of illumination.
- White Light Diagnostic Aid Mode - Visual aid for whenever the dental professional may need a source of natural light, e.g., determining shade. Gives 60 seconds of illumination.

Use – Interproximal Lens (VALO X)

- White Light Diagnostic Aid Mode - Transillumination of teeth: use to aid in visualization of fractures, cracks, or defects in teeth by illuminating the teeth from the lingual and observing shadows.

Use – Interproximal Lens (VALO XG and XS)

- This accessory converts the blue light from the Curing Mode to white light for the transillumination of teeth: use to aid in visualization of fractures, cracks, or defects in teeth by illuminating the teeth from the lingual and observing shadow. Gives 10 seconds of illumination when using Standard Power Mode.

Use – White Light Accessory Lens (VALO XG)


- This accessory converts the blue light from the Curing Mode to white light for visual aid for whenever the dental professional may need a source of natural light, e.g., determining shade. Gives 10 seconds of illumination when using Standard Power Mode.

Use – Translume Lens

- White Light Diagnostic Aid Mode (VALO X) or Curing Mode - Neon-colored lens, which looks like the PointCure Lens, provides longer wavelength light to transilluminate the teeth to aid in visualization of cracks, fractures, defects, etc. The White Light Diagnostic Aid Mode gives 60 seconds of illumination. For curing, use the Standard Power Mode which gives 10 seconds of illumination.

Powering the VALO X Family of Curing Lights

- **Charging Batteries/Cells** - Remove cap at the base of the curing light with a counter-clockwise turn. Insert dead battery/cells in the provided Ultradent VALO X light battery charger with VALO X light power supply in the orientation shown in the charger. Battery/cells will be fully charged in 1–3 hours.
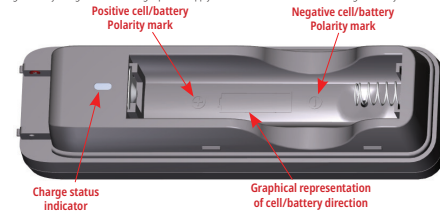
| Charge State Description | Charge Status Indicator | Voltage Condition    | Charge Current |
|--------------------------|-------------------------|----------------------|----------------|
| Pre-Charge               | Solid Amber             | 2.86 Vdc to 3.15 Vdc | 95 mA          |
| Bulk Charge              | Solid Amber             | 3.15 Vdc to 4.20 Vdc | 950 mA         |
| Charge Complete          | Solid Green             | >4.20Vdc             | N/A            |
| Standby                  | Both off                | N/A                  | N/A            |

- **Cord Adapter** - Remove cap at the base of the curing light with a counter-clockwise turn. Remove battery/cell. Insert cord adapter into battery/cell compartment, fully seat adapter and turn clockwise to lock into place. Position VALO X light power supply so it can be disconnected without difficulty from MAINS power.

**Mounting Bracket Instructions**

- 1) Bracket should be mounted to a flat, oil-free surface.
- 2) Clean surface with rubbing alcohol.
- 3) Peel backing off the bracket's adhesive tape.
- 4) Position bracket so the curing light lifts upward when removed. Press firmly into place.
- 5) If running the curing light on battery/cell power, the curing light should be oriented downward in the bracket. If running the curing light on the VALO X light cord adapter power, the curing light should be oriented upward in the bracket.

**Maintenance**

**General Cleaning for Curing Lights**

After each use, moisten a gauze or soft cloth with an approved surface disinfectant and wipe the surface and lens. Unauthorized cleaners may cause damage to the curing lights.

ACCEPTABLE CLEANERS:

- 70% Isopropyl alcohol
  - 70% Ethanol
- Use-Performed Maintenance
- 1) End of life is determined by wear and damage due to use. Inspect all components for damage prior to use and contact authorized service personnel to repair damage.
  - 2) Routinely check the lens for cured dental resins. If necessary, use a plastic or stainless-steel dental instrument to carefully remove any adhered resin.
  - 3) Routinely check the output in Standard Power mode. Light meters differ greatly and are designed for specific light guide tips and lenses. NOTE: the true numeric output will be skewed due to the inaccuracy of common light meters and the custom LED pack in the curing light; however, contact authorized service personnel if a decrease in output is observed when measuring with the same meter.

Manufacturer Repair

- 1) Repairs are only to be performed by authorized service personnel. Ultradent will provide service personnel with documentation to perform repairs.

**Warranty**

Ultradent Products, Inc. ("Ultradent") warrants that these products (including the original cord adapter that accompanied certain products when purchased, if applicable) shall, for a period of 5 years from the date of purchase, when operated according to the operating instructions included with the product, (i) conform in all material respects to the specifications set forth in Ultradent's documentation accompanying the product; and (ii) be free from defects in material and workmanship.

This limited warranty is nontransferable and applies solely to the original purchaser and does not extend to subsequent owners of the product. This limited warranty does not cover any other accessory components such as, but not limited to, batteries, chargers, or adaptive lenses. This limited warranty is void if the product fails or is damaged due to negligence, abuse, misuse, accident, modification, tampering, alteration, or failure to follow the applicable instructions for use. For purposes of example only, a product that is dropped and damaged is not covered under this warranty. To qualify under this limited warranty, proof of purchase (e.g., sales receipt or similar documentation) must be submitted to Ultradent along with the defective product.

A defective product meeting the warranty conditions set forth herein will, at Ultradent's sole discretion, either be repaired or replaced. In no event shall Ultradent's liability for the product exceed the purchase price paid by the purchaser. Under no circumstances shall Ultradent be liable for any indirect, incidental, foreseen, unforeseen, special, or consequential damages arising out of or in connection with the use of this product.

**Processing**

**Disinfection Reprocessing Instructions**

**VALO X Family of Curing Lights Processing Instructions**

| Disinfection Reprocessing Instructions |                                                                                                                                                                                                                                                                                                                                                                                                              |              |               |                        |                  |
|----------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------|---------------|------------------------|------------------|
| Limitations to Reprocessing            | A new barrier sleeve must be used for each patient. The use of a barrier sleeve is the primary means of contamination control for the device. In the event that a barrier sleeve is compromised or other extenuating circumstance, and the device is contaminated, follow the instructions included in this table to disinfect the device.<br><br>Note: The device is NOT intended for regular reprocessing. |              |               |                        |                  |
| Preparation Before Disinfection        | Use a 70% isopropyl alcohol (IPA) or 70% ethanol wipe to thoroughly wipe the entire surface. Discard soiled wipe                                                                                                                                                                                                                                                                                             |              |               |                        |                  |
| Rinsing                                | None                                                                                                                                                                                                                                                                                                                                                                                                         |              |               |                        |                  |
| Drying                                 | Allow to completely air dry at ambient temperature.                                                                                                                                                                                                                                                                                                                                                          |              |               |                        |                  |
| Maintenance, inspection, and testing   | Visually inspect that the device is not damaged and is free of debris. If debris is present, repeat preparation step with a new wipe.                                                                                                                                                                                                                                                                        |              |               |                        |                  |
| Packaging                              | No packaging required for disinfection.                                                                                                                                                                                                                                                                                                                                                                      |              |               |                        |                  |
| Disinfection                           | Wipe the entire surface and keep wet for the exposure time while paying special attention to cracks, crevices, seams, and hard-to-reach areas. Use additional wipes as necessary to keep the surface wet. <table><tr><td>Disinfectant</td><td>Exposure Time</td></tr><tr><td>70% IPA<br/>70% Ethanol</td><td>4 min.<br/>5 min.</td></tr></table>                                                             | Disinfectant | Exposure Time | 70% IPA<br>70% Ethanol | 4 min.<br>5 min. |
| Disinfectant                           | Exposure Time                                                                                                                                                                                                                                                                                                                                                                                                |              |               |                        |                  |
| 70% IPA<br>70% Ethanol                 | 4 min.<br>5 min.                                                                                                                                                                                                                                                                                                                                                                                             |              |               |                        |                  |
| Drying                                 | Allow to completely air dry at ambient temperature.                                                                                                                                                                                                                                                                                                                                                          |              |               |                        |                  |
| Storage                                | Store in a clean and dry location. To prevent moisture buildup, do not store in a barrier sleeve.                                                                                                                                                                                                                                                                                                            |              |               |                        |                  |

|                        |                                                                                                                                                                                                                                                                                                                                                                     |
|------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Additional information | This procedure was validated for intermediate-level disinfection by an independent and accredited laboratory. Note: Use of cleaning agents other than what is prescribed above, such as bleach or cleaners containing phosphoric acid, may cause damage to equipment. DO NOT attempt to clean the electrical contacts, or any part of the battery/cell compartment. |
|------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|

**VALO X Family of Curing Lights Accessories Processing Instructions**

| Disinfection Reprocessing Instructions                     |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                         |              |               |                                                            |           |
|------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------|---------------|------------------------------------------------------------|-----------|
| Initial Treatment at the Time of Use                       | Remove accessory prior to disinfection between each use.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                |              |               |                                                            |           |
| Preparation Before Disinfection                            | Remove any cured dental composites from accessory lenses prior to cleaning.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                             |              |               |                                                            |           |
| Cleaning                                                   | Place in an ultrasonic cleaner with 1 part detergent to 9 parts warm water using the product listed below: <table><tr><td>Detergent</td><td>Cleaning Time</td></tr><tr><td>Henry Schein General Purpose Cleaner or equivalent product</td><td>2-10 min.</td></tr></table>                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                               | Detergent    | Cleaning Time | Henry Schein General Purpose Cleaner or equivalent product | 2-10 min. |
| Detergent                                                  | Cleaning Time                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                           |              |               |                                                            |           |
| Henry Schein General Purpose Cleaner or equivalent product | 2-10 min.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                               |              |               |                                                            |           |
| Rinsing                                                    | Rinse accessory with warm water to remove detergent for 1–2 min..                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                       |              |               |                                                            |           |
| Drying                                                     | Dry with gauze. Air dry 30 min.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                         |              |               |                                                            |           |
| Maintenance, Inspection, and Testing                       | Visually inspect that the accessory is not damaged and free of debris. If cracked or damaged, discontinue use. If debris is present, repeat cleaning process.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                           |              |               |                                                            |           |
| Packaging                                                  | No packaging required for disinfection (sterilization packages can be used for storage after disinfection).                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                             |              |               |                                                            |           |
| Disinfection                                               | Remove from storage packaging (if applicable). Submerge accessory in disinfection solution as per manufacturer's instructions listed below: <table><tr><td>Disinfectant</td><td>Exposure Time</td></tr><tr><td>Cidex® OPA Solution or equivalent product</td><td>12 min.</td></tr></table>                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                              | Disinfectant | Exposure Time | Cidex® OPA Solution or equivalent product                  | 12 min.   |
| Disinfectant                                               | Exposure Time                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                           |              |               |                                                            |           |
| Cidex® OPA Solution or equivalent product                  | 12 min.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                 |              |               |                                                            |           |
| Rinsing                                                    | Rinse per disinfectant manufacturer's instructions listed below:<br>1. Following removal from CIDEX® OPA Solution, thoroughly rinse the accessory by immersing it completely in a large volume (e.g., 2 gal.) of water.<br>2. Keep the accessory totally immersed for a minimum of 1 min. in duration.<br>3. Remove the accessory and discard the rinse water. Always use fresh volumes of water for each rinse. Do not reuse the water for rinsing or any other purpose.<br>4. Repeat the procedure two (2) additional times, for a total of three (3) rinses, with large volumes of fresh water to remove CIDEX® OPA Solution residues. Residues may cause serious side effects. SEE WARNINGS. THREE (3) SEPARATE, LARGE VOLUME WATER IMMERSION RINSES ARE REQUIRED.<br>Additionally, rinse in water in ultrasonic cleaner for 5 min. |              |               |                                                            |           |
| Drying                                                     | Dry with gauze. Air dry 30 min.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                         |              |               |                                                            |           |
| Storage                                                    | Store in a clean and dry location.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      |              |               |                                                            |           |
| Additional information                                     | This procedure was validated by an independent and accredited laboratory.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                               |              |               |                                                            |           |

**Storage and Disposal**

If storing the curing light for periods longer than 2 weeks, or packing it for travel, always remove the batteries. If batteries are left in the unit for long periods of time without recharging, they may become nonfunctional or un-chargeable. Do not store batteries in temperatures over 40°C (104°F) or in direct sunlight.

Storage and Transport Conditions:

- Temperature: +10°C to +40°C (+50°F to +104°F)
- Relative Humidity: 10% to 95%
- Ambient Pressure: 500 hPa to 1060 hPa

Dispose of waste according to local rules, guidelines, and regulations.

When disposing of electronic waste, (i.e. devices chargers, batteries and power supplies), follow local waste and recycling guidelines.


**Technical Considerations**

Accessories

| Item                                                              |                                                                                     | CE Information                                                                                   |                                                                                                                    |
|-------------------------------------------------------------------|-------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------|
| VALO X Batteries<br>VALO X Power Supply<br>VALO X Battery Charger |  | Manufactured by:<br>GlobTek, Inc.<br>186 Veterans Drive,<br>Northvale, NJ 07647<br>Made in China | Distributed by:<br>Ultradent Products Inc.<br>505 West Ultradent Drive (10200 South)<br>South Jordan, UT 84095 USA |

**Technical Information/Data**


VALO X Curing Light:
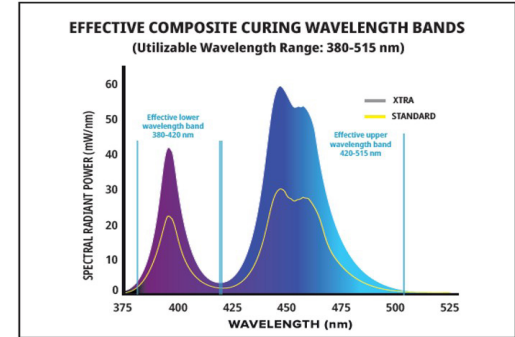
| Attribute             | Information/Specification                                                                                                                                                                                                                                                                                                                                                      |
|-----------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Lens                  | Diameter 12.5 mm = 1.2271 cm² (122.71 mm²)<br>PointCure Lens: Diameter 2.5 mm<br>ProxiCure Ball Lens: Diameter 2.1 mm                                                                                                                                                                                                                                                          |
| Wavelength range      | Utilizable wavelength range: 380–515 nm<br>Peak wavelengths: 380–420 nm and 420–515 nm                                                                                                                                                                                                                                                                                         |
| Light Intensity Table | Standard Power - 1,100 mW/cm² (±10%)<br>Xtra Power - 2,200 mW/cm² (±10%)<br>(Conforms to ISO 10650 when measured with a Gigahertz spectrum analyzer)<br>Black Light - 200-648 mW/cm²<br>White Light - 160-300 mW/cm²<br>NOTE: PointCure and ProxiCure Ball lenses focus the light energy into a smaller aperture, thus increasing the light intensity above calibrated values. |

|                            |                                                                                                                                                                                                                                                                                                                   |                                                                                                                                                                                                                       |
|----------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| VALO X Curing Light        | The curing light with a barrier sleeve is an applied part.<br>Ratings: IEC 60601-1 (Safety), IEC 60601-1-2 (EMC), Ingress Protection Rating: IP54                                                                                                                                                                 | Weight:<br>With battery: 4.8 oz. (136 grams)<br>Without battery: 3.8 oz. (108 grams)<br>With cord adapter: 5.6 oz. (158 grams)<br><br>Dimension (battery configuration): 8.9 x 0.83 x 0.83 inches, (226 x 21 x 21 mm) |
| Operating Conditions       | Temperature: +10°C to +32°C (+50°F to +90°F)<br>Relative Humidity: 10% to 95%<br>Ambient Pressure: 700 hPa to 1060 hPa                                                                                                                                                                                            |                                                                                                                                                                                                                       |
| Curing Mode<br>Duty Cycle: | At maximum ambient temperature (32°C) Xtra Power Mode has 5 seconds ON and 60 seconds OFF.<br>Note - Duty cycle is defined by primary function of device (polymerize resin) using Xtra Power mode, which creates highest amount of thermal energy.<br>Maximum Device temperature when following duty cycle: ≤45°C |                                                                                                                                                                                                                       |


VALO XG Curing Light:
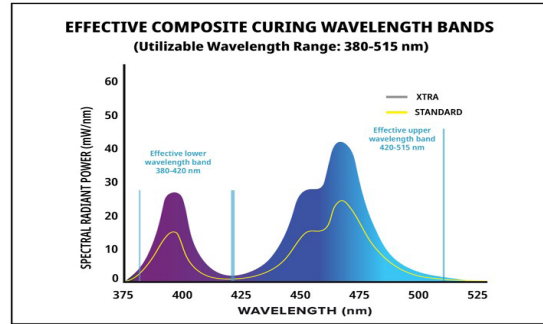
| Attribute             | Information/Specification                                                                                                                                                                                                                                                                                                                                                          |
|-----------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Lens                  | Diameter 11.71 mm = 1.0770 cm <sup>2</sup> (107.70 mm <sup>2</sup> )<br>PointCure Lens: Diameter 2.5 mm<br>ProxiCure Ball Lens: Diameter 2.1 mm                                                                                                                                                                                                                                    |
| Wavelength range      | Utilizable wavelength range: 380-515 nm<br>Peak wavelengths: 380-420 nm and 420-515 nm                                                                                                                                                                                                                                                                                             |
| Light Intensity Table | Standard Power - 1,000 mW/cm <sup>2</sup> (±10%)<br>Xtra Power - 2,000 mW/cm <sup>2</sup> (±10%)<br>(Conforms to ISO 10650 when measured with a Gigahertz spectrum analyzer)<br>Black Light - 180-594 mW/cm <sup>2</sup><br>NOTE: PointCure and ProxiCure Ball lenses focus the light energy into a smaller aperture, thus increasing the light intensity above calibrated values. |

|                            |                                                                                                                                                                                                                                                                                                                   |                                                                                                                                                                                                           |
|----------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| VALO XG Curing Light       | The curing light with a barrier sleeve is an applied part.<br>Ratings: IEC 60601-1 (Safety), IEC 60601-1-2 (EMC), Ingress Protection Rating: IP54                                                                                                                                                                 | Weight:<br>With battery: 4.8 oz. (134 g)<br>Without battery: 3.8 oz. (105 g)<br>With cord adapter: 5.5 oz. (156 g)<br><br>Dimension (battery configuration): 8.9 x 0.83 x 0.83 inches, (226 x 21 x 21 mm) |
| Operating Conditions       | Temperature: +10°C to +32°C (+50°F to +90°F)<br>Relative Humidity: 10% to 95%<br>Ambient Pressure: 700 hPa to 1060 hPa                                                                                                                                                                                            |                                                                                                                                                                                                           |
| Curing Mode<br>Duty Cycle: | At maximum ambient temperature (32°C) Xtra Power Mode has 5 seconds ON and 60 seconds OFF.<br>Note - Duty cycle is defined by primary function of device (polymerize resin) using Xtra Power mode, which creates highest amount of thermal energy.<br>Maximum Device temperature when following duty cycle: ≤45°C |                                                                                                                                                                                                           |

VALO XS Curing Light:
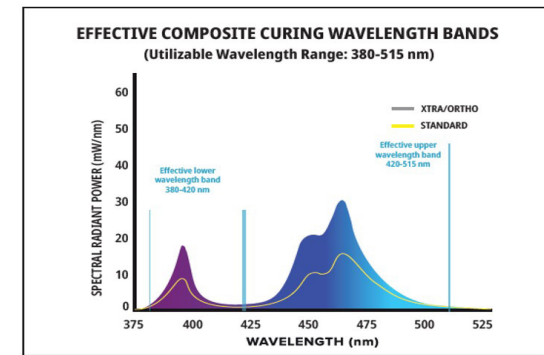
| Attribute        | Information/Specification                                                                                                                       |
|------------------|-------------------------------------------------------------------------------------------------------------------------------------------------|
| Lens             | Diameter 9.75 mm = 0.74662 cm <sup>2</sup> (74.662 mm <sup>2</sup> )<br>PointCure Lens: Diameter 2.5 mm<br>ProxiCure Ball Lens: Diameter 2.1 mm |
| Wavelength range | Utilizable wavelength range: 380-515 nm<br>Peak wavelengths: 380-420 nm and 420-515 nm                                                          |

|                            |                                                                                                                                                                                                                                                                                                                                                                                                                                          |                                                                                                                                                                                                           |
|----------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Light Intensity Table      | Standard Power - 1,000 mW/cm <sup>2</sup> (±10%)<br>Xtra Power - 2,000 mW/cm <sup>2</sup> (±10%)<br>Ortho Mode Power - 2,000 mW/cm <sup>2</sup> (±10%)<br>(Conforms to ISO 10650 when measured with a Gigahertz spectrum analyzer)<br>Black Light - 180-594 mW/cm <sup>2</sup><br>NOTE: PointCure and ProxiCure Ball lenses focus the light energy into a smaller aperture, thus increasing the light intensity above calibrated values. |                                                                                                                                                                                                           |
| VALO XS Curing Light       | The curing light with a barrier sleeve is an applied part.<br>Ratings: IEC 60601-1 (Safety), IEC 60601-1-2 (EMC), Ingress Protection Rating: IP54                                                                                                                                                                                                                                                                                        | Weight:<br>With battery: 4.8 oz. (133 g)<br>Without battery: 3.7 oz. (105 g)<br>With cord adapter: 5.6 oz. (156 g)<br><br>Dimension (battery configuration): 8.8 x 0.83 x 0.83 inches, (224 x 21 x 21 mm) |
| Operating Conditions       | Temperature: +10°C to +32°C (+50°F to +90°F)<br>Relative Humidity: 10% to 95%<br>Ambient Pressure: 700 hPa to 1060 hPa                                                                                                                                                                                                                                                                                                                   |                                                                                                                                                                                                           |
| Curing Mode<br>Duty Cycle: | At maximum ambient temperature (32°C) 5 seconds ON and 60 seconds OFF.<br>Ortho Mode has 15 seconds ON and 60 seconds OFF.<br>Note - Duty cycle is defined by primary function of device (polymerize resin) using Xtra Power Mode and Ortho Mode, which creates highest amount of thermal energy.<br>Maximum Device temperature when following duty cycle: ≤45°C                                                                         |                                                                                                                                                                                                           |

VALO X Component Information (applies to all curing lights in VALO X family):

| Attribute                     | Information/Specification                                                                                                                                                                                                                                                                                                                                                                              |                                                                                                                                                                      |
|-------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| VALO X Power Supply           | Output - 9VDC at 2A<br>Input - 100VAC to 240VAC, 50-60 Hz<br>Ultradent P/N 4952 - VALO X Power Supply with international plug inserts                                                                                                                                                                                                                                                                  | Ratings: IEC 60601-1 (Safety)<br>Cord Length - 6 ft. (1.8 m)<br>VALO X Power Supply is a Medical Grade Class II power supply and provides isolation from MAINS power |
| VALO X Battery Charger        | VALO X Lithium Ion (Li-Ion) Battery Charger:<br>Output - 4.2VDC at 950mA bulk charge current<br>Input - 9VDC at 2A (700mA minimum)<br>Automatic shut off when fully charged<br>Auto-detection of defective cell<br>Protections: Overcharge, Short-circuit, Reverse polarity<br>Amber LED - Charging<br>Green LED - Fully charged<br>Charging time: 1-3 hrs.<br>Ratings: IEC 60601-1 (Safety), CE, WEEE |                                                                                                                                                                      |
| VALO X Li-Ion 14/65 Batteries | Protected, Rechargeable, Li-Ion Battery<br>11MR14/65 3.7V, 900mAh 3.33Wh<br>3.7V, 950mAh Max. 3.52Wh<br>Ratings: WEEE, CE, RoHS, UL 62133-2<br>Cells are rated for >500 charge/discharge cycles before hitting 70% of its capacity. It is advisable to replace batteries every 2-3 years.                                                                                                              |                                                                                                                                                                      |
| Charging Conditions           | Temperature: +10°C to +32°C (+50°F to +90°F)<br>Relative Humidity: 10% to 95%<br>Ambient Pressure: 700 hPa to 1060 hPa                                                                                                                                                                                                                                                                                 |                                                                                                                                                                      |

## Troubleshooting

| If the solutions suggested below do not rectify the problem, please call Ultradent at 800.552.5512. Outside the United States, call your Ultradent distributor or dental dealer. |                    |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------|
| Problem                                                                                                                                                                          | Possible Solutions |

|                                         |                                                                                                                                                                                                                                                                                                                                                                                                             |
|-----------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Light will not turn on                  | •Press the Power Button to wake from Sleep Mode<br>•Confirm battery/cell is fully charged<br>•Check that cell is correctly inserted into the unit<br>•If Amber LED alert appears, the curing light has reached its internal temperature safety limit. Allow the curing light to cool down.<br>•If Red or Amber LED alert flashes continuously, call Ultradent Customer Service for repair                   |
| Light does not stay on for desired time | •Check mode indicator for correct device setting<br>•Check the Low Power indicator for cell charge status<br>•Check that a fresh cell is properly inserted into the unit                                                                                                                                                                                                                                    |
| Light is not curing resins properly     | •Check lens for residual cured resins/composites<br>•Using proper amber UV eye protection, verify the LED lights are working<br>•Check power level with light meter. Refer to maintenance section for specific instructions on using light meter.<br>•Check expiration date on curing resin<br>•Ensure proper technique is being followed according to manufacturer recommendations                         |
| Battery will not charge                 | •Make sure the cell is inserted into the charger in the correct orientation and allow cell to fully charge for up to 3 hours<br>•If Amber lights on the charger do not change to Green, call Ultradent Customer Service to order replacement cells and/or charger<br>•If neither Green nor Amber lights on the charger are visible when a battery is inserted into charger, call Ultradent Customer Service |
| Charger does not charge battery         | •Make sure charger is plugged in and AC power supply is plugged into a working power outlet<br>•If neither Green nor Amber lights on the charger are visible, call Ultradent Customer Service                                                                                                                                                                                                               |

|                                                                                                                                                                                                                                                          |                                              |                       |          |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------|-----------------------|----------|
| Miscellaneous Information                                                                                                                                                                                                                                |                                              |                       |          |
| <b>Electromagnetic Environment Guidance –</b>                                                                                                                                                                                                            |                                              |                       |          |
| a. The VALO X family of curing lights are suitable for use in all establishments, including domestic establishments, that are directly connected to the public low-voltage power supply network (AC MAINS).                                              |                                              |                       |          |
| b. The VALO X family of curing lights can operate from two different power sources.                                                                                                                                                                      |                                              |                       |          |
| i. Battery powered: When the VALO X family of curing lights use the battery, they are not affected by EMI, RF, electric transients, surges, voltage dips, shorts, interruptions, or variations from AC MAINS power.                                      |                                              |                       |          |
| ii. Adaptor powered: the VALO X family of curing lights can use a GlobTek medical grade 9VDC adaptor to operate from AC MAINS power ranging from 100–240 VAC. It provides limited EMI, RF, brown out, and surge suppression.                             |                                              |                       |          |
| c. The VALO X family of curing lights use electrical and electromagnetic energy only for their internal functions. However, any RF emissions are very low and are not likely to cause interference in nearby electronic equipment.                       |                                              |                       |          |
| d. The VALO X family of curing lights are not equipped with a DC power port, Patient Coupling Ports, or Signals Input/Output ports, therefore testing to any of these ports is not applicable. Only compliance with applicable configurations is listed. |                                              |                       |          |
| <b>Table 5.1 Summary of Emission Tests:</b> the VALO X curing light complies with CISPR 11 Group 1, Class B for radiated emissions.                                                                                                                      |                                              |                       |          |
| Port                                                                                                                                                                                                                                                     | Test Description                             | Frequency Range (MHz) | Result   |
| AC Power                                                                                                                                                                                                                                                 | Conducted Emissions (Hot Lead to Ground)     | 0.15 to 30            | Complied |
| AC Power                                                                                                                                                                                                                                                 | Conducted Emissions (Neutral Lead to Ground) | 0.15 to 30            | Complied |
| Enclosure                                                                                                                                                                                                                                                | Radiated Emissions (Vertical Polarity)       | 30 to 1000            | Complied |
| Enclosure                                                                                                                                                                                                                                                | Radiated Emissions (Horizontal Polarity)     | 30 to 1000            | Complied |

|                                                                                |                                                           |          |
|--------------------------------------------------------------------------------|-----------------------------------------------------------|----------|
| <b>Table 5.2 Summary of Immunity Tests (Enclosure Port and AC Power Port):</b> |                                                           |          |
| Basic Standard                                                                 | Environmental Phenomena                                   | Result   |
| EN 61000-4-2                                                                   | Electrostatic Discharge                                   | Complied |
| EN 61000-4-3                                                                   | Radio Frequency Electromagnetic Field                     | Complied |
| EN 61000-4-3                                                                   | Proximity fields from RF wireless communication equipment | Complied |
| EN 61000-4-4                                                                   | Electrical fast transient/burst                           | Complied |
| EN 61000-4-5                                                                   | Lightning Surge (line-line)                               | Complied |
| EN 61000-4-6                                                                   | Radio-frequency continuous conducted                      | Complied |
| EN 61000-4-11                                                                  | Voltage Dips and Interruptions                            | Complied |

|                                                                                  |                              |          |
|----------------------------------------------------------------------------------|------------------------------|----------|
| <b>Table 5.3 Summary of Tests (EN 60601-3-2 and EN 61000-3-3; AC Power Port)</b> |                              |          |
| Basic Standard                                                                   | Environmental Phenomena      | Result   |
| EN 61000-4-2                                                                     | Harmonic Current Measurement | Complied |
| EN 61000-4-3                                                                     | Voltage Flicker Measurement  | Complied |

|                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                       |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>FCC Regulatory Statement</b>                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                       |
| This device complies with Part 15 of the FCC rules. Operation is subject to the following two conditions:<br>1) This device may not cause harmful interference, and<br>2) This device must accept any interference received, including interference that may cause undesired operation.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                               |
| <b>Warning</b>                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                        |
| Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |
| <b>Interference</b>                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                   |
| This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.<br>This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.<br>If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one of the following measures: <ul style="list-style-type: none"><li>• Re-orient or relocate the receiving antenna.</li><li>• Increase the separation between the equipment and receiver.</li><li>• Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.</li><li>• Consult the dealer or an experienced radio/TV technician for help.</li></ul> |
| <b>ISED Regulatory Statement (includes French)</b>                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                    |
| This radio transmitter has been approved by Industry Canada to operate with the antenna types listed with the maximum permissible gain indicated. Antenna types not included in this list, having a gain greater than the maximum gain indicated for that type, are strictly prohibited for use with this device. This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions:<br>1) This device may not cause interference, and<br>2) This device must accept any interference, including interference that may cause undesired operation of the device.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                          |
| This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to ICES-003. These limits are designed to provide reasonable protection against harmful interference in a residential installation.<br>This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.<br>If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one of the following measures: <ul style="list-style-type: none"><li>• Re-orient or relocate the receiving antenna.</li><li>• Increase the separation between the equipment and receiver.</li><li>• Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.</li><li>• Consult the dealer or an experienced radio/TV technician for help.</li></ul>                 |

Report any serious incident to the manufacturer and the competent authority.  
ultradent.com / 1.800.552.5512 / 801.572.4200

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|                                                                                   |                                                                                                                                      |
|-----------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------|
|  | Warning to read and following instructions and pay attention to specific warnings and cautions related to the use of a curing light. |
|  | Refer to Instruction manual/booklet.                                                                                                 |
|  | Compliance with EU regulations for safety, electrical and electronic interference and immunity.                                      |
|  | Authorized representative.                                                                                                           |
|  | Authorized representative.                                                                                                           |
|  | Compliance with Medical Class B electrical safety.                                                                                   |
|  | WEEE compliant: Recycle; do not discard improperly.                                                                                  |

|                                                                                     |                                                                                                                                        |
|-------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------|
|  | Class II Equipment. Wall power supply for the charger is a Class II medical grade, double insulated supply.                            |
|  | Manufactured by Ultradent in accordance with GMP and ISO 13485 practices.                                                              |
|  | Date of Manufacture                                                                                                                    |
|  | Do not re-use to avoid cross contamination                                                                                             |
|  | Upper and lower temperature limitations                                                                                                |
|  | Upper and lower relative humidity limitations: Store or transport the device in areas within the relative humidity range of 10 to 95%. |
|  | Follow local waste and recycling guidelines when disposing of Lithium Ion batteries.                                                   |

|                                                                                                               |                                                                                                                                                               |
|---------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------|
|  <p>500 hPa<br/>1060 hPa</p> | EN - Upper and lower atmospheric pressure limitations: Store or transport the device in areas within the atmospheric pressure range of 1060 to 500 hPa.       |
|                              | Keep Dry: Cargoes bearing this symbol must be protected from excessive humidity and must be stored under cover.                                               |
|                              | Serial Number                                                                                                                                                 |
|                              | Nemko Electrical Safety and EMC                                                                                                                               |
|                              | Ingress Protection Rating                                                                                                                                     |
|                              | Home symbol indicates that the equipment is only for indoor use and not for industrial purposes.                                                              |
|                              | Indicates that the product contains certain toxic or hazardous substances or elements, and can be used safely during its environmental protection use period. |

|                                                                                     |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                        |
|-------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|  | EN -For professional use only DE - Nur zur Anwendung durch den Zahnarzt FR - Pour usage professionnel seulement NL - Alleen voor professioneel gebruik IT - Solo per uso professionale ES - Solo para uso profesional PT - Apenas para utilização profissional SV-Endast för professionell användning DA - Kun til professionel brug FI - Vain ammattikäyttöön EL - Για επαγγελματική χρήση μόνο CS - Pouze pro odborné použití. BG - Само за професионално ползване HR - Samo za profesionalnu uporabu ET - Ainult professionaalseks kasutuseks HU - Kizárólag professzionális felhasználásra GA - Le haghaidh úsáide gairmiúla amháin LV -Tikai profesionālai lietošanai LT - Tik profesionaliam naudojimui MT - Għall-użu professjonali biss RO - Numai pentru uz profesional SK - Len na profesionálne použitie SL - Samo za profesionalno uporabo NO - Kun for profesjonell bruk PL - Wyrób wyłącznie do użytku profesjonalnego. RU - Только для профессионального использования TR - Sadece profesyonel kullanım için AR - للاستخدام من قبل متخصصين فقط TH - สำหรับใช้งาน: เฉพาะผู้เชี่ยวชาญเท่านั้น KO - 전문가 전용 ZH - 仅供专业人士使用 JA - 専門家のみが使用できます |
|-------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|

Manufactured by:

Ultradent Products, Inc.  
505 West Ultradent Drive (10200 South)  
South Jordan, Utah 84095 USA  
Manufactured in the USA from globally sourced materials

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