Dentatus CLASSIC SURTEX® POSTS


PREMIUM DENTAL POSTS WITH PROVEN RELIABILITY

Whatever your crown retention requirements...

DENTATUS CLASSIC SURTEX® POSTS HAS THE SOLUTION:

- Surface treatment for superior retention
- · Innovative and practical assortment box for ease-of-use
- Swedish high quality posts with 100% quality control

The original Dentatus Classic Surtex® Posts have become a global standard in dental anchors, and offer dentists a reliable and simple solution for post and core build-up. Our high-quality, pre-fabricated posts are dependable and economic for most endodontic treatment needs, available in a range of materials including pure titanium, stainless steel and gold-plated.

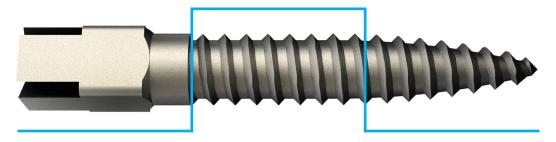

Square crosscut head for core material adaptation and retention.

Parallel sides for passive insertion with thread pitch for secure retention to all cements.

Anatomical geometry for safe extra length of posts and for tapered roots.



What ever your crown retention requirements, Dentatus Classic Surtex Post system has the solution. The system is thoroughly tested for retention of core materials in endodontically treated teeth. The posts are cemented permanently or temporarily in the prepared root canal and will offer improved retention and a non-glare appearance due to its texturized surface. The system consists of reamers and specially fitted posts in 6 different diameters and 5 different lengths.

The anatomically shaped reamers and posts correspond to the natural morphology of the root canal to minimize risk for root fractures.

SUPERIOR MATERIALS FOR UNPARALLELED QUALITY

<u>The Titanium posts</u> are made of pure titanium and are recommended for most retentive purposes. Titanium is biologically highly compatible with surrounding tissue.

<u>The Stainless Steel posts</u> are recommended for all retentive purposes where extra high shear strength is required.

<u>The Gold Plated posts</u> are suitable for most retentive purposes, especially where an economical alternative is required (e.g. in temporary crowns).



SIZES Available in Titanium, Gold Plated and Stainless Steel

	SHORT (7,8 MM)	MEDIUM (9,3 MM)	LONG (11,8 MM)	EXTRA LONG (14,2 MM)	SUPER LONG (17 MM)
1	Ø 1,05 mm	Ø 1,05 mm	Ø 1,05 mm		
2	Ø 1,20 mm	Ø 1,20 mm	Ø 1,20 mm		
3	Ø 1,35 mm	Ø 1,35 mm	Ø 1,35 mm	Ø 1,35 mm	Ø 1,35 mm
4	Ø 1,50 mm	Ø 1,50 mm	Ø 1,50 mm	Ø 1,50 mm	Ø 1,50 mm
5	Ø 1,65 mm	Ø 1,65 mm	Ø 1,65 mm	Ø 1,65 mm	Ø 1,65 mm
6	Ø 1,80 mm	Ø 1,80 mm	Ø 1,80 mm	Ø 1,80 mm	Ø 1,80 mm

OUTSTANDING ASSORTMENT IN AN INNOVATIVE BOX

Surtex posts are delivered as complete assortments of posts, reamers and accessories in the One-Stop Dispenser. Equiped with built-in gauges, factory mounted safety seal and anti-slip bottom to prevent unwanted sliding on table top. All sizes of posts and reamers are also available in refills.

Safety Seal

Tamper proof

from factory

Full range of refills available


Measuring of length



Easy diameter measuring (beneath safety seal)

Click-stop, see-through lid

Compartment floor designed for easy catch of indicator rings

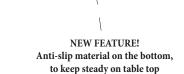

Scan and read more about Dentatus Post Systems

