

The complete system for tooth isolation.



Please call 1-800-343-5336 for quarterly specials

IVORY[®] Rubber Dam



Made from natural latex, thoroughly washed in the manufacturing stage GREEN

to reduce surface proteins. Ivory dam has extremely high tear strength, reducing waste from tearing during application. The dam is thin and pliable, allowing easy placement in especially difficult areas. The bright blue or green tint increases visibility by its contrast.

Available in 6x6 or 5x5 (medium weight).



IVORY® Punch

This precision instrument is made from high quality stainless steel. The unique floating plunger strikes all edges of the cutting disc evenly, guaranteeing a perfect, tear-resistant hole every time. The cutting disc rotates through a range of 6 hole sizes, providing the correct opening for the smallest anterior tooth to the largest molar.

IVORY® Forceps

The lvory forceps are made from high grade stainless steel, which will assure many years of use. The forceps will reach around the bow of any clamp without tilting, which is especially important when reversing a clamp. The straight head design allows easy gripping and placement of any clamp size.

IVORY® Rubber Dam Kits

The Ivory Clamp Starter Kit offers a selection of the most widely used clamps. These popular clamps provide a variety of general use clamps that fulfill most needs. The Rubber Dam Complete Kit is a self-contained unit providing all the necessary equipment for the Rubber Dam Technique. The Complete Kit is especially useful for the practitioner just beginning their practice or for one who is implementing the technique for the first time.

Kit components include:

8 Clamp Starter Kit

One each of the following clamps:

0 2A 8A 14 2 7 9 14A

Available in Shiny or Matte Finish

Complete Kit

(1) Box 6x6 Blue Rubber Dam
(1) 6x6 Marking Template
(1) 6 Inch Frame
(8) Most Popular Clamps (O 2A 8A 14 2 7 9 14A)
Clamp Forceps
Rubber Dam Punch

Ivory Rubber Dam Clamps by Heraeus

Heraeus manufactures many styles of uniquely designed IVORY[®] Rubber Dam clamps to meet the most demanding situations. You will find one of these clamps suitable for practically every application. The tension of each clamp is meticulously hand-set in our factory and is appropriate for the tooth on which it is designed to be used. Ivory Clamps are available in shiny and matte finish as noted in the charts.

Listed below is the complete Ivory Clamp line and suggestions for routine procedures.

Most Popular Clamps

PREMOLARS – WINGED

- 00 Small lower, anterior
- 0 Elongated small neck
- 1 General purpose upper
- 2 General purpose lower
- 2A Large neck bicuspid

PREMOLARS – WINGLESS

W2 General purpose lower

W2A Large neck bicuspid

MOLAR - WINGED

- 3 Small lower
- 4 Small upper for round teeth
- 7 General purpose lower
- 8 General purpose upper
- 8A Upper and lower molar root
- 12A Right partially erupted, 3rd molar
- 13A Left partially erupted, 3rd molar
- 14 Partially erupted, most popular clamp

MOLAR – WINGED (CONTINUED) 14A Large, partially erupted 56 All large molars

MOLAR – WINGLESS

W3 Small lowers W8A Upper and lower molar root W14A Large partially erupted

Specialty Clamps

PREMOLARS – WINGED 1A Root clamp

PREMOLARS – WINGLESS

- W0 Elongated small neck
- W1 General purpose upper
- 27N Universal

MOLAR - WINGED

- 5 Large upper
- 7A Large lower
- 10 Right, small short crown

MOLAR – WINGED (CONTINUED)

- 11 Left, small short crown
- 26 Small lower molar
- 27 Distal access, large molar

MOLAR - WINGLESS

- W4 Small upper
- W5 Large upper
- W7 General purpose lower, square teeth
- W8 General purpose upper, round teeth
- W14 Partially erupted

MOLAR – WINGLESS (CONTINUED)

- W56 Large molar
- 26N Universal upper
- 28N General purpose lower molar

LABIAL - WINGED

- 6 Central and incisor
- 9 Labial prep, all anteriors

LABIAL - WINGLESS

212 Labial prep, cervical

Labial



For labial cavities on centrals and cuspids. For second bicuspids and second upper molars.





For labial cavities on all anterior teeth.




For labial cavities. Anterior cervical clamp.

Tooth Diagram

Tooth Clamp

- **1** 4, 5, 8, 8A, 12A, 14, 26N
- 4, 5, 8, 8A, 12A, 14, 26N
 4, 5, 8, 8A, 12A, 14, 14A,
- 26N, 56
- **4** 0, 1, 1A, 2A, 27N
- **5** 0, 1, 1A, 2A, 27N
- 6 6, 9, 212, 212SA
- **7** 6, 9, 212, 212SA
- **8** 6, 9, 212, 212SA
- **9** 6, 9, 212, 212SA
- **10** 6, 9, 212, 212SA
- 11 6, 9, 212, 212SA
- **12** 0, 1, 1A, 2A, 27N
- **13** 0, 1, 1A, 2A, 27N
- 14 4, 5, 8, 8A, 13A, 14, 14A, 26N
- **15** 4, 5, 8, 8A, 13A, 14, 14A, 26N
- **16** 4, 5, 8, 8A, 12A, 14, W14, 14A, 26N, 56

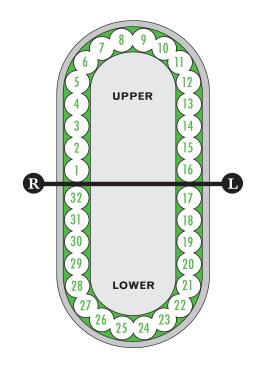

Tooth Clamp

- **17** 3, 7, 7A, 8A, 12A, 14, 14A, 27, 56
- **18** 3, 7, 7A, 8A, 12A, 14, 14A, 27, 56
- **19** 3, 7, 7A, 8A, 12A, 14, 14A, 27, 56
- **20** 0, 1A, 2, 2A, 6, 27N
- **21** 0, 1A, 2, 2A, 6, 27N
- **22** 00, 6, 9, 212, 212SA
- **23** 00, 6, 9, 212, 212SA
- **24** 00, 6, 9, 212, 212SA
- **25** 00, 6, 9, 212, 212SA
- **26** 00, 6, 9, 212, 212SA
- **27** 00, 6, 9, 212, 212SA
- **28** 0, 1A, 2, 2A, 6, 27N
- **29** 0, 1A, 2, 2A, 6, 27N
- **30** 3, 7, 7A, 8A, 12, 12A, 14, 14A, 27, 56
- **31** 3, 7, 7A, 8A, 12, 12A, 14, 14A, 27, 56
- **32** 3, 7, 7A, 8A, 14, 14A, 27, 56

Bicuspid



Universal bicuspid clamp. Jaws conform to curvature at gum line. Available wingless only



For small lower anteriors and teeth set irregularly in arch. High bow, narrow jaws.



General purpose lower. Same as #1, only jaws are flat so as not to impinge on tissue. Also available wingless - W2





General purpose upper. Jaws fit curvature at gum line. Also available wingless - W1





For bicuspids with large necks. Broad jaws. Also available wingless - W2A

* Available wingless only

- Also available wingless
- T Also available with serrated jaws (Tiger clamps)
- A Slightly larger
- N No wing

Tiger Clamps

- 1T Bicuspid, general purpose upper
- 2T Premolar, general purpose lower
- 2AT Premolar, large neck bicuspid
- 9T Anterior labial
- 14T Molar, partially erupted
- 14DT Molar, distal access, serrated
- 56T For large molars

Accessories

IVORY® RUBBER DAM CLAMP STARTER KIT See page 5 for contents


IVORY® RUBBER DAM COMPLETE KIT See page 5 for contents

IVORY® RUBBER DAM Clamp Board

25 Peg, unnumbered, plastic or stainless steel

IVORY® RUBBER DAM TEMPLATE 6X6

Molar

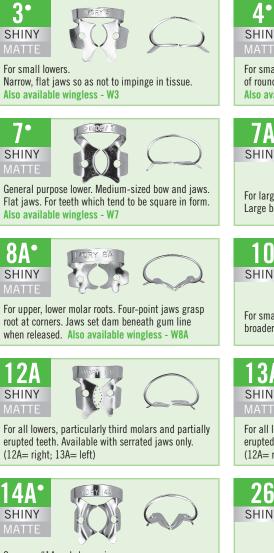





For large molar distal cavities.



For small uppers. Narrow jaws designed to fit necks of round-shaped teeth. Also available wingless - W4



For large lowers. Same as #7, only larger. Large bow, broad jaws.



For small molars with short crowns. Buccal jaw broader than lingual. (10=right; 11=left)



For all lowers, particularly third molars and partially erupted teeth. Available with serrated jaws only. (12A= right; 13A= left)



For small lower molars.



General purpose lower molar clamp. Available wingless only



For large uppers. Jaws conform to curvature at gum line. Also available wingless - W5





General purpose upper. For round upper teeth. Jaws fit curvature at gum line, to prevent rotating. Also available wingless - W8



For small molars with short crowns. Buccal jaw broader than lingual. (10=right; 11=left)



For partially erupted, undersized, round teeth or for teeth sloping from gum line. Small bow. Also available wingless - W14





Universal upper with festooned jaws. Available wingless only





For all large molars. Also available wingless - W56 56 and W56 are SHINY, 56T is MATTE

