




INFORMATIONS


Definition: Syringe Sleeve

Indication: Diminishing of anxiety caused by needle fear during anesthesia.

Presentation: Package with 2 (two) green alligators, 1 (one) pink and 3 (Three) adhesive labels.

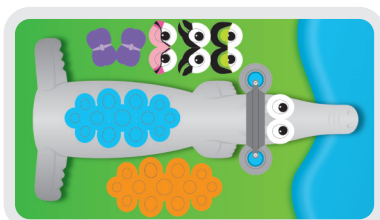

Adhesive labels: For the customization of the alligator

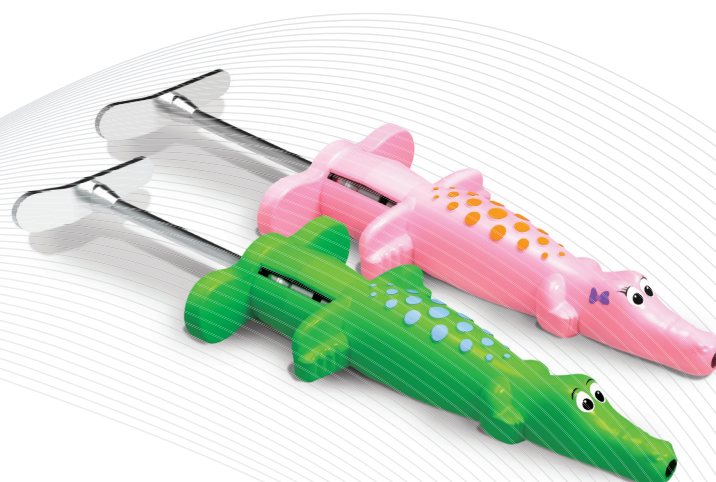


ADVANTAGES

- Facilitates the work of the dentist during anesthesia application.
- It may be given to the patient as a nice gift.
- Sterilizable – Guarantee of Biosafety.



INSTRUCTIONS FOR USE

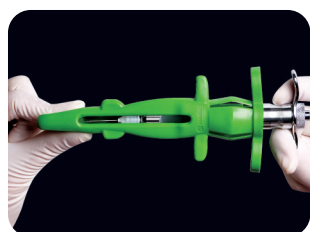

Prepare the syringe and insert it in the sleeve, the needle should not touch the internal wall of the sleeve. Ensure that the cartridge is clearly visible at the lower opening of sleeve.



Position index and middle fingers on the rear feet of the Alligator Sleeve and place thumb on the plunger. At this moment, the needle is hidden inside the "mouth" of the alligator.



Place the mounted assembly into the area of the mouth to receive anesthesia. The lower part of the Alligator Sleeve must be positioned upwards; this allows the cartridge to be seen by the dentist so that the flow of the anesthetic solution and any resultant bleeding can be checked.



Press the rear feet of the Alligator Sleeve with the index and middle fingers to expose the needle. Proceed with the injection of the anesthetic solution. After completing the injection, release the rear feet of the sleeve, which will hide the needle again.

See the video of the reaction of the child that saw the syringe with and without the sleeve: www.angelus.ind.br

The Angelus Syringe Sleeve, tranquility during the application of anesthetics.




SIMPLIFYING PROCEDURES
THROUGH INNOVATION.

WWW.ANGELUS.IND.BR
PHONE +55 43 2101.3200
FAX +55 43 2101.3201

DO YOU WANT TRANQUILITY DURING ANESTHESIA?

START USING AN ANGELUS
ALLIGATOR SYRINGE SLEEVE

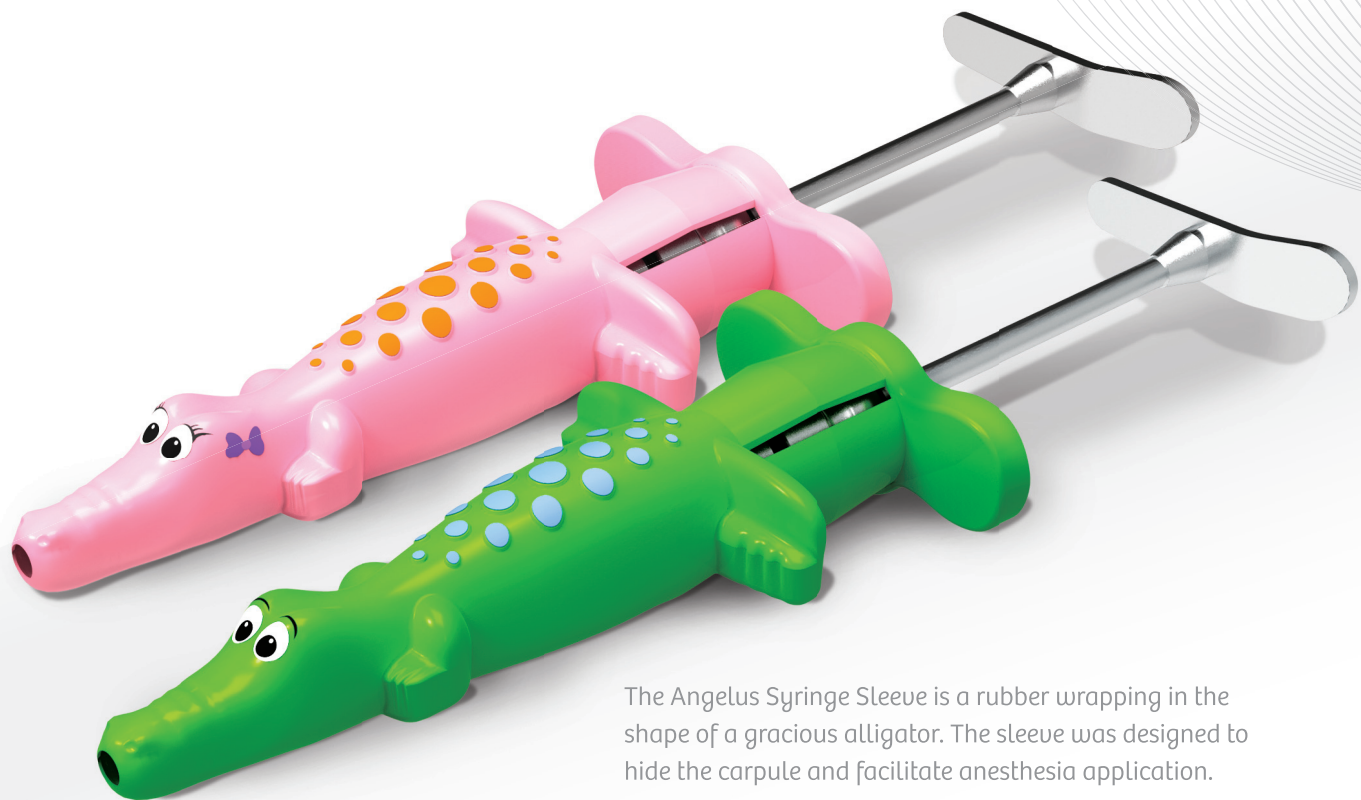

The Angelus Syringe Sleeve is a rubber wrapping in the shape of a gracious alligator. The sleeve was designed to hide the carpule and facilitate anesthesia application.

 **angelus**[®]
science and technology

