



# **Carbopost<sup>®</sup> / Snowlight<sup>®</sup> / Snowpost<sup>®</sup>**

COMPOSITE POSTS - Strength and Flexibility

# Carbopost®

## REPLACES METALS WHERE IT MATTERS.

The fiber reinforced composite post is as strong as metal and reduces shear stress by 50% by transmitting stress progressively and regularly to the entire coated surface.

Made of 60% T700 carbon fibers and 40% epoxy resin Carbopost has proven to have the same mechanical behavior as metal posts.

- Shock absorbing
- Radiopaque
- Silanated for superior adhesive qualities
- Autoclavable up to 135°C/ 275°F
- Bio-compatible and corrosion free
- Anatomical shape and 5 different diameters (0.8mm – 1.6mm)

### MECHANICAL PROPERTIES


	Carbopost
Elastic Modulus	141 GPa
Shear Strength	23.9 MPa

### Carbopost Kits

C801001	Basic Kit : 50 assorted posts (1.0, 1.2, 1.4, 1.6) and 4 matching finishing drills
C801002	Intro Kit: 20 assorted posts (1.0, 1.2, 1.4, 1.6) and 4 matching finishing drills

### Carbopost Refills

10 posts/pack		10 posts/pack with matching finishing drill	
C802013	Ø 0.8		
C802001	Ø 1.0	C802006	Ø 1.0
C802002	Ø 1.2	C802007	Ø 1.2
C802003	Ø 1.4	C802008	Ø 1.4
C802004	Ø 1.6	C802009	Ø 1.6

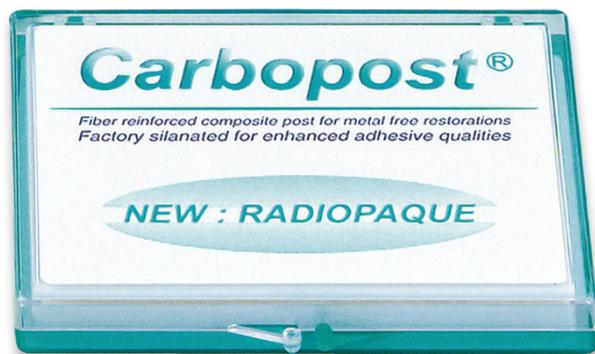

# Snowpost® Snowlight®

## STRENGTH AND FLEXIBILITY

With their elastic modulus and flexural strength close to dentin, Snowpost and Snowlight follow the movement of the tooth, preventing tensions, fractures, and other problems.


Made of zircon-rich glass fiber embedded in a resin matrix the posts ensure the perfect balance of strength and flexibility.

## EXCELLENT AESTHETICS

No grey shadows. No show through. The white root canal posts answer all demands of aesthetic dentistry.

## RADIOPAQUE AND LIGHT TRANSMITTING

Best-in-class radiopacity to view post and placement in tooth. Snowlight posts allow transmission of light, thus are ideal for all dual or light cured materials.

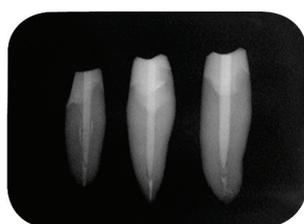

- Superior fatigue resistance
- Radiopaque and light transmitting
- Silanated for superior adhesive qualities
- Autoclavable up to 135°C/ 275°F
- Bio-compatible and corrosion free
- Anatomical shape and 5\* different diameters (0.8mm\* – 1.6mm) \*Snowlight posts only

### MECHANICAL PROPERTIES

	Snowpost	Snowlight
Elastic Modulus	45 GPa	52 GPa
Shear Strength	38 MPa	64 MPa



### Snowpost Kits

C810001	Basic Kit : 50 assorted posts (1.0, 1.2, 1.4, 1.6) and 4 matching finishing drills
C810002	Intro Kit: 20 assorted posts (1.0, 1.2, 1.4, 1.6) and 4 matching finishing drills

### Snowpost Refills

10 posts/pack		10 posts/pack with matching finishing drill	
C811001	Ø 1.0	C811005	Ø 1.0
C811002	Ø 1.2	C811006	Ø 1.2
C811003	Ø 1.4	C811007	Ø 1.4
C811004	Ø 1.6	C811008	Ø 1.6



### Snowlight Kits

C830001	Basic Kit : 50 assorted posts (1.0, 1.2, 1.4, 1.6) and 4 matching finishing drills
C830002	Intro Kit: 20 assorted posts (1.0, 1.2, 1.4, 1.6) and 4 matching finishing drills

### Snowlight Refills

10 posts/pack		10 posts/pack with matching finishing drill	
C831019	Ø 0.8		
C831001	Ø 1.0	C831006	Ø 1.0
C831002	Ø 1.2	C831007	Ø 1.2
C831003	Ø 1.4	C831008	Ø 1.4
C831004	Ø 1.6	C831009	Ø 1.6



# STAINBUSTER®

## THE INNOVATIVE, PAINLESS COMPOSITE BUR

The innovative composite bur made of zircon-rich glass fiber features permanent abrasive power for minimal invasive dentistry. Stainbuster removes cement and stains from tooth surface with no damage to enamel, ceramics or soft tissues - it simply slides over the soft tissue without cutting or grinding.

The self-sharpening burs continuously maintain their abrasive power.

- Minimally invasive
- No damage to enamel, ceramics, cement, or soft tissue
- Self-Sharpener
- Pain free

**DAILY PRACTICE:** Stain removal, cement or composite surplus grinding, temporary cement removal, light maintenance scaling

**ORTHODONTICS:** Elimination of adhesive cement remnants after bracket removal

**PERIODONTICS:** Scaling, root surfacing, easy access to furcation

**IMPLANTOLOGY:** Sealing cement surplus removal even in narrow and difficult-to-reach places



Effective root planing and surfacing




Easy calculus removal



Gentle and efficient implant cleaning

### AVAILABLE SHAPES

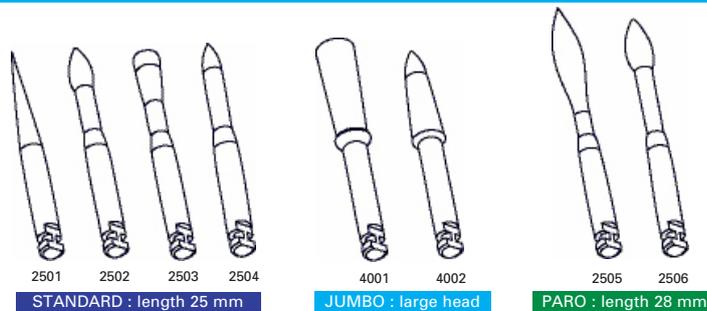

#### Stainbuster Assorted Pack - Contains 6 burs

2540 One each of 2501, 2502, 2504, 2505, 2506, 4002

#### Stainbuster Individual Packs - Contains 6 burs

2501	2503	4001	2505
2502	2504	4002	2506



**abrasive technology**

[www.abrasive-tech.com](http://www.abrasive-tech.com)

8400 Green Meadows Dr. • P.O. Box 545  
Lewis Center, OH 43035 USA  
740.548.4100 • Fax: 740.548.7617

Two Striper® diamond dental burs are manufactured under U.S. Pat. No: 5,511,718  
©Abrasive Technology 2/09