

WORKING LENGTH DETERMINATION & ROOT CANAL TREATMENT



**MicroMega**

**Dual Move**

Cordless endo motor



**Dual Pex**

Electronic apex locator




 **COLTENE**

## CORDLESS ENDO MOTOR WITH CONTROL OF MOVEMENT, SPEED & TORQUE

### GREAT SHAPING EXPERIENCE

- > Well-balanced in hand
- > Excellent visibility
- > Simple to operate
- > Adaptative, ideal for both left and right-handed practitioners

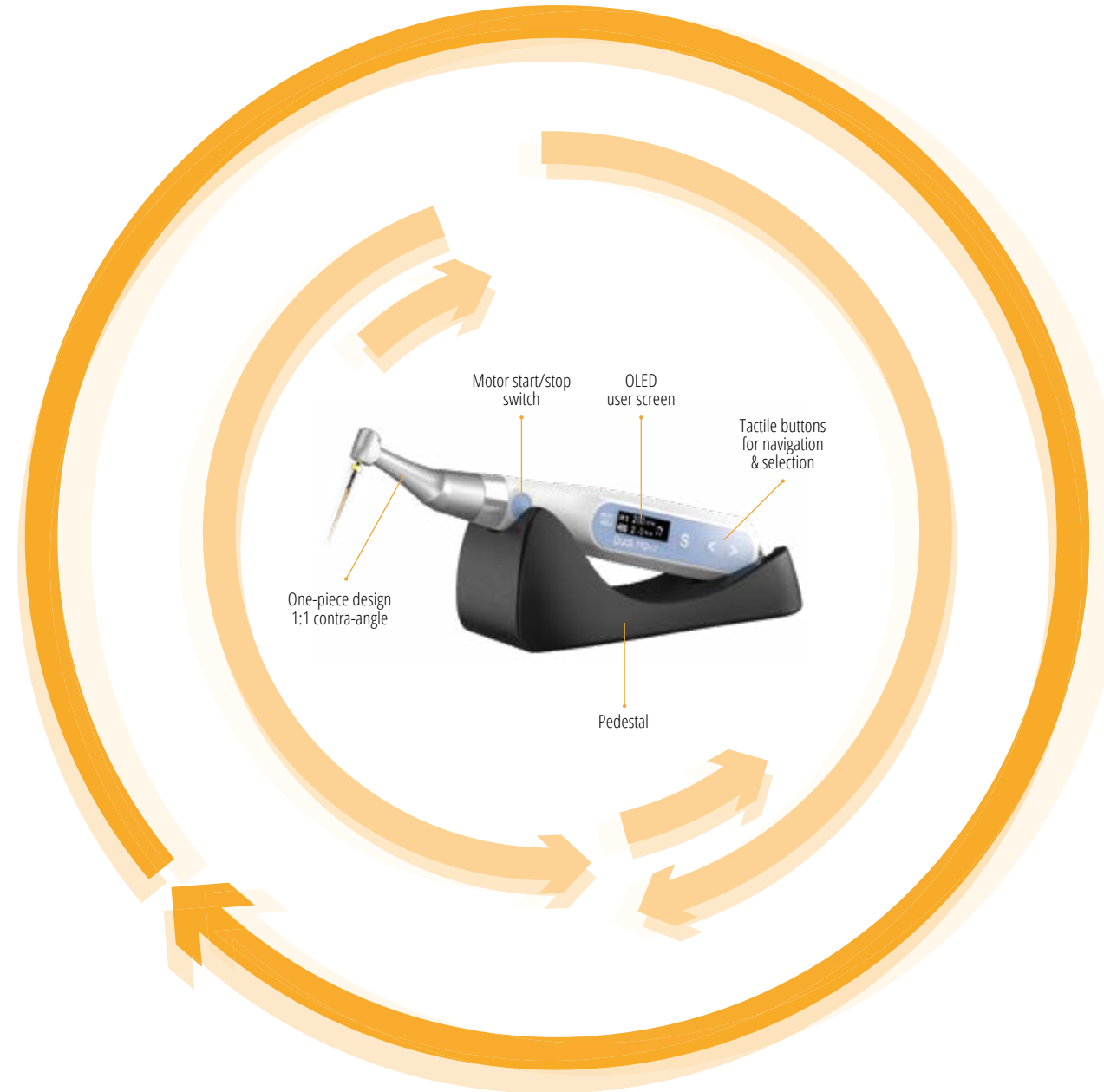

### SAFE ENDO TREATMENT

- > Prevention of file separation
- > Speed and torque control

### ADAPTED TO INDIVIDUAL NEEDS

- > Flexible settings with large number of parameters
- > Suitable for all NiTi files
- > Versatile, stand-alone use or combined with Dual Pex
- > Dual movement : continuous / reciprocating\*

\* Not available in some countries.

# Dual·Move

## GREAT SHAPING EXPERIENCE

### > FREEDOM WHEN SHAPING!

- Slim and compact design
- Light and well-balanced handpiece



## ADAPTED TO INDIVIDUAL NEEDS

### > ALL LEVELS OF PRACTICE, ALL NiTi FILES!

- 6 customisable memory programmes with speed, rotation direction, torque limit
- Continuous rotation
- Reciprocating movement\*




### > TURN THE CONTRA-ANGLE TO THE PERFECT POSITION AND GET EXCELLENT VISIBILITY!

- 340° rotative contra-angle
- Miniature head and thin neck



## SAFE ENDO TREATMENT

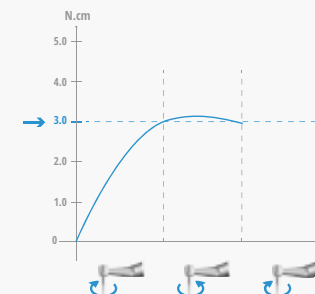

### > PREVENT FILE SEPARATION!

- Auto-reverse on torque limit
- Alarms with adjustable volume level

### > CONTROL YOUR PREPARATION WITH VISION AND SOUND!

- Ambidextrous display (rotation of 180°)
- Automatic locking of the file on the contra-angle
- Perfect position of start/stop switch for single hand operation
- User interface with short learning curve



### > CONTROL YOUR TREATMENT!

- Reliable speed and torque control
- On-demand calibration function
- Instantaneous torque value on the screen and audible sounds



\* Not available in some countries.

## STAND ALONE OR COMBINED WITH DUAL MOVE


### FOCUSED ON THE ESSENTIALS

- > User friendly screen
- > Accurate localisation even in wet and dry root canals, owing to multi-frequency positioning technology
- > Stable delivered voltage whatever the battery level
- > Real-time visual and sound indications
- > Auto power-off to save the battery
- > Plug & Play self tester




### YOUR PERSONAL COMPANION

- > Compact and robust design
- > Easy-to-read 3.5" colour LCD display
- > Stand-alone or combined with Dual Move
- > Rechargeable battery
- > Adjustable apical limit (Reference Point) from 06 (apical constriction) to 00 (apical foramen)

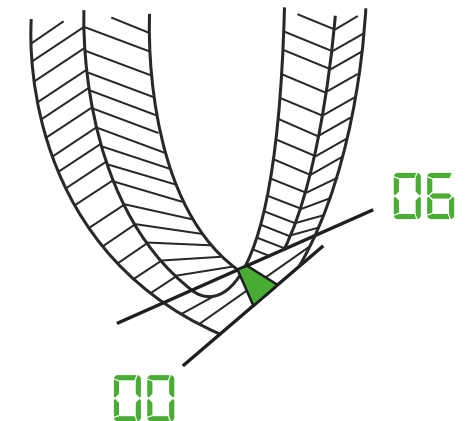

# THE PERFECT DUO WITH DUAL MOVE

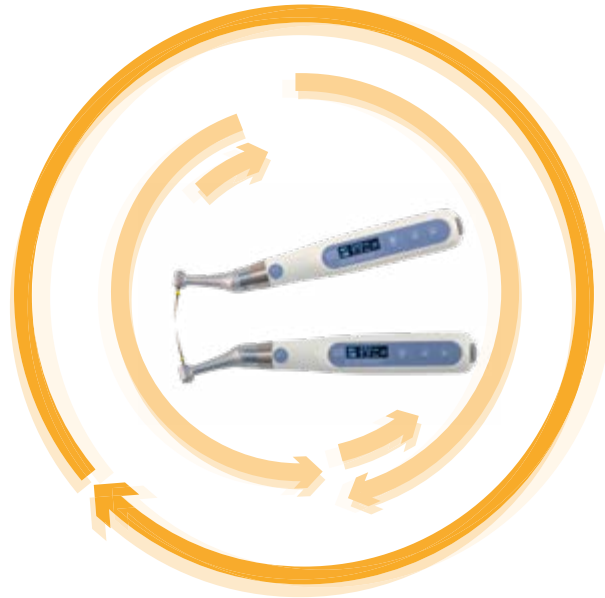

Move

Dual



Pex



ALL READINGS AT A GLANCE AND CLOSE TO YOU!

> Screen mirroring



> Duplication of the file location on Dual Move's screen

CHOOSE THE WAY TO WORK WITH ADDITIONAL FEATURES!

> Auto Start & Stop



> Auto-start when the file enters the root canal

> Auto-stop when it leaves the canal

> Apical Reverse



> Auto-reverse when the files reaches the specified apical limit

> Automatic restart as soon as the files comes back above the apical limit

> Apical Slow Down



> Automatic speed reduction when the file reaches your preset apical limit



« Dual Move, Dual Pex, 2Shape TS1 and 2Shape TS2: your 4 cardinal points for a safe navigation of the root canal system »

*Dr. Filippo CARDINALI - Italy*



« Dual Move is well-balanced and easy-to-use. Dual Pex is really accurate. Combine this cordless motor with a prebendable, efficient, single file like One Curve and endo becomes so easy! »

*Dr. Tara Mc MAHON - Belgium*

#### Dual Move and Dual Pex Presentation

**CE** 0197 Year of CE marking : 2017

Medical devices class IIA according to directive 93/42/EEC.

Notified body: TÜV Rheinland.

Medical device for dental care, meant for professional dental use only.

See product labelling and instructions for use.

 Manufactured by  
Changzhou Sifary Medical  
Technology Co., Ltd.  
No. 99, Qingyang Road,  
Xuejia County, Xinbei District  
Changzhou City  
213000 Jiangsu  
China

**EC REP** Landlink GmbH  
Dorfstrasse 2/4, 79312,  
Emmendingen - Germany

Distributed by Micro-Mega SA  
12, rue du Tunnel  
25000 BESANCON  
France

#### One Curve Presentation


**CE** 0459 Year of CE marking : 2017

Medical device class IIA according to directive 93/42/EEC.

Notified body: LNE/G-MED.

Medical device for dental care, meant for professional dental use only.

See product labelling and, where applicable, instructions for use.

 Micro-Mega SA  
12, rue du Tunnel  
25000 BESANCON / France

#### 2Shape Presentation


**CE** 0459 Year of CE marking : 2017

Medical device class IIA according to directive 93/42/EEC.


Notified body: LNE/G-MED.

Medical device for dental care, meant for professional dental use only.

See product labelling and, where applicable, instructions for use.

 Micro-Mega SA  
12, rue du Tunnel  
25000 BESANCON / France

© COLTENE – MicroMega.Coltene.com



**Micro-Mega SA**  
12, rue du Tunnel  
25000 BESANCON / France  
T +33 3 81 54 42 42  
F +33 3 81 54 42 30  
commercial@micro-mega.com

