



THE ORIGINAL KNITTED CORD


with Ultrapak® CleanCut

Easy packing. Excellent absorption. Exceptional retraction.




Ultrapak® Knitted Cord


Available in six versatile sizes with and without epinephrine

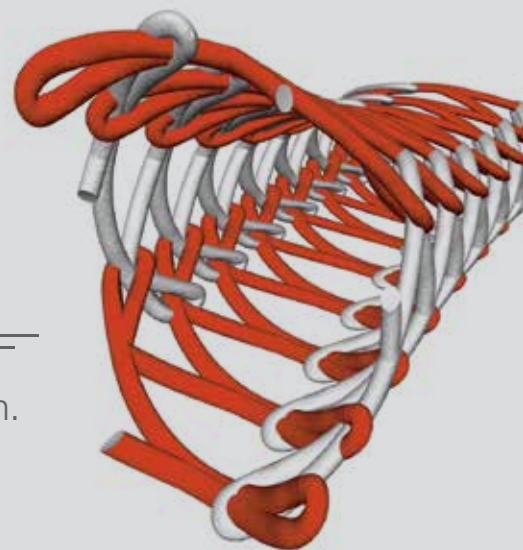

When tissue management is key, reach for Ultrapak knitted displacement cord. Its proprietary, knitted design works better than braided or twisted cord at **packing** into the sulcus, **absorbing** fluids, and **expanding** for adequate retraction. As an important part of Ultradent's proven tissue management system, Ultrapak ensures maximum control with minimum effort.

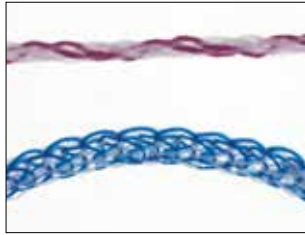

Knitted Ultrapak cord (shown at bottom) is composed of thousands of tiny, interlocking loops so it compresses and expands easier than other cords. 100% cotton fibers provide high absorption of hemostatic agents and sulcular fluids.

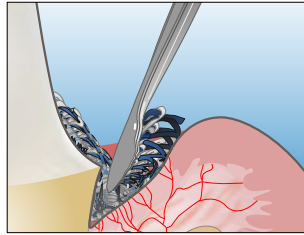

Ultrapak cord packs effortlessly into the sulcus, then exerts a gentle, continuous outward force as knitted loops seek to open. Optimal tissue displacement occurs in 1-3 minutes.



Ultrapak packing cord does not entangle in burs or packers. Bright color facilitates easy location and removal.



Ultrapak CleanCut design features a blade in the cap for efficient cutting. A redesigned dispensing orifice prevents cord from falling into the bottle.

#137 - Ultrapak Cord #000 Refill



#132 - Ultrapak Cord #1 Refill

#199 - Ultrapak Cord #1 Refill - Epinephrine



#136 - Ultrapak Cord #00 Refill

#197 - Ultrapak Cord #00 Refill - Epinephrine



#133 - Ultrapak Cord #2 Refill

#200 - Ultrapak Cord #2 Refill - Epinephrine



#131 - Ultrapak Cord #0 Refill

#198 - Ultrapak Cord #0 Refill - Epinephrine



#134 - Ultrapak Cord #3 Refill



#130 - Ultrapak Kit

- 1 - each #00, 0, 1, and 2 Cord
- 1 - Ultrapak organizer



Each bottle contains 8fl/244cm of cord.

Ultrapak E is impregnated with racemic epinephrine hydrochloride and has a shelf life of 24 months.



Recommended Product Fischer's Ultrapak Packers

Specially designed packers ease the packing of Ultrapak knitted cord. Available with thin, serrated or non-serrated edges, they allow practitioners to place the cord in one circular, sliding motion, without lifting the instrument from the cord.



Ultradent Products, Inc.

TISSUE MANAGEMENT

— Since 1978 —

800.552.5512 ULTRADENT.COM

©2012 Ultradent Products, Inc. All Rights Reserved.

87577.4 121912

ULTRADENT
PRODUCTS, INC.

Improving Oral Health Globally