

ORTHODONTIC



ORTHO LINGUAL COLLECTION

REDEFINING THE ART
OF ORTHODONTICS

How the best perform



Introducing the Hu-Friedy ORTHO LINGUAL COLLECTION

14% TO 29% THINNER
ORBIT FORMED JOINT
THAN STANDARD

The Ortho Lingual Collection consists of ten instruments that feature:

- Optimal **intraoral access**
- Ideal **visibility** for clinician
- Enhanced **comfort** for patient
- More **time savings** allowing for additional patient education

Every ortho instrument from Hu-Friedy is created using the highest quality materials, processes and quality assurance.

DID YOU KNOW?

Hu-Friedy Orthodontic instruments are designed to withstand the rigors of steam sterilization while maintaining their beauty and functionality.

1.2" TO 2" LONGER
THAN STANDARD
HANDLE & 3/4" TO
1.2" LONGER THAN
LONG HANDLE

LINGUAL CUTTERS

Lingual Distal End Cutter

| 678-700

Wire size: .012" to .018"




Occlusal view of the Lingual Distal End Cutter cutting the excess archwire distal to the buccal tube.

Lingual Pin & Lig Cutter

| 678-701

Wire size: < = .012"

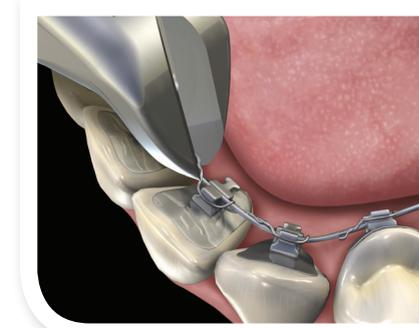

Occlusal view of the Lingual Pin & Lig Cutter cutting the excess steel ligature.

LINGUAL PLIERS



Lingual Weingart Pliers

| 678-702

Wire size: < = .030"



Occlusal view of the Lingual Weingart Pliers guiding the archwire into the bracket slot.



Lingual Weingart Pliers,
90° Bend

| 678-702-90

Wire size: < = .030"

LINGUAL PLIERS



Lingual Mathieu Pliers

| 678-705




Occlusal view of the Lingual Mathieu Pliers tying steel ligature to hold the arch wire into bracket slot.



Lingual Bracket
Removing Pliers

| 678-703

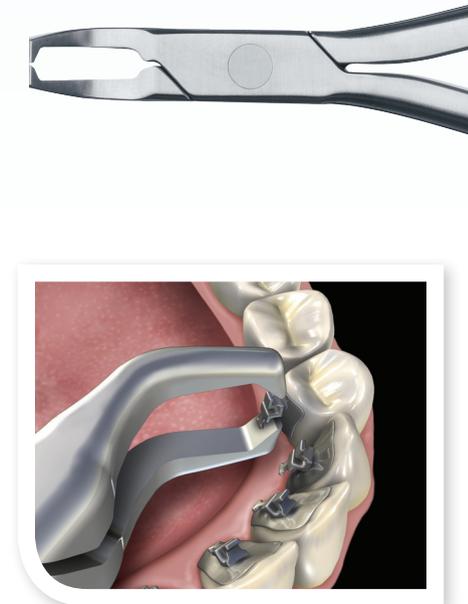

Lingual view of the Lingual Bracket Removing
Pliers debonding the bracket.

LINGUAL PLIERS



Lingual Hammerhead
NiTi Pliers

| 678-704

Wire size: $\leq .025''$



Occlusal view of the Lingual Hammerhead NiTi Pliers bending NiTi wire distal to the buccal tube.



Lingual Hammerhead NiTi
Pliers with Serrations

| 678-704SER

Wire size: $\leq .025''$




LINGUAL HAND INSTRUMENTS

Lingual Ligature
Director

| 678-706

Wire size: $\leq .021'' \times .025''$

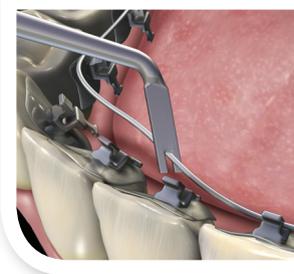

Ligature Side: Occlusal view of the
Lingual Ligature Director guiding the
archwire into the bracket slot.



Picker Side: Occlusal view of the
Lingual Ligature Director tucking
the excess steel ligature under-
neath the bracket.



Lingual Ligature
Director, U Shape

| 678-707

Wire size: $\leq .021'' \times .025''$




Occlusal view of the U Shaped
Lingual Ligature Director guiding
the archwire into the bracket slot.



Lingual Torquing
Key with Picker

| 678-708

Slot size .018''

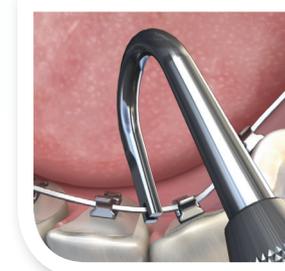

Occlusal view of the Torquing
Key guiding the archwire into
the bracket slot.



RECYCLE BROKEN & OLD INSTRUMENTS

HU-FRIEDY
PROGRAMS



When your instrument blades cannot be sharpened anymore, take advantage of Hu-Friedy's Envirodent® Instrument Recycling Program. Your old, broken instruments will be recycled to help create new car bumpers, bridges and street lights. For every 12 you recycle, we'll give you 1 FREE Hu-Friedy instrument up to \$42!



To start recycling visit
Envirodent.com

THE HU-FRIEDY ORTHO LINGUAL COLLECTION

CUTTERS

Lingual Distal End Cutter	678-700
Lingual Pin & Lig Cutter	678-701

PLIERS

Lingual Weingart Pliers	678-702
Lingual Weingart Pliers, 90°	678-702-90
Lingual Bracket Removing Pliers	678-703
Lingual Hammerhead NiTi Pliers	678-704
Lingual Hammerhead NiTi Pliers with Serrations	678-704SER
Lingual Mathieu Pliers	678-705

HAND INSTRUMENTS

Lingual Ligature Director	678-706
Lingual Ligature Director, U Shape	678-707
Lingual Torquing Key with Picker	678-708

For more information about
Hu-Friedy's Total Solutions visit
HU-FRIEDY.COM

PERIODONTAL



RESTORATIVE



ORTHODONTIC



INSTRUMENT
MANAGEMENT



ENDODONTIC



SURGICAL



ULTRASONIC
SCALING



Hu-Friedy Mfg. Co., LLC: 3232 N. Rockwell St. | Chicago, IL 60618 | USA
Telephone: 1-800-Hu-Friedy (1-773-975-6100) | Website: Hu-Friedy.com

©2014 Hu-Friedy Mfg. Co., LLC. All rights reserved. HF-899/0614

Hu-Friedy