

# Order Information

EDM NITI FILE SYSTEM

## Shaping

HyFlex EDM OGSF Set, 4 pcs	21 mm	25 mm
OGSF Sequence	60036141	60036142
HyFlex EDM OGSF Refills, 3 pcs	18 mm	21 mm 25 mm
Orifice Opener 18/.11	60036140	-
Glidepath File 15/.03	-	60023493 60023492
Shaping File 18/.045	-	60036136 60036137
Finishing File 30/.04	-	60036138 60036139
HyFlex EDM Finishing File, 3 pcs	21 mm	25 mm
Finishing File 40/.04	60022271	60019597
Finishing File 50/.03	60022272	60036598
Finishing File 60/.02	60022274	60036599

## Drying and Obturation

Size / Taper	HyFlex EDM Paper Points, 100 pcs	HyFlex EDM Guttapercha Points, 60 pcs
30/.04	60036503	60036504
40/.04	60022676	60022708
50/.03	60022741	60022709
60/.02	60022740	60022710



## Your Endo Guide

Practical. Innovative. Leading the way.

COLTENE provides simple, safe, and efficient solutions for every step of the endo process, carefully guiding the dentist through every endodontic situation with dedicated trainings and experienced endo consultants.

Our portfolio, which has been refined over decades, includes innovative consumables and sophisticated equipment that reliably cover all five treatment steps.



© COLTENE – www.coltene.com



**Imprint / Publisher**  
**Coltene/Whaledent AG**  
 Feldwiesenstrasse 20  
 CH-9450 Altstätten  
 +41 71 757 5300  
 info.ch@coltene.com

**Sales Contact**  
**Coltene/Whaledent Ltd.**  
 The President Suite  
 Kendal House, Victoria Way  
 Burgess Hill, West Sussex  
 RH15 9NF/UK  
 +44 1444 235486  
 info.uk@coltene.com

**Service Center**  
 service@coltene.com  
 00 800 11 0000 12 (UIFN)  
 Mon.-Fri. 8:00 to 17:00 MESZ

# HyFlex EDM – OGSF Sequence

Take control in your hands

**COLTENE**

**COLTENE**

# HyFlex EDM

Outstanding fracture resistant files for safer treatment

HyFlex EDM files have a unique surface and stand out by their durability and high fracture resistance. This is achieved by an innovative manufacturing process using "Electric Discharge Machining" (EDM). In addition, HyFlex files are manufactured from Nickel Titanium alloy with Controlled Memory heat treatment (CM).

The entirely unique combination of flexibility and fracture resistance makes it possible to reduce the number of HyFlex EDM files required for cleaning and shaping during root canal treatments without having to compromise preservation of the root canal anatomy. Like stainless steel files, the HyFlex EDM files can be pre-bent. This is a crucial advantage in curved root canals to avoid the creation of steps.

## Safely regenerable multiple times

Resetting behavior allows visual control after autoclaving and enables multiple use



✓ Ready for safe reuse




✗ Untwisted (discard)

HyFlex EDM NiTi have a built-in shape memory. They prevent stress during preparation by changing their spiral shape. This change can be quickly reversed by autoclaving to return the instruments to their original shape and resistance to fatigue.

Should the file fail to regain its shape after heating, it is plastically deformed, thus increasing the risk of fracture. It should not be used after visual inspection.


 **Pre-Sterile**  
Immediately usable

 **Excellent fracture resistance**  
Safer treatment thanks to the unique hardened surface

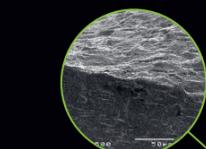

 **Powerful cutting performance**  
Allows a smaller number of files for quick preparation

 **Highly flexible and pre-bendable**  
Evenly, centered preparation reduces risk of ledging or perforation due to the controlled memory effect

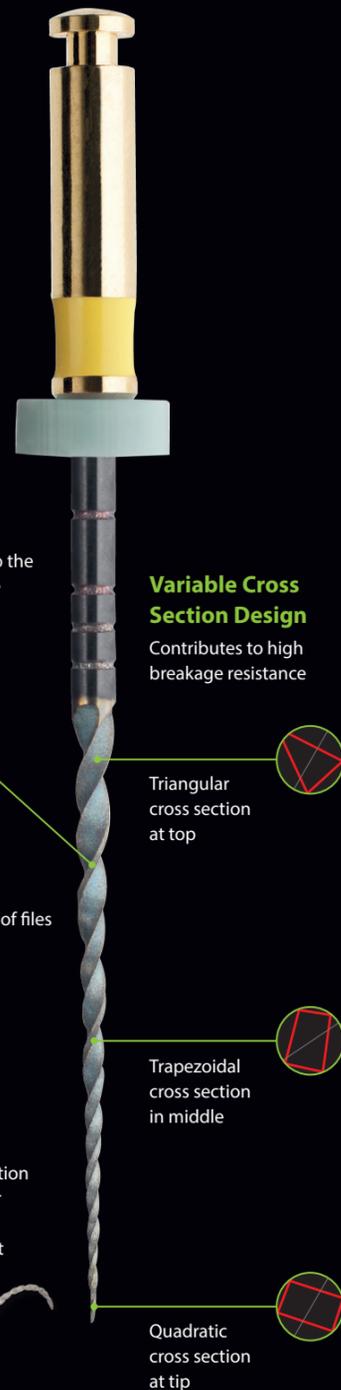

**Variable Cross Section Design**  
Contributes to high breakage resistance

Triangular cross section at top

Trapezoidal cross section in middle

Quadratic cross section at tip

HyFlex EDM Shaping File 18/.045

## OGSF Sequence

The new HyFlex EDM files – Designed to be used as a sequence

Did you say fast and safe?

Discover the new sequence solution!

The HyFlex EDM file system offers a tailor-made solution for fast and reliable preparation of simple as well as complex root canal anatomies. Evidence shows that by applying checklists and following standard procedures, reproducibility and safety are increased. Based on this awareness COLTENE offers new HyFlex EDM files designed to be used as a sequence with one, straight forward protocol.

The new HyFlex EDM OGSF sequence includes a new Orifice Opener and Glidepath file for opening the glide path, as well as a Shaping file for shaping the root canal over the full length. The defined treatment procedure makes the preparation easy and keeps the learning curve short. Due to the optimized match of tip sizes and tapers, designed in cooperation with Prof. Pedullà, the files provide a smooth feeling when shifting from one file to another and help you to keep the control over the treatment.



## Features and Benefits

- Efficient treatment with high cutting performance
- High fracture resistance and flexibility
- Centered preparation of the canal
- Only one shaping file
- Sterile blister
- 2 lengths: 21 and 25 mm