
YOTUEL® PATIENT

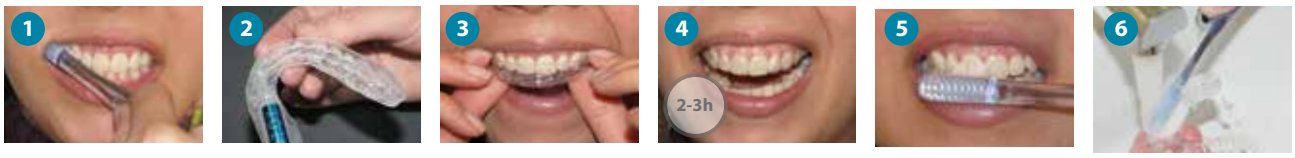ENGLISH INSTRUCTIONS

1. Brush your teeth with YOTUEL whitening toothpaste for three minutes.
2. Place a fine layer of YOTUEL whitening gel in the mouthguard.
3. Fix the mouthguard comfortably onto the teeth and wipe away any excess gel with a tissue. For the first few minutes you may notice a slight increase of saliva.
4. Keep the mouthguard on for at least 1 hour a day. For quicker results you may keep the mouthguard on for 2-3 hours or as recommended by your dentist.
5. After the whitening session brush your teeth again with YOTUEL whitening toothpaste.
6. After every session wash and dry the mouthguard.

For whiter teeth and faster results, ask your dentist for the in office whitening techniques with any of the YOTUEL dental office options.

YOTUEL WHITENING PROGRAMME

YOTUEL is a complete and advanced whitening programme for those who want to achieve perfect oral care.
YOTUEL WHITENING TOOTHPASTE With the lowest abrasion ratio worldwide (36 RDA) and a high concentration of Xylitol to fight tooth decay. It is non-foaming and detergent free.
YOTUEL WHITENING MOUTHWASH To regenerate soft tissues and enamel. Alcohol free for your total safety.

WARNINGS

Not for use by pregnant women or children.
It is not advisable to use after dental restorations.
It is not advisable to use when there are previous lesions on dental tissues.

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