



Sani-Tip®

Separating the Elements

Disposable 3-in-1 syringe tip

- ◆ Guaranteed dry air
- ◆ Separate water channel
- ◆ Improved infection control

For better dentistry

DENTSPLY

Why Sani-Tip®?



Sani-Tip®

All you need to start using the Sani-Tip® system is a small, simple-to-fit adapter, supplied FREE from DENTSPLY.

Adapters available for all major chair types, including those from Adec, Belmont, Castellini, Kavo and Siemens.

Sani-Tip is the leading disposable 3-in-1 syringe tip in the UK.

Guaranteed dry air

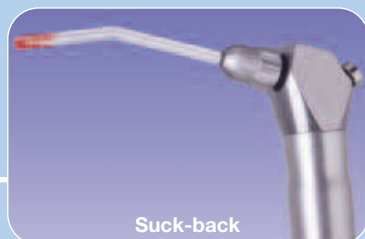
- The patented adapter design guarantees dry air every time.
- Separate air and water channels also help to keep the air dry.

Excellent cross-infection control


- There are no suck-back concerns with a disposable tip.
- The free Sani-Shields® included in every pack, give additional barrier protection for the syringe. This helps extend the syringe life by reducing exposure to the harsh disinfectants used for sterilisation.

Universal practicality


- The clear plastic tip doesn't obstruct your view.
- They are quick to install and change.
- They can be bent to over 140° without loss of function.
- An adapter is available for the majority of three-in-one syringes, past and present.
- Millions of SaniTips® have already been used by thousands of dentists across the UK.



Suck-back

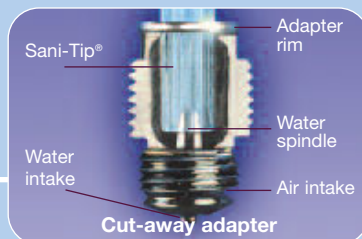

Cut-away adapter

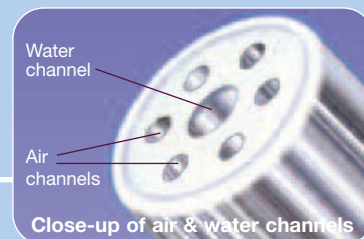

Close-up of air & water channels

Metal tips

Metal tips can cause problems with cross-infection

Published papers have shown that 'suck-back' of contaminated aerosols and oral fluids into the air and water syringe, occurs when the syringe buttons are released.¹

Although the exterior of the tip may be smoothly polished, the internal surface is rougher, making cleaning the internal surface impossible. As patient debris and body fluids must be removed before sterilising, a metal tip can never be effectively sterilised, increasing the risk of patient-to-patient cross-contamination.^{2,3}

Metal tips can become clogged

Metal tips can become clogged, as the internal surfaces corrode during use.

Metal tips do not guarantee dry air

The design of metal tips, and their repeated insertion and removal, can destroy the O-ring inside the syringe, which separates the air and water. This can lead to water leaking into the air supply, compromising the success of restorative procedures.

**To convert to SaniTip®, call DENTSPLY
on 01932 837317 or complete the attached
freepost reply card.**

REFERENCES

1. Covington J et al Quintessence International 12/89, University of Tenn
2. Miller, C. (1993) Cleaning Sterilization and Disinfection: Basics of Microbial Killing for Infection Control JADA Vol 124 1/93
3. Miller C. et al. (1994) Infection Control and Management of Hazardous Material for the Dental Team, 1st Ed
4. SDM market share data Q3 2003

Introducing Disposa-Shield® - have you got it covered?

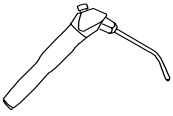
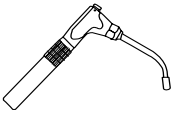
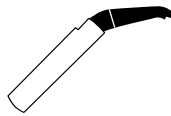
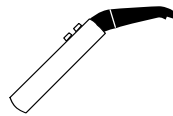
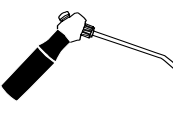
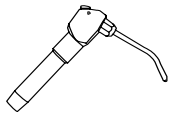
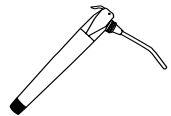
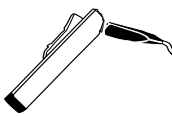
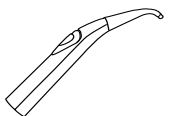
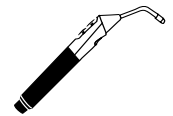
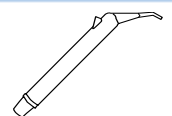
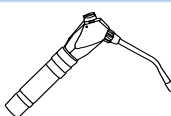
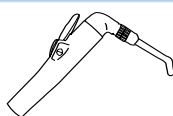
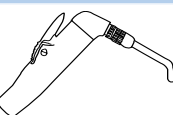
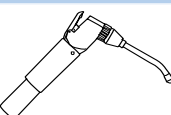
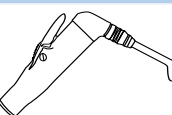
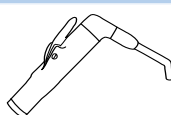
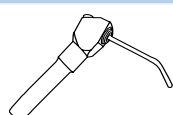
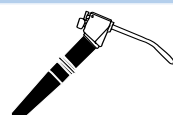
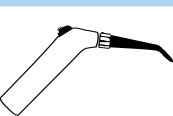
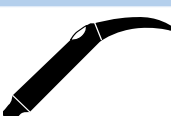
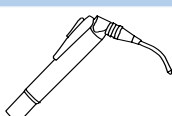
DENTSPLY's Disposa-Shield is designed to fit over surfaces to prevent cross-contamination in every practice. These protective barriers have been designed to cover a variety of shapes from instruments to chairs, including hoses, light switches, handles and handpieces.

Furthermore, they are easy to dispense and easy to change.



Sani-Tip® conversion kit reference guide

For the most popular syringes*

 A-1 122220N S1 ADEC Style Syringes	 B-1 122290N S1 BELMONT Takara Syringes	 C-1 122370N S1 CASTELLINI Cleanthree and Cleanwarm Syringes	 C-2 122371N S1 CASTELLINI Multisteril Syringes	 C-3 122372N S1 CASTELLINI Multisteril II (Puma Unit) Syringes	 D-1 122240N S1 DEN-TAL-EZ Syringes
 E-1 122230N S1 ENGLE/MARCO Syringes	 F-1 122260N S1 FARO Syringes Pre-1991	 F-2 122262N S1 FARO Electronique Syringes	 F-3 122263N S1 FARO SR903 / SR906 Syringes	 K-2 122274N S1 KAVO 573 0600 / 573 0610 Syringes (with buttons)	 K-3 122275N S1 KAVO 573 0620 / 573 0630 Syringes (Systematica no optics)
 K-4 122276N S1 KAVO 773 0500 / 773 0510 Estetica Syringes	 L-2 122325N S1 LUZZANI Minimatic Syringes (Black Plastic)	 L-3 122326N S1 LUZZANI Minilight, Minimatic, Minimate Anthos A4 Syringes	 MC-1 122340N S1 MICRON / NSK Syringes	 MC-1 122340N S1 YOSHIDA 250 Model Syringes	 M-1 122330N S1 MORITA WS9 Syringes (silver)
 M-1 122330N S1 MORITA WS10 Syringes (reddish brown)	 M-2 122336N S1 MORITA WS66 / WS88 / WS89 / WS90 Syringes	 M-3 122331N S1 MORITA WS11 Syringes (grey plastic body)	 M-4 122332N S1 MORITA WS12 Syringes (dark brown autoclavable plastic body)	 P-1 122280N S1 PRESS RING DCI Syringes	 SA-1 122350N S1 SANDRI Steritip Syringes
 S-1 122250N S1 SIEMANS Sprayvit 4000 / 4000L / 4000E Syringes	 SL-2 122251N S1 SIEMENS Sirona Syringes	 T-1 122360N S1 TRIDAC Mark 3 (MK3) Syringes (blue/silver)			

*Based on SDM Market Share Data 2006