


Sani-Tip

By using the Sani-Tip System,
you will achieve optimal bonding
results. Each Sani-Tip disposable
tip features a central water
channel encircled by six separate
air lumens. This innovative design
when matched with the appropriate
Sani-Tip adaptor, provides a
self-sealing water channel
which eliminates moisture
contamination and the chances
of adhesive failures.

DENTSPLY
PROFESSIONAL
Better products. Better practice.

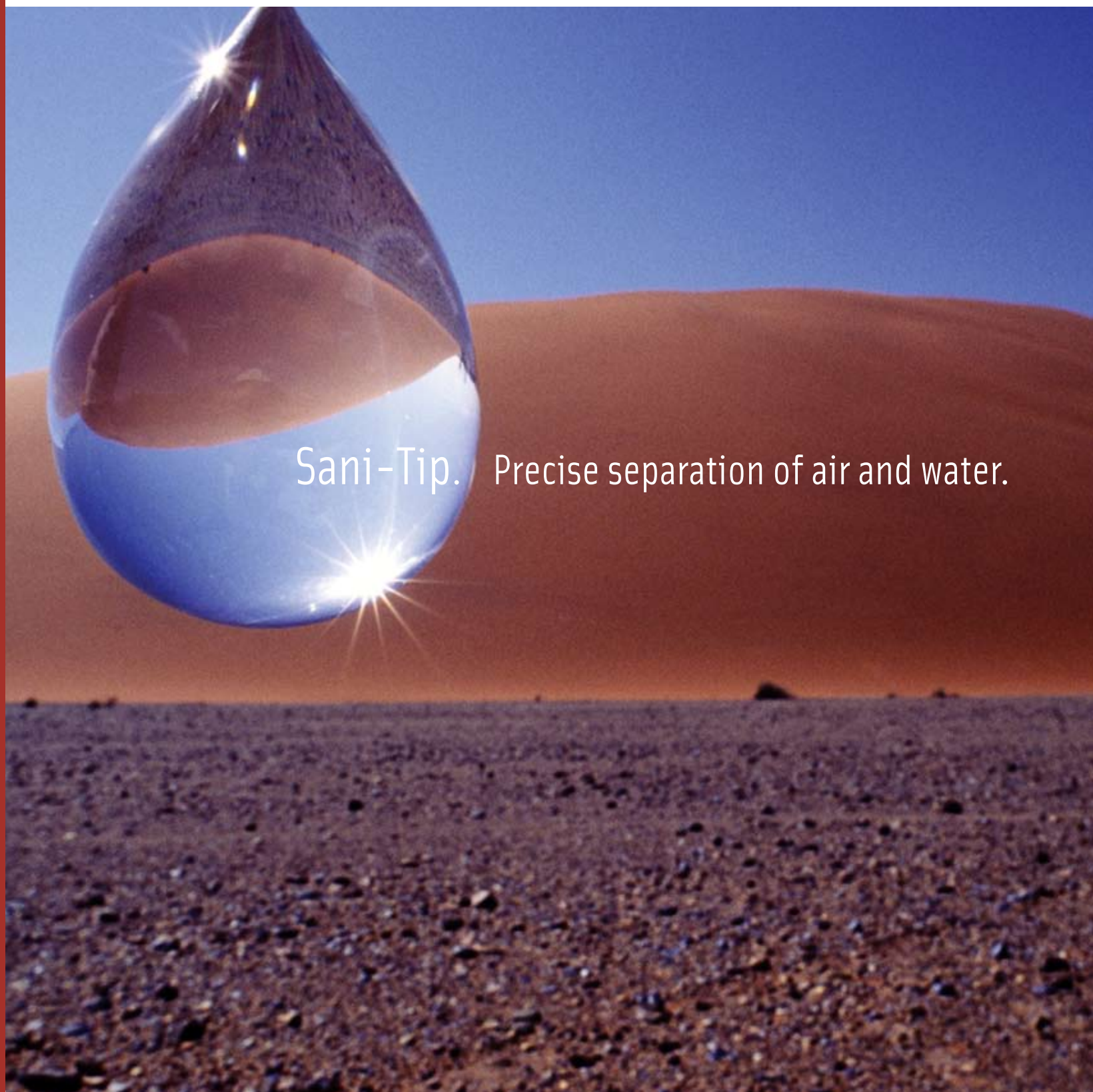

Sani-Tip. Precise separation of air and water.

Sani-Tip® Disposable Air / Water Syringe Tips. Dry air guaranteed.


DENTSPLY
PROFESSIONAL
Better products. Better practice.

Sani-Tip® Disposable Air / Water Syringe Tips

The Sani-Tip System is the surest, safest, and most cost effective way to eliminate the air/water syringe tip as a means of cross contamination.

Innovative Design

Sani-Tip's central water channel and six separate air lumens offer precise separation of air and water. This innovative design, when matched with the appropriate Sani-Tip adaptor, provides a self-sealing water channel which guarantees tips will not leak.

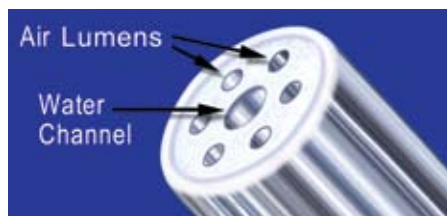

Dry Air Guaranteed

The key to the Sani-Tip System is a self-sealing water channel which guarantees dry air when a tip is fully seated over the patented tapered water spindle.



Benefits of the Sani-Tip System

Disposable

- Prevents cross-contamination.
- No clogging from internal corrosion.
- Cost effective.

Transparent

- Visually verify sanitary condition.
- Does not block view of operating field.

Superior Performance

- Self-sealing water channel guarantees dry air under optimal pressure.
- Delivers precise separation or atomization of air and water.
- Reduces adhesive restoration failure caused by moisture crossover into air lines.
- Consistent results with bonding.
- Provides secure retraction.

METAL TIP CROSS-SECTIONS



New Tip

Two Months

Six Months


One Year

Metal tip cross-sections showing a progression of corrosion and buildup on the internal surfaces of each lumen. Each used tip was processed in an autoclave four times per day. The top row shows the water channels which fit inside the air lumens in the bottom row.

From the makers of

Midwest | Nupro | Cavitron | Delton | Sani-Tip | ALLSolutions

ORDERING INFORMATION

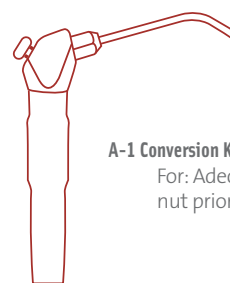

A-1 Conversion Kit 122220N5


For: Adec (all Adec styles), DC1 (thumb nut prior of 1992), Marus, Forest Medical



D-1 Conversion Kit 122240N5


For: Belmont (American), Den-Tal-Ez (Some older models use A-1)

Important: If Den-Tal-Ez unit is manufactured after 1/97, the right angle syringe requires an A-1 Conversion Kit and the straight Continental syringe requires a P-1 Conversion Kit.

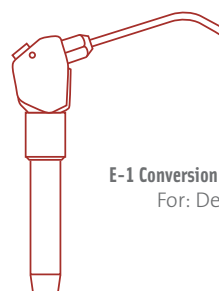

E-1 Conversion Kit 122230N5

For: Dentsply, Marco, Engle

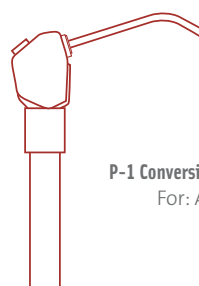

P-1 Conversion Kit 122280N5

For: All Press Ring Collar, DCI (press ring)



S-1 Conversion Kit 122250N5

For: Siemens Sprayvit 4000

REGULAR PACKS OF 250 TIPS

1222105

SANI-TIP Disposable AWS Tip-Regular Pack (250)
Includes 1 pack of 250 Tips & 250 Sani-Shields

122218-5

SANI-TIP MINI Disposable AWS Tip-Regular Pack (250)
Includes 1 pack of 250 Tips & 250 Sani-Shields

BULK PACKS OF 1500 TIPS

122224-5

SANI-TIP Disposable AWS Tip-Bulk Pack (1500)
Includes 1 pack of 1500 Tips & 1500 Sani-Shields

1222155

SANI-TIP Disposable AWS Tip-Bulk Pack (1500)
Includes 1 pack of 1500 Tips & 250 Sani-Shields



For air/water syringes not shown, contact DENTSPLY Professional Customer Service at 1.800.989.8826.



© 2004 Dentsply International. Printed in U.S.A. Form #79873 Rev. A 11/04
All rights reserved. U.S. and International Patents Pending